

Contents

1. Description
 - 1.1 Principle of the MACS® Separation
 - 1.2 Technical specifications of LS Columns
 - 1.3 Background information
 - 1.4 Applications
 - 1.5 Reagent and instrument requirements
2. Protocol
 - 2.1 Preparation of LRSCs
 - 2.2 Magnetic labeling
 - 2.3 Magnetic separation
3. Example of a separation using the StraightFrom™ LRSC CD14 MicroBead Kit

Warnings

Reagents contain sodium azide. Under acidic conditions sodium azide yields hydrazoic acid, which is extremely toxic. Azide compounds should be diluted with running water before discarding. These precautions are recommended to avoid deposits in plumbing where explosive conditions may develop.

1. Description

This product is for research use only.

Components	2 mL StraightFrom™ LRSC CD14 MicroBeads, human: MicroBeads conjugated to monoclonal anti-human CD14 antibodies (isotype: mouse IgG1). 50 mL Whole Blood Column Elution Buffer 12 LS Columns columns and plungers, sterile packed.
Capacity	For one LRSC.
Product format	StraightFrom LRSC CD14 MicroBeads are supplied in buffer containing stabilizer and 0.05% sodium azide. Whole Blood Column Elution Buffer contains stabilizer and 0.09% sodium azide.
Storage	Store StraightFrom LRSC CD14 MicroBeads and Whole Blood Column Elution Buffer protected from light at 2–8 °C. Do not freeze. Store LS Columns dry at 10–35 °C and protected from light. The expiration date is indicated on the vial or box label.

1.1 Principle of the MACS® Separation

First, the CD14⁺ cells in an LRSC sample are magnetically labeled with StraightFrom LRSC CD14 MicroBeads. Then, the cell suspension is loaded onto an LS Column, which is placed in the magnetic field of a MACS® Separator. The magnetically labeled CD14⁺ cells are retained within the column. The unlabeled cells run through; this cell fraction is thus depleted of CD14⁺ cells. After removing the column from the magnetic field, the magnetically retained CD14⁺ cells can be eluted as the positively selected cell fraction.

1.2 Technical specifications of LS Columns

- Columns are “flow stop” and do not run dry.
- Void volume: 400 µL. Reservoir volume: 8 mL.
- LS Columns are for single use only.
- Use Whole Blood Column Elution Buffer to elute positive cell fraction from column.

▲ To remove clumps and to prevent aggregates in the sample, pass cells through 30 µm nylon mesh (MACS SmartStrainers (30 µm), # 130-098-458) before separation.

1.3 Background information

Leukoreduction system chambers (LRSCs), also known as buffy cones, arise after standard plateletpheresis procedures, where leukocytes are filtered out into the LRSCs. They contain highly concentrated leukocytes in a small volume.

The StraightFrom LRSC CD14 MicroBead Kit has been developed for the positive selection of CD14⁺ cells directly from LRSCs by using the MultiMACS™ Cell24 Separator Plus or the MidiMACS™ Separator or QuadroMACS™ Separator. No sample preparation is required, including density gradient centrifugation or erythrocyte lysis. The CD14 antigen belongs to the LPS receptor complex. It is strongly expressed on most monocytes and macrophages and weakly on neutrophils.

1.4 Applications

- Isolation of CD14⁺ cells from LRSCs. The purified CD14⁺ cells are well suited for further flow cytometric, functional, or molecular analysis.

1.5 Reagent and instrument requirements

- Separation buffer: Prepare a solution containing phosphate-buffered saline (PBS), pH 7.2, 0.5% bovine serum albumin (BSA), and 2 mM EDTA by diluting MACS BSA Stock Solution (# 130-091-376) 1:20 with autoMACS® Rinsing Solution (# 130-091-222). Keep buffer cold (2–8 °C). Alternatively, use autoMACS Running Buffer. Degas buffer before use, as air bubbles could block the column.

▲ **Note:** BSA can be replaced by other proteins such as human serum albumin, human serum, or fetal bovine serum (FBS). Buffers or media containing Ca²⁺ or Mg²⁺ are not recommended for use.

- Semi-automated separation:
 - MultiMACS™ Cell24 Separator Plus (# 130-098-637) with the MultiMACS 12× Single-Column Adapter LS (# 130-108-816)
 - Single-well Deep Well Plates (# 130-114-966) or 24-well Deep Well Plates (# 130-110-500)
- Manual separation:
 - MidiMACS™ Separator (# 130-042-302) or QuadroMACS™ Separator (# 130-090-976)
 - MACS® Acrylic Tube Rack (# 130-041-406) or MACS 15 mL Tube Rack (# 130-091-052)
- (Optional) Fluorochrome-conjugated antibodies for flow cytometric analysis, e.g., CD14-PE, CD15-APC, and CD45-VioBlue® or the MC CD14 Monocyte Cocktail, human. For more information about antibodies refer to www.miltenyibiotec.com/antibodies.
- (Optional) Propidium Iodide Solution (# 130-093-233) or 7-AAD Staining Solution (# 130-111-568) for flow cytometric exclusion of dead cells.
- (Optional) MACS SmartStrainers (30 µm) (# 130-098-458) to remove cell clumps.

2. Protocol

2.1 Preparation of LRSCs

1. Transfer content of one LRSC into a collection tube. Fill up to 40 mL with separation buffer.

▲ **Note:** LRSC can be rinsed with separation buffer to maximize cell recovery from the LRSC. Do not exceed 40 mL sample volume in total.
2. Proceed to magnetic labeling (2.2).



2.2 Magnetic labeling

▲ Work fast, keep cells cold, and use pre-cooled solutions. This will prevent capping of antibodies on the cell surface and non-specific cell labeling.

▲ Volumes for magnetic labeling given below are for one entire LRSC (filled up to 40 mL sample volume). When working with fractions of LRSC, scale down all reagents and total volumes accordingly (e.g. for half an LRSC, fill up to 20 mL and use 1 mL StraightFrom™ LRSC MicroBeads as well as 6 LS Columns).

▲ For optimal performance it is important to obtain a single-cell suspension before magnetic labeling. Pass cells through 30 µm nylon mesh (MACS SmartStrainers (30 µm), # 130-098-458) to remove cell clumps which may clog the column. Moisten filter with separation buffer before use.

▲ The recommended incubation temperature is 2–8 °C. Higher temperatures and/or longer incubation times may lead to non-specific cell labeling. Working on ice may require increased incubation times.

1. Add 2 mL StraightFrom LRSC CD14 MicroBeads per 40 mL sample volume.
2. Mix well by inverting the tube and incubate for 15 minutes in the refrigerator (2–8 °C).
3. Proceed directly to magnetic separation (2.3).



2.3 Magnetic separation

2.3.1 Semi-automated magnetic separation with the MultiMACS™ Cell24 Separator Plus

▲ For more detailed instructions on how to use the MultiMACS™ Cell24 Separator Plus, please refer to the user manual.

▲ The MultiMACS Cell24 Separator Plus, including the MACS Elution Station, has to be used with the Single-Column Adapter, LS Columns, and Deep Well Plates for magnetic separation with StraightFrom LRSC MicroBeads.

▲ **Note:** To reach maximum cell recovery, rinse the Deep Well Plate of the positive fraction and combine with the positive fraction.

▲ Buffer volumes are as follows:

Equilibration: 3 mL (separation buffer)

Wash: 2×2 mL (separation buffer)

Elution: 4 mL (Whole Blood Column Elution Buffer)

▲ Divide the sample equally between the 12 LS Columns, e.g., when starting volume is 42 mL (40 mL sample and 2 mL StraightFrom LRSC MicroBeads), add 3.5 mL onto each column.

▲ If equilibration solution and negative fraction should be collected in the same Deep Well Plate as waste, select the program **POSSEL_SCA** and follow the on-screen instructions of the MultiMACS Cell24 Separator Plus.

▲ If equilibration solution and negative fraction with wash fractions should be collected in separate Deep Well Plates, select the program **POSSEL2_SCA** and follow the on-screen instructions of the MultiMACS Cell24 Separator Plus.

▲ After the separation, centrifuge positive fraction at 300×g for 10 minutes. Aspirate supernatant carefully. Resuspend cell pellet in a suitable amount of buffer or medium for subsequent analysis.

2.3.2 Manual magnetic separation with the MidiMACS™ Separator or QuadroMACS™ Separator

▲ The protocol given below is for one LS Column. Repeat steps for each of the 12 LS Columns.

▲ Divide the sample equally between the LS Columns, e.g., when starting volume is 42 mL, add 3.5 mL onto each column.


1. Place LS Column in the magnetic field of a MidiMACS™ Separator or QuadroMACS™ Separator. Insert LS Column with the column wings to the front into the separator.
2. Prepare column by rinsing with 3 mL separation buffer.
3. Apply magnetically labeled cell suspension onto the prepared LS Column. Collect flow-through containing unlabeled cells.
4. Wash LS Column with 2×2 mL separation buffer. Collect unlabeled cells that pass through and combine with the flow-through from step 3.

▲ **Note:** Perform washing steps by adding buffer aliquots as soon as the column reservoir is empty.

5. Remove LS Column from the separator and place it on a new collection tube.
6. Pipette 4 mL **Whole Blood Column Elution Buffer** onto the LS Column. Immediately flush out the magnetically labeled cells by firmly pushing the plunger into the column.
7. Centrifuge cells at 300×g for 10 minutes. Aspirate supernatant carefully. Resuspend cell pellet in a suitable amount of buffer or medium for subsequent analysis.

3. Example of a separation using the StraightFrom™ LRSC CD14 MicroBead Kit

Separation of an LRSC sample using the StraightFrom™ LRSC CD14 MicroBead Kit and the MultiMACS™ Cell24 Separator Plus with the Single-Column Adapter and LS Columns. Cells were fluorescently stained with CD14-PE, CD15-APC, as well as CD45-VioBlue® and analyzed by flow cytometry using the MACSQuant® Analyzer. Cell debris and dead cells were excluded from the analysis based on scatter signals and propidium iodide fluorescence.

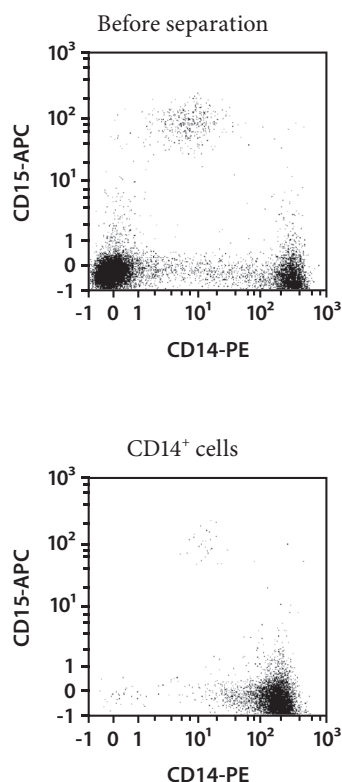

Refer to www.miltenyibiotec.com for all data sheets and protocols. Miltenyi Biotec provides technical support worldwide. Visit www.miltenyibiotec.com/local to find your nearest Miltenyi Biotec contact.

Legal notices

Limited product warranty

Miltenyi Biotec B.V. & Co. KG and/or its affiliate(s) warrant this product to be free from material defects in workmanship and materials and to conform substantially with Miltenyi Biotec's published specifications for the product at the time of order, under normal use and conditions in accordance with its applicable documentation, for a period beginning on the date of delivery of the product by Miltenyi Biotec or its authorized distributor and ending on the expiration date of the product's applicable shelf life stated on the product label, packaging or documentation (as applicable) or, in the absence thereof, ONE (1) YEAR from date of delivery ("Product Warranty"). Miltenyi Biotec's Product Warranty is provided subject to the warranty terms as set forth in Miltenyi Biotec's General Terms and Conditions for the Sale of Products and Services available on Miltenyi Biotec's website at www.miltenyibiotec.com, as in effect at the time of order ("Product Warranty"). Additional terms may apply. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING IF A PRODUCT IS SUITABLE FOR CUSTOMER'S PARTICULAR PURPOSE AND APPLICATION METHODS.

Technical information

The technical information, data, protocols, and other statements provided by Miltenyi Biotec in this document are based on information, tests, or experience which Miltenyi Biotec believes to be reliable, but the accuracy or completeness of such information is not guaranteed. Such technical information and data are intended for persons with knowledge and technical skills sufficient to assess and apply their own informed judgment to the information. Miltenyi Biotec shall not be liable for any technical or editorial errors or omissions contained herein.

All information and specifications are subject to change without prior notice. Please contact Miltenyi Biotec Technical Support or visit www.miltenyibiotec.com for the most up-to-date information on Miltenyi Biotec products.

Licenses

This product and/or its use may be covered by one or more pending or issued patents and/or may have certain limitations. Certain uses may be excluded by separate terms and conditions. Please contact your local Miltenyi Biotec representative or visit Miltenyi Biotec's website at www.miltenyibiotec.com for more information.

The purchase of this product conveys to the customer the non-transferable right to use the purchased amount of the product in research conducted by the customer (whether the customer is an academic or for-profit entity). This product may not be further sold. Additional terms and conditions (including the terms of a Limited Use Label License) may apply.

CUSTOMER'S USE OF THIS PRODUCT MAY REQUIRE ADDITIONAL LICENSES DEPENDING ON THE SPECIFIC APPLICATION. THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING FOR ITSELF WHETHER IT HAS ALL APPROPRIATE LICENSES IN PLACE. Miltenyi Biotec provides no warranty that customer's use of this product does not and will not infringe intellectual property rights owned by a third party. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

Trademarks

autoMACS, MACS, MACSQuant, MidiMACS, the Miltenyi Biotec logo, MultiMACS, QuadroMACS, StraightFrom, and VioBlue are registered trademarks or trademarks of Miltenyi Biotec and/or its affiliates in various countries worldwide. All other trademarks mentioned in this publication are the property of their respective owners and are used for identification purposes only.

Leukopak is a registered trademark of StemExpress, LCC.

Copyright © 2020 Miltenyi Biotec and/or its affiliates. All rights reserved.