



**Miltenyi Biotec**

# ImSpector 8.0

User manual



Copyright © 2025 Miltenyi Biotec and/or its affiliates. All rights reserved.

No part of this user manual may be reproduced, stored in a retrieval system, transmitted, published, or distributed in any form or by any means, electronically, mechanically, by photocopying, microfilming, recording, or otherwise, without the prior written consent of Miltenyi Biotec; however, notwithstanding the foregoing, the owners of the Imaging System may make copies solely for purposes of training personnel in the use and servicing of the unit within their business or organization. Maximal care has been taken by Miltenyi Biotec in the preparation of this user manual. However, Miltenyi Biotec shall not be liable for any technical or editorial errors or omissions contained herein, or for incidental or consequential damages in connection with the furnishing, performance, or use of this document. The information in this document is provided "as is" without warranty of any kind, and is subject to change without notice. Unless otherwise specifically indicated, Miltenyi Biotec products and services are for research use only and not for therapeutic or diagnostic use.



# ImSpector 8.0

## User manual

Original instructions

Issued: 2025-07

140-006-592-007.07

Miltenyi Biotec B.V. & Co. KG

Friedrich-Ebert-Straße 68

51429 Bergisch Gladbach

Germany

Phone +49 2204 8306-0

Fax +49 2204 85197

[macsde@miltenyi.com](mailto:macsde@miltenyi.com)

[www.miltenyibiotec.com](http://www.miltenyibiotec.com)





# Contents

|                                      |           |
|--------------------------------------|-----------|
| <b>1 About this user manual</b>      | <b>7</b>  |
| 1.1 Read the user manual first       | 7         |
| <b>2 Introduction</b>                | <b>9</b>  |
| 2.1 The ImSpector software           | 9         |
| 2.2 Intended use                     | 9         |
| <b>3 Basic functions</b>             | <b>11</b> |
| 3.1 Switching on the instrument      | 11        |
| 3.2 Switching off the instrument     | 11        |
| <b>4 Graphical user interface</b>    | <b>13</b> |
| 4.1 Menu bar                         | 13        |
| 4.1.1 File menu                      | 14        |
| 4.1.2 View menu                      | 14        |
| 4.1.3 Analysis and Processing menu   | 15        |
| 4.1.4 Settings menu                  | 15        |
| 4.1.5 Help menu                      | 16        |
| 4.1.6 User menu                      | 16        |
| 4.2 Toolbar                          | 17        |
| <b>5 Live Dialogs</b>                | <b>19</b> |
| 5.1 Rearranging the dialogs          | 19        |
| 5.2 Using the dialogs                | 20        |
| 5.2.1 Channel Setup                  | 20        |
| 5.2.2 Light Sheet Properties         | 21        |
| 5.2.3 Objective Lens Focus           | 21        |
| 5.2.4 Optical Setup                  | 22        |
| 5.2.5 Light Sheet Focus              | 22        |
| 5.2.6 Light Sheet Selection          | 22        |
| 5.2.7 Calibration tab for excitation | 23        |
| 5.2.8 Advanced tab                   | 23        |

|                                                      |           |
|------------------------------------------------------|-----------|
| 5.2.9 Acquisition Mode tab .....                     | 24        |
| 5.2.10 xyz-Table Z view .....                        | 26        |
| 5.2.11 xyz-Table Visual XY view .....                | 27        |
| 5.2.12 Camera .....                                  | 30        |
| <b>6 Setting up an experiment .....</b>              | <b>33</b> |
| <b>6.1 Measurement Wizard .....</b>                  | <b>33</b> |
| <b>6.2 Instrument Mode .....</b>                     | <b>34</b> |
| <b>6.3 Measurement Mode .....</b>                    | <b>34</b> |
| 6.3.1 Acquiring an image .....                       | 35        |
| 6.3.2 Acquiring a multicolor z-stack .....           | 36        |
| 6.3.3 Acquiring a 3D mosaic image .....              | 37        |
| 6.3.4 Acquiring a multicolor mosaic image .....      | 38        |
| 6.3.5 Acquiring a 3D multiposition dataset .....     | 40        |
| <b>6.4 Experiment Manager .....</b>                  | <b>42</b> |
| <b>6.5 Setting up a batch measurement .....</b>      | <b>43</b> |
| <b>7 Data visualization and analysis .....</b>       | <b>45</b> |
| <b>7.1 Stack window .....</b>                        | <b>45</b> |
| <b>7.2 Image Series Viewer .....</b>                 | <b>49</b> |
| 7.2.1 Image Series Browser .....                     | 49        |
| 7.2.2 Image Series Manager .....                     | 50        |
| 7.2.3 Image Series Axes .....                        | 50        |
| <b>7.3 Stack Properties .....</b>                    | <b>50</b> |
| <b>7.4 Changing the view .....</b>                   | <b>50</b> |
| 7.4.1 Colormap dialog .....                          | 51        |
| 7.4.2 Simple view .....                              | 51        |
| 7.4.3 Advanced view .....                            | 54        |
| <b>7.5 Merging and splitting .....</b>               | <b>55</b> |
| <b>7.6 Defining a region of interest (ROI) .....</b> | <b>55</b> |
| <b>7.7 Showing measurement information .....</b>     | <b>56</b> |
| <b>8 Data management .....</b>                       | <b>57</b> |

|                                        |           |
|----------------------------------------|-----------|
| <b>8.1 Types of windows</b>            | <b>57</b> |
| <b>8.2 Managing files</b>              | <b>57</b> |
| 8.2.1 New documents                    | 59        |
| 8.2.2 New graph or stack windows       | 59        |
| 8.2.3 Open documents                   | 60        |
| 8.2.4 Save documents                   | 60        |
| <b>9 Settings</b>                      | <b>63</b> |
| <b>9.1 User settings</b>               | <b>63</b> |
| <b>9.2 Hardware settings</b>           | <b>63</b> |
| <b>10 Important keyboard shortcuts</b> | <b>65</b> |
| <b>11 Technical support</b>            | <b>67</b> |
| <b>12 Legal notes</b>                  | <b>69</b> |
| <b>12.1 Limited warranty</b>           | <b>69</b> |
| <b>12.2 Trademarks</b>                 | <b>70</b> |





# 1

## About this user manual

### 1.1 Read the user manual first

#### **WARNING**

Before using the ImSpector software, read the chapter **Important safety information** in the instrument user manual and all other information contained in this software user manual, including any safety and operating instructions. Pay special attention to all warnings displayed on the instrument. Failure to read and follow these guidelines could lead to improper or incorrect usage and result in damage to the instrument. Improper usage could also cause severe personal injury, death, unpredictable results, instrument malfunction, and premature wear to components shortening the lifetime of the instrument. Such actions may void your warranty. Keep the user manual and any other safety and operating instructions provided with the instrument in a safe place accessible to all users for future reference.

If you have a serious concern regarding the safe use of your instrument, contact your authorized Miltenyi Biotec service provider or call Miltenyi Biotec Technical Support.



# 2

## Introduction

### 2.1 The ImSpector software

The ImSpector software is a Windows program for setting up the image acquisition with the UltraMicroscope Blaze Imaging System and the UltraMicroscope Choros Imaging System. The software also provides data viewing and basic analytic features. The integrated Python script package allows access to imaging data and microscope devices. Python can be used for advanced measurement setups or analysis of imaging data.

### 2.2 Intended use

The ImSpector software in this version is configured for the use with the UltraMicroscope Blaze Imaging System and the UltraMicroscope Choros Imaging System. The software controls the instrument and generates high-resolution images.

The ImSpector software is for research use only in a laboratory environment. Only laboratory professional users may use the software.



# 3

## Basic functions

This software user manual applies to the UltraMicroscope Blaze Imaging System and the UltraMicroscope Choros Imaging System running ImSpector software version 8.0 or later.

### 3.1 Switching on the instrument

- 1 Ensure that all components are connected to the power supply.
- 2 Switch on the external laser module. The key switch is located on the laser module. Refer to the manufacturer documentation.
- 3 Switch on the instrument. The power switch is located on the rear side of the instrument. Refer to the instrument user manual.
- 4 Switch on the acquisition workstation, and sign in to the Windows system.
- 5 Open the ImSpector software. The switched on devices are registered and checked.

### 3.2 Switching off the instrument

- 1 Close the ImSpector software.
- 2 Shut down the Windows system.
- 3 Switch off the instrument.
- 4 Switch off the external laser module.






# 4

## Graphical user interface

This chapter describes the graphical user interface. Most functions can be accessed via the menu bar (1). Frequently used functions can also be accessed directly from the toolbar (2).

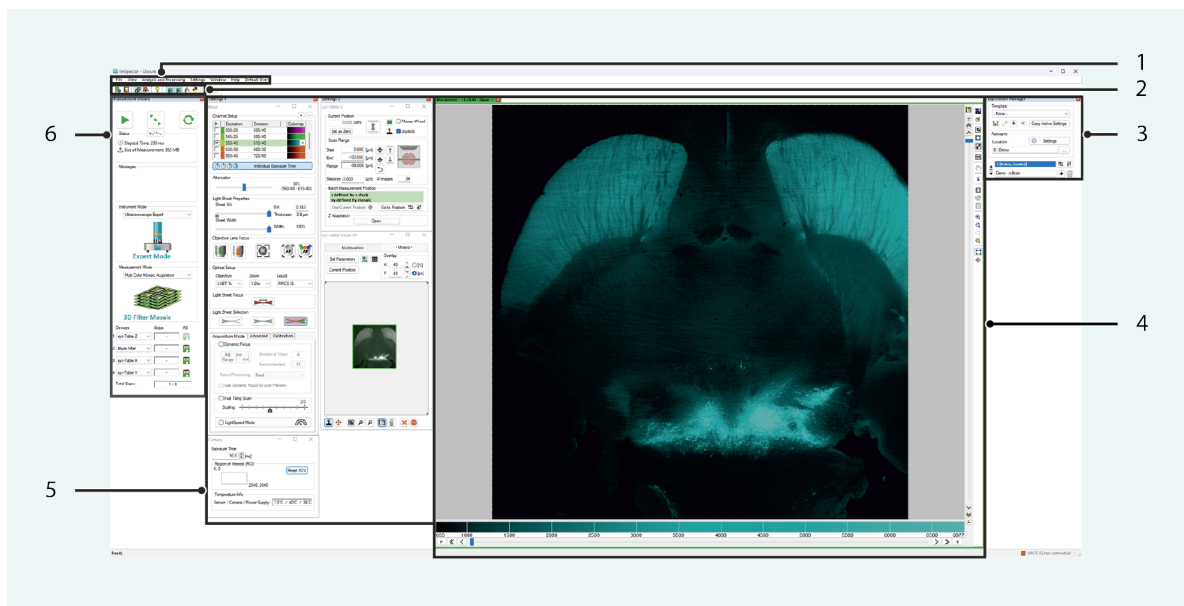


**Figure 4.1:** Graphical user interface

- |   |                    |   |                    |
|---|--------------------|---|--------------------|
| 1 | Menu bar           | 4 | Stack window       |
| 2 | Toolbar            | 5 | Settings pane      |
| 3 | Experiment Manager | 6 | Measurement Wizard |

### 4.1

#### Menu bar

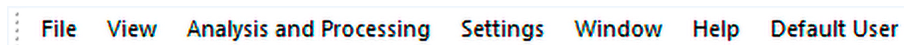

**Figure 4.2:** Menu bar

The menu bar provides access to the following menus:

- File
- View
- Analysis and Processing
- Settings
- Window
- Help
- Default User

4.1.1

File menu

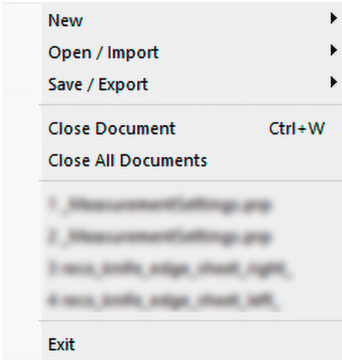


Figure 4.3: File menu

|                     |                                                                                        |
|---------------------|----------------------------------------------------------------------------------------|
| New                 | Creates a new document, a new stack window, or a graph window for the current document |
| Open / Import       | Opens or imports documents or templates                                                |
| Save / Export       | Saves or exports documents or templates                                                |
| Close Document      | Closes the currently opened document window                                            |
| Close All Documents | Closes all document windows                                                            |
| Exit                | Closes the software                                                                    |

Table 4.1: File menu

All other previously listed opened files are also shown in this menu.

4.1.2

View menu

All dialogs are listed in the **View** menu.


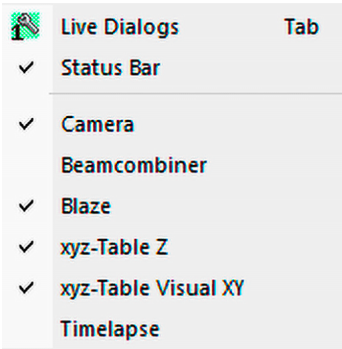


Figure 4.4: View menu

|                     |                               |
|---------------------|-------------------------------|
| <b>Live Dialogs</b> | Shows all live dialogs        |
| <b>Status Bar</b>   | Shows or hides the status bar |

Table 4.2: View menu

### 4.1.3 Analysis and Processing menu

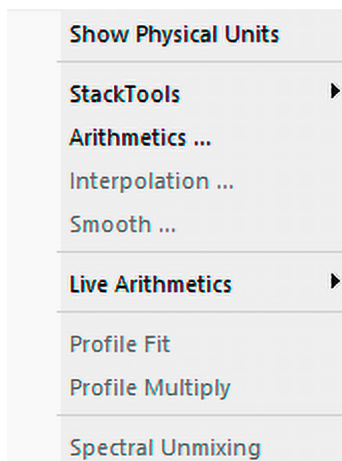

|                            |                                               |
|----------------------------|-----------------------------------------------|
| <b>Show Physical Units</b> | Shows physical units as unit of length        |
| <b>StackTools</b>          | Selects the stack tools                       |
| <b>Arithmetics</b>         | Opens the <b>Arithmetic Operations</b> dialog |
| <b>Live Arithmetics</b>    | Selects the arithmetic tools                  |
| <b>Profile Fit</b>         | <i>Not used for this instrument</i>           |
| <b>Profile Multiply</b>    | <i>Not used for this instrument</i>           |
| <b>Spectral Unmixing</b>   | <i>Not used for this instrument</i>           |

Table 4.3: Analysis and Processing menu

### 4.1.4 Settings menu

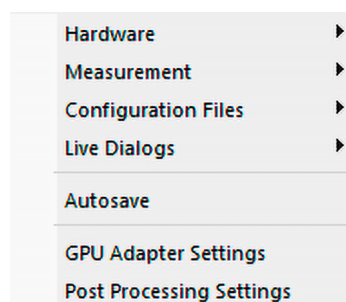



Figure 4.5: Settings menu

|                                 |                                                                                                        |
|---------------------------------|--------------------------------------------------------------------------------------------------------|
| <b>Hardware</b>                 | Opens a menu to change the hardware configuration; some changes may require support by Miltenyi Biotec |
| <b>Measurement</b>              | Opens a menu to change advanced measurement configuration settings and show or hide buttons or boxes   |
| <b>Configuration Files</b>      | Opens a menu to change configuration files; some changes may require support by Miltenyi Biotec        |
| <b>Live Dialogs</b>             | Changes the settings for visibility and position of dialogs                                            |
| <b>Autosave Settings</b>        | Changes the settings for autosave, for example, the storage location                                   |
| <b>User Settings</b>            | Changes the user settings                                                                              |
| <b>Remote Control Settings</b>  | Changes the remote control settings ( <i>only used if two microscopes need to be synchronized</i> )    |
| <b>GPU Adapter Settings</b>     | <i>Not used for this instrument</i>                                                                    |
| <b>Post Processing Settings</b> | <i>Not used for this instrument</i>                                                                    |

**Table 4.4: Settings** menu

#### 4.1.5

#### Help menu

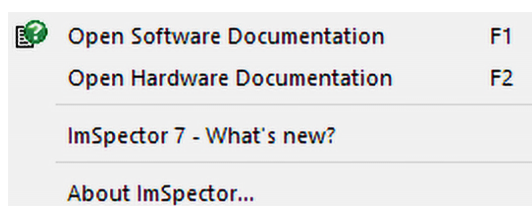


**Figure 4.6: Help** menu

|                                    |                                                  |
|------------------------------------|--------------------------------------------------|
| <b>Open Software Documentation</b> | Opens the software documentation as an HTML file |
| <b>Open Hardware Documentation</b> | Opens the hardware documentation as a PDF file   |
| <b>What's new?</b>                 | Opens the Release Notes                          |
| <b>About ImSpector...</b>          | Shows information about the software             |

**Table 4.5: Help** menu

#### 4.1.6

#### User menu

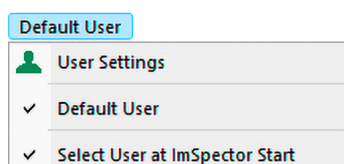

**Figure 4.7: User** menu

|                                       |                                                               |
|---------------------------------------|---------------------------------------------------------------|
| <b>User Settings</b>                  | Opens a dialog to set up user profiles                        |
| <b>Default User</b>                   | Selects the default user or changes the active user           |
| <b>Select User at ImSpector Start</b> | Activates selecting the active user when opening the software |

**Table 4.6:** User menu

## 4.2 Toolbar

The toolbar provides quick access to often used functions and features.



**Figure 4.8:** Toolbar

|                                                                                     |                            |                                                                            |
|-------------------------------------------------------------------------------------|----------------------------|----------------------------------------------------------------------------|
|    | <b>New Document</b>        | Creates a new document                                                     |
|    | <b>Image Series Viewer</b> | Shows or hides the <b>Image Series Viewer</b>                              |
|    | <b>Save as</b>             | Saves a document                                                           |
|   | <b>Dark Mode</b>           | Turns on the <b>Dark Mode</b> , which softens the screen with a dark theme |
|  | <b>Settings 1</b>          | Shows or hides the Settings 1 pane. In this case, the pane is shown.       |
|  | <b>Settings 2</b>          | Shows or hides the Settings 2 pane. In this case, the pane is shown.       |
|  | <b>Settings 3</b>          | Shows or hides the Settings 3 pane. In this case, the pane is hidden.      |
|  | <b>Python</b>              | Opens the Python script editor                                             |
|  | <b>User Settings</b>       | Opens the <b>User Settings</b> dialog                                      |
|  | <b>Save Settings</b>       | Saves the settings as default                                              |

**Table 4.7:** Buttons of the toolbar



# 5

## Live Dialogs

The **Live Dialogs** are used to configure and control various components of the instrument. Each component has its own dialog that asks for input.


### 5.1 Rearranging the dialogs

You can arrange the individual dialogs in up to three settings panes on the left-hand side of the stack window.

- 1 To show or hide each pane, click **Settings > Live Dialogs > Show / Hide**.
- 2 To dock and undock panes, double-click the corresponding title bar.
- 3 To show or hide all panes except the **Measurement Wizard**, press the Tab key.


The availability of the dialogs depends on the **Instrument Mode** and **Measurement Mode**.

- 1 Select the dialog layout you want to use:
  - To save the dialog layout for the current **Instrument Mode** and for each **Measurement Mode**, select **Settings > Live Dialogs > Keep Dialog Positions (All Measurement Modes)**.

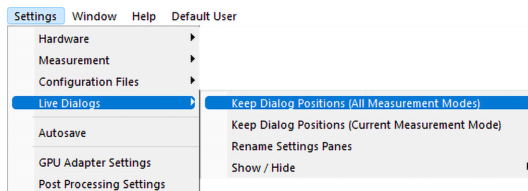

or

- To save the dialog layout for the current **Instrument Mode** and for the current **Measurement Mode**, select **Settings > Live Dialogs > Keep Dialog Positions (Current Measurement Mode)**.
- 2 Optional: To rename each pane, select **Settings > Live Dialogs > Rename Settings Panes** and enter the new name.

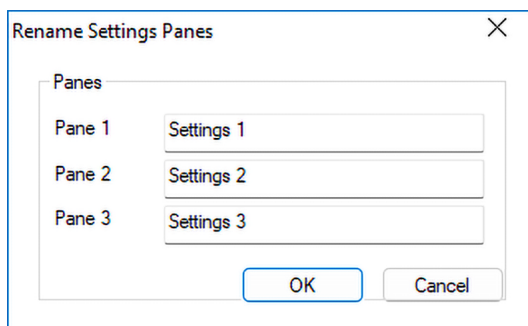

- 3 To show or hide each pane, select **Settings > Live Dialogs > Show / Hide**.

## 5.2 Using the dialogs

You can use the dialogs to control various aspects of the instrument. For use cases, refer to the instrument user manual.

### 5.2.1 Channel Setup

In the **Channel Setup**, you can find all available excitation filters and emission filters for the instrument. The list shows a checkbox followed by **Excitation** wavelength, **Emission** wavelength, and a configurable **Colormap**.

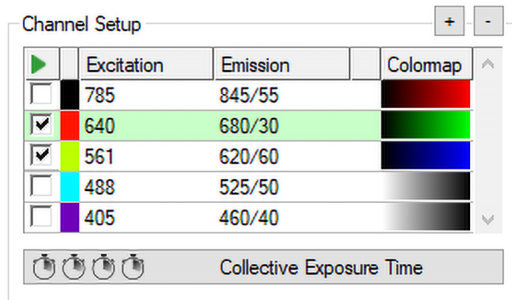


Figure 5.1: Channel Setup


To change the colormap:

- 1 Click the colormap for the respective excitation wavelength and emission wavelength.




- 2 Select the desired color.

The stack is automatically adjusted.

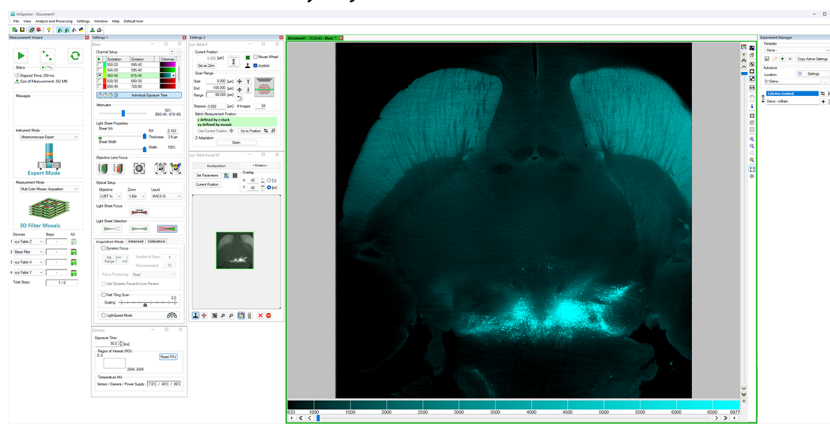

To set individual exposure times for the selected channels:


- 1 Click the **Collective Exposure Time** button. The button changes to **Individual Exposure Time**.

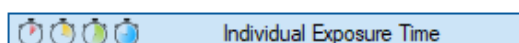



To set collective exposure times for the selected channels:

- 1 Click the **Individual Exposure Time** button. The button changes to **Collective Exposure Time**.

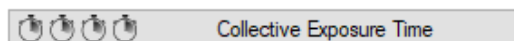

### 5.2.2 Light Sheet Properties

To change the numerical aperture (NA) of the light sheet and its width, you can use **Light Sheet Properties**.

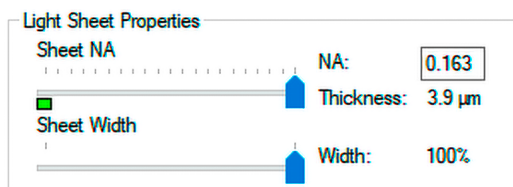


Figure 5.2: Light Sheet Properties

- 1 Move the slider under **Sheet NA** to adjust the thickness of the light sheet. A higher NA creates a thinner light sheet near the light sheet focus. A lower NA creates a more homogeneous light sheet, but reduces the resolution in the z-plane.

The shown thickness of the light sheet is only a reference value. The quality of the sample and the alignment of the instrument can also have an effect on the thickness of the light sheet.

- 2 Move the slider under **Sheet Width** to adjust sheet width of the light sheet. A lower sheet width results in more light being focused on a smaller region of the sample.

If the sheet width is reduced too much, graduated or incomplete illumination of the field of view may occur.

### 5.2.3 Objective Lens Focus



#### Move Detection Unit Down

Moves the detection unit to the measurement position; must be held clicked for safety reasons



#### Move Detection Unit Up

Moves the detection unit to the park position



#### Manual Focus

Activates slider to focus manually



#### Single Channel Autofocus

Focuses automatically on current channel. When clicked, the detection unit goes through a short path in z-direction. The software selects the z-position with the strongest contrast. Right-click causes to toggle autofocus on the ROI.



#### Multichannel Autofocus

Focuses automatically on all selected channels in the **Channel Setup**. Right-click causes to focus only on the ROI.

Table 5.1: Buttons of the **Objective Lens Focus** dialog

Autofocus is applied to the reference channel first, as it affects the focus position of all other channels. When focusing manually, the reference channel must be focused first.

## 5.2.4 Optical Setup

To select the correct settings for **Objective**, **Zoom**, and **Liquid** of the experiment, you can use **Optical Setup**.

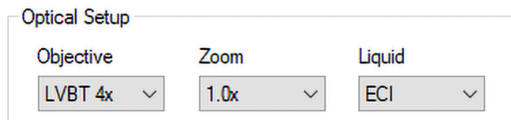


Figure 5.3: Optical Setup

The selected **Liquid** must match the liquid in the cuvette because the light sheet focus depends on the refractive index of the liquid used. Additional liquids can be added to the preset ones.

## 5.2.5 Light Sheet Focus

To toggle the visibility of the light sheet focus indicator in the stack window, you can click the **Show Light Sheet focus adjustment** button under **Light Sheet Focus**.

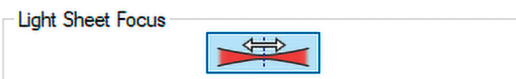


Figure 5.4: Light Sheet Focus


## 5.2.6 Light Sheet Selection

To select the optimal light sheet setting for the experiment, you can use **Light Sheet Selection**.



Figure 5.5: Light Sheet Selection

To select the left light sheet only:

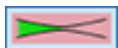

1 Click **Select left light sheet**.

To select the right light sheet only:



1 Click **Select right light sheet**.

To select both light sheets:



1 Click **Select both light sheets**.

### 5.2.7 Calibration tab for excitation

You can use the **Calibration** tab to set a reference point for the software focus position to the actual position of the light sheet focus. Use the marker at the top of the slider to set the global offset for the reference wavelength. You can fine-tune the chromatic correction for all other excitation wavelengths with the marker at the bottom. You can use the **Set as Reference** button to set the current excitation wavelength as the new reference wavelength. If the **Coupled Chromatic Correction** checkbox is selected, the chromatic correction is applied symmetrically to the left and right light sheet. Otherwise, the chromatic correction can be set individually for each light sheet.

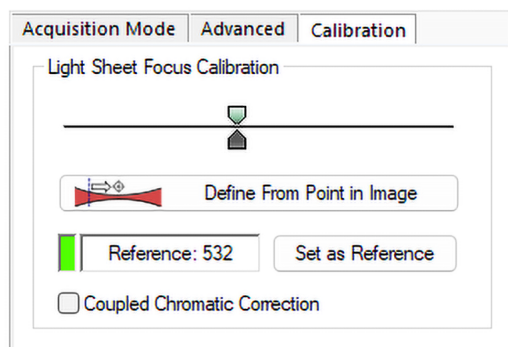


Figure 5.6: Calibration tab

### 5.2.8 Advanced tab

To change the **Light Sheet Processing** and the **Light Sheet Configuration**, you can use the **Advanced** tab.


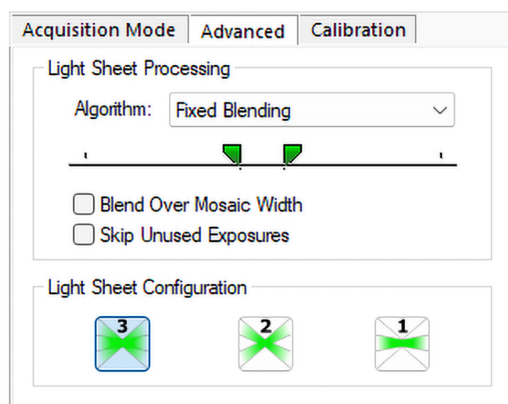


Figure 5.7: Advanced tab

To use the image algorithm **Fixed Blending**:

- 1 Select **Skip Unused Exposures**. By this, all exposures that are not going to contribute to the processed image are skipped.

- 2 Select **Blend Over Mosaic Width** to blend over the whole mosaic width. The following example shows a  $3 \times 3$  mosaic, in which the left-hand side is not covered by the right light sheet (orange) and the right-hand side is not covered by the left light sheet (blue). For the middle, both light sheets are cross-faded.

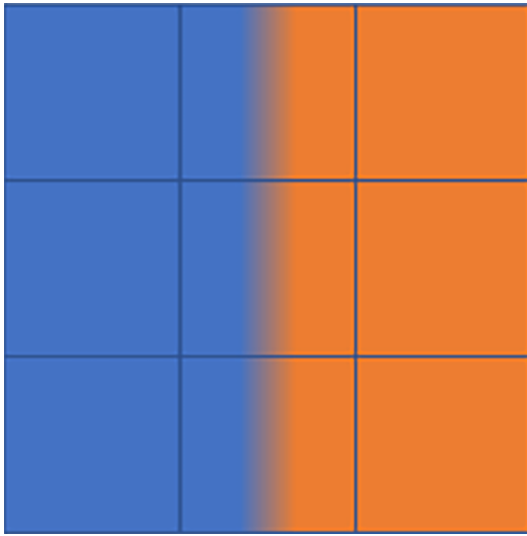


The **Blend Over Mosaic Width** combo box is not available if either the left or right light sheet is selected or if the x-axis is not one of the active devices.

If **Blend Over Mosaic Width** is selected, the acquisition speed can be increased and a lower value under **Laser Transmission Control** can be used to decrease photobleaching.

With **Dynamic Focus** activated, the number of steps for the dynamic focus is limited to parts of the ROI that will be visible in the image after blending the two sides.

To let the ImSpector software determine the **Light Sheet Processing** settings, which can be used to compensate errors from a poor sheet focus calibration:

- 1 Select **Adaptive Blending** under **Algorithm**.

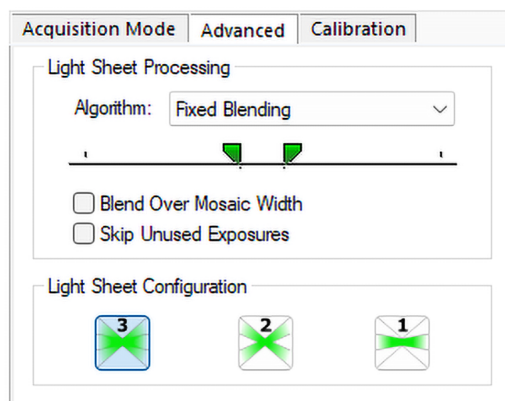

Under **Light Sheet Configuration**, you can select which laser beams are active.

### 5.2.9 Acquisition Mode tab

In the **Acquisition Mode** tab, you have two options:

- Dynamic Focus
- Fast Tiling Scan

Depending on your ordered configuration, **LightSpeed Mode** as a third option might be visible.

## Dynamic Focus

To activate the dynamic focus in the **Acquisition Mode** tab, you have to select the **Dynamic Focus** checkbox.


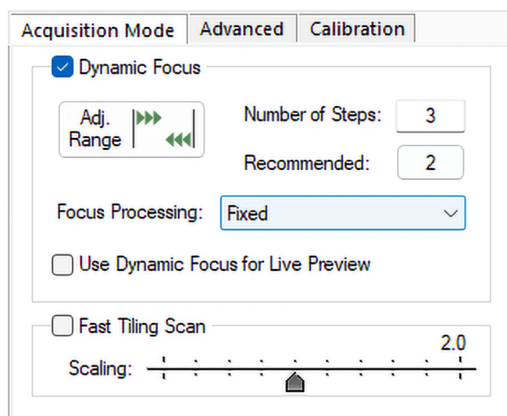


Figure 5.8: Dynamic Focus mode

With dynamic focus activated, the exposure time of the camera is slightly adjusted to the travel time of the focus to synchronize the camera acquisition with the movement of the light sheet focus.

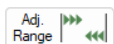

- 1 To toggle the visibility of the horizontal focus range indicators and adjust the focus range in the stack window, click the **Show dynamic focus range adjustment** button. During a measurement with enabled dynamic focus, only the area between the indicators will be covered by the light sheet focus.
- 2 Adjust the range of the dynamic focus to cover the region of interest.
- 3 Enter the **Number of Steps** or click **Recommended** to select the number of steps taken.

Consider that the **Number of Steps** corresponds to the number of recorded raw images. A large number of steps increases the image acquisition time and the amount of data significantly. The recommended number of steps approximates the number of images necessary to achieve a homogenous z-resolution.

To blend between the theoretical focus positions:

- 1 Select **Fixed** under **Focus Processing**. The results are best with proper sheet focus calibration in poor samples.

To estimate the sheet focus position during the measurement:

- 1 Select **Adaptive** under **Focus Processing**. This setting can compensate for bad sheet focus calibration.

To use the dynamic focus during live preview:

- 1 Select the **Use Dynamic Focus for Live Preview** checkbox.

With dynamic focus activated during live preview, the acquisition time is longer.

### 5.2.10 xyz-Table Z view

To move the sample in live preview in z-direction, you can use the **xyz-Table Z** view. The sample can be moved by the jog wheel, the mouse wheel, or the arrows on the keyboard. For the use of the mouse wheel or the arrows, the brain icon has to be clicked and the cursor kept on this position.

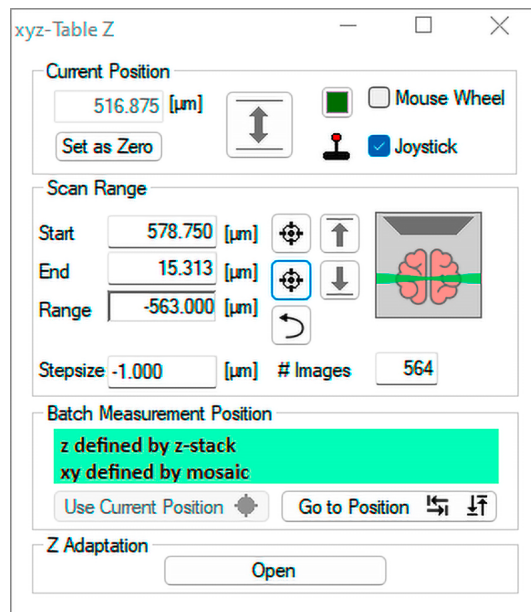


Figure 5.9: xyz-Table Z view

#### Current Position

After switching on the instrument, the sample arm moves to its start position. The start position in z-direction can be changed by entering the **Current Position**, or selecting **Mouse Wheel** or **Joystick**. By clicking **Set as Zero**, you can set a new reference point of the sample.



Figure 5.10: Current Position

#### Scan Range

Within this group box, you can configure the scan range together with the step size and the number of images. The start and end values can be typed into the fields, as well as the step size or the number of images. The range is automatically calculated.

Alternatively, you can move the stage to the start position and the end position by clicking the corresponding arrow and click the **Set current position as start** button and the **Set current position as end** button.

By entering **Stepsize** or **# Images**, you can define the number of images taken.

Scan Range

Start 1000.000 [μm]

End 2000.000 [μm]

Range 1000.000 [μm]

Stepsize 5.000 [μm] # Images 201

Figure 5.11: Scan Range

### Batch Measurement Position

Within this group box, you can configure a batch measurement.

### Z Adaptation

To adjust the laser power relative to the z-position:

- 1 Click **Z Adaptation > Open**.

Z Adaptation

Close

Select channel to adapt

None None None

Function

Points

Start End + - X

☐ Follow in Live Preview

- 2 Define a z-stack (see **Acquiring an image on page 35**) by entering a start and end position and **Stepsize** or **# Images**.
- 3 Select the channel that will be used.
- 4 Move to the start position.
- 5 Define the laser power at the start position and click **Start**.
- 6 Move to the end position.
- 7 Define the laser power at the end position and click **End**.
- 8 Select the **Function** of choice: **Points** or **Linear** are recommended.
- 9 Optional: To adapt the laser power during live preview, select the **Follow in Live Preview** checkbox. By this, the detection unit and/or sample can be moved in z, and the laser power is adapted in real time.

Start

End

#### 5.2.11 xyz-Table Visual XY view

To move the sample in live preview horizontally, you can use the **xyz-Table Visual XY** view. The sample can be moved by the jog wheel or the arrows on the keyboard. The field of view is represented by the green tile.

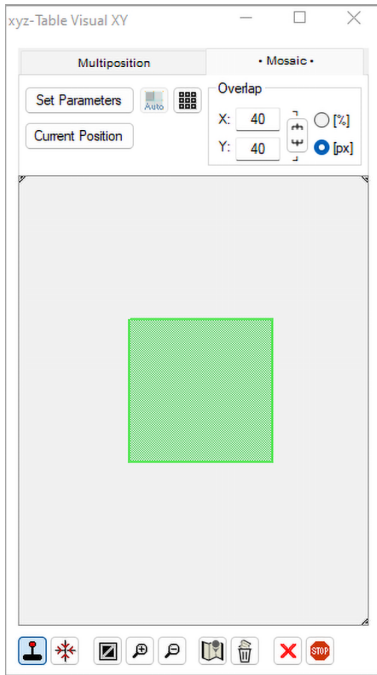


Figure 5.12: xyz-Table Visual XY Mosaic measurement

The window has got two tabs for **Mosaic** and **Multiposition** measurements. Both tabs use the buttons described in the following table.

|                                                                                     |                                                           |                                                                                                               |
|-------------------------------------------------------------------------------------|-----------------------------------------------------------|---------------------------------------------------------------------------------------------------------------|
|  | <b>Allow the joystick to move to the current position</b> | Activates or deactivates the jog wheel to move to the current position                                        |
|  | <b>Set current position as origin</b>                     | Sets the current position as origin                                                                           |
|  | <b>Zoom to the range of positions</b>                     | Zooms to the range of positions                                                                               |
|  | <b>Delete all positions of this mode</b>                  | Deletes all positions of this measurement mode                                                                |
|  | <b>Stop</b>                                               | Stops the xy-movement immediately and stops the measurement                                                   |
|  | <b>Zoom in</b>                                            | Zooms into the view                                                                                           |
|  | <b>Zoom out</b>                                           | Zooms out of the view                                                                                         |
|  | <b>Capture overview MIP at current position</b>           | Acquires overview maximum intensity projection (MIP) at current position (see <b>Simple view on page 51</b> ) |
|  | <b>Clear overview MIP</b>                                 | Clears overview MIP                                                                                           |
|  | <b>Define mosaic automatically</b>                        | Defines mosaic automatically when overview MIP is used (only available in the <b>Mosaic</b> tab)              |
|  | <b>Shift Mosaic positions to Multiposition mode</b>       | Converts <b>Mosaic</b> mode positions to <b>Multiposition</b> mode (only available in the <b>Mosaic</b> tab)  |




|                                                                                   |                                                                      |                                                                                                                                    |
|-----------------------------------------------------------------------------------|----------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------|
|  | <b>Set overlap in x and y to the same value (OFF state)</b>          | Activates or deactivates the same overlap value in x and y (in percent or pixel) (only available in the <b>Mosaic</b> tab)         |
|  | <b>Show ascending measurement order of positions</b>                 | Shows the ascending order of the measurement positions (only available in the <b>Multiposition</b> tab)                            |
|  | <b>Define ascending measurement order to shorten table movements</b> | Defines the descending order of measurements to reduce the movement of the sample (only available in the <b>Multiposition</b> tab) |
|  | <b>Set Z offsets for different positions</b>                         | Activates different z-offsets for different positions (only available in the <b>Multiposition</b> tab)                             |

**Table 5.2:** Buttons of the **xyz-Table Visual XY** view

### Mosaic measurements

Observe the instructions for creating a mosaic measurement under **Acquiring a 3D mosaic image on page 37** or **Acquiring a multicolor mosaic image on page 38**.


You can click the **Set Parameters** button to fine-tune the experiment by setting the **Size of One Image** and the **Number of Images**.

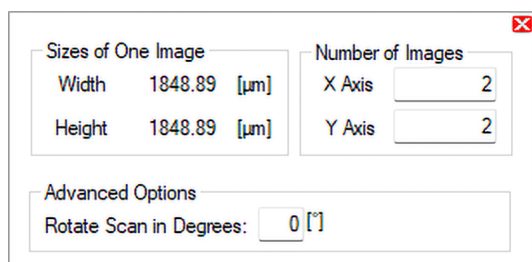

**Figure 5.13:** **Set Parameters** dialog

**Rotate Scan in Degrees** must be set to zero.

You can click the **Current Position** button to edit the current position in x- and y-direction by entering numerical values.

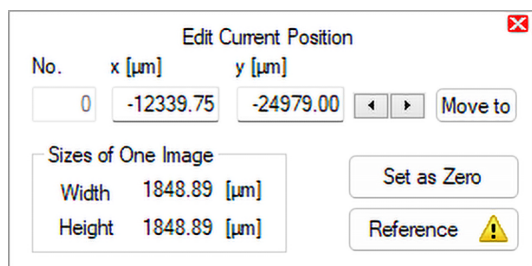

**Figure 5.14:** **Current Position** dialog

### Multiposition measurements

Observe the instructions for creating a multiposition measurement under **Acquiring a 3D multiposition dataset on page 40**.

To retrieve more details on the individual positions, click the **Edit Position** button. You can move through all active positions by using the arrows or the **Move to** button. Alternatively, you can enter numerical values for the current position. Select an active position to show its position and highlight it in dark gray color.

To set the coordinates of one position as reference point, you can click **Set as Zero**. Clicking **Reference** moves the stage to this reference point.

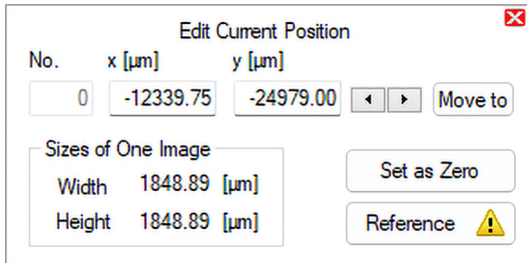


Figure 5.15: Current Position dialog

To show the ascending order of the measuring positions, you can click the **Show ascending measurement order of positions** button.

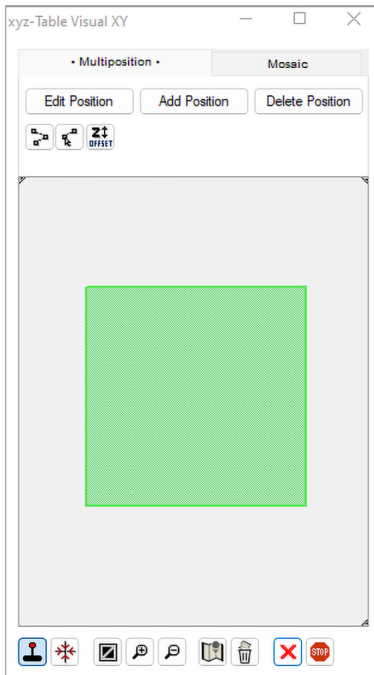



Figure 5.16: xyz-Table Visual XY Multiposition measurement

To set the ascending order of measuring positions, you can click the **Define ascending measurement order to shorten table movements** button.

To activate setting different z-offsets when clicking the **Edit Position** button, you can click the **Set Z offsets for different positions** button.

## 5.2.12 Camera

To control the camera settings of the instrument, you can use the **Camera** view. **Temperature Info** shows the temperature of the sensor, the camera, and the power supply.

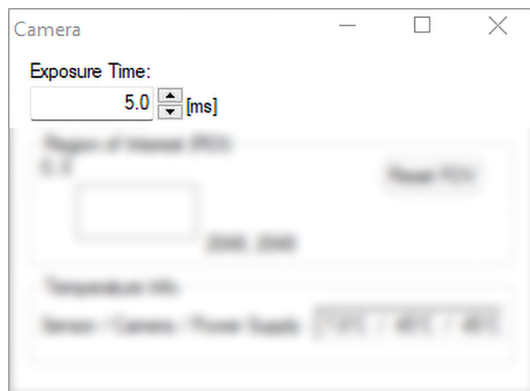

**Figure 5.17: Camera view**

- 1 Set an **Exposure Time** depending on the sample used. A suitable starting value can be anything between 100 and 200 ms.



# 6

## Setting up an experiment

During the installation of the instrument, the ImSpector software has been installed with the standard settings and validated for the use with the instrument. For more detailed information on the operation of the instrument, including light sheet alignment, refer to the instrument user manual.

### 6.1 Measurement Wizard

You can start and stop a measurement or start a live preview in the document window using the **Measurement Wizard**. If activated under **Settings > Measurement**, you can also set up a batch measurement.

|                                                                                     |                                |                                                           |
|-------------------------------------------------------------------------------------|--------------------------------|-----------------------------------------------------------|
|   | <b>Start Measurement</b>       | Starts the defined measurement                            |
|  | <b>Start Measurement</b>       | Starts the defined measurement in LightSpeed Mode         |
|  | <b>Stop Measurement</b>        | Stops the defined measurement                             |
|  | <b>Start Batch Measurement</b> | Starts the defined batch measurement                      |
|  | <b>Start Live Preview</b>      | Starts a live preview                                     |
|  | <b>Toggle Estimation</b>       | Toggles to estimate parameters for the active measurement |
|  | <b>Autosave</b>                | Activates autosave for the selected device                |

**Table 6.1:** Buttons of the **Measurement Wizard**

After setting up an experiment, the amount of necessary storage space is shown under **Status** as well as the estimated time. Under **Messages**, information and warnings for the measurement setup are shown.

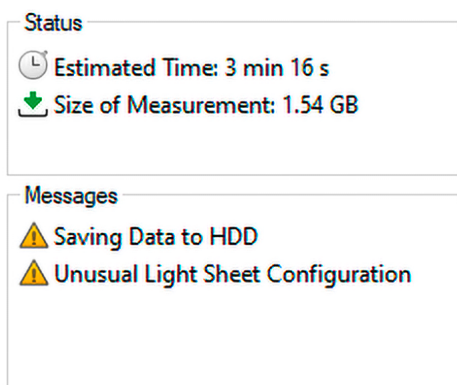


Figure 6.1: Status and Messages area

## 6.2 Instrument Mode

In the **Measurement Wizard**, you can select the **Instrument Mode** for the ordered configuration.

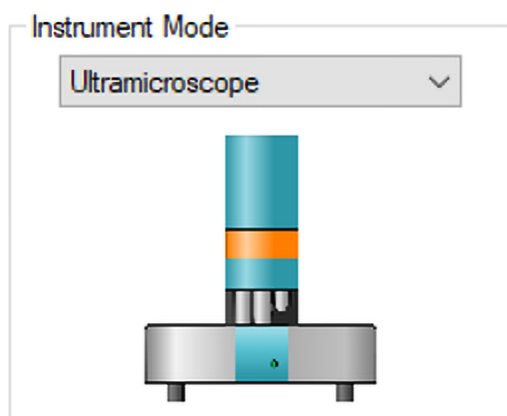


Figure 6.2: Instrument Mode list

All instrument modes are predefined by Miltenyi Biotec. If the settings need to be changed, contact Miltenyi Biotec Technical Support. See **Technical support on page 67**.

## 6.3 Measurement Mode

To acquire data, you must define the order of the devices to be used during acquisition. You can freely define these devices by selecting them in the **Devices** lists or use predefined settings, called measurement modes. The available measurement modes can be selected in the **Measurement Mode** list.

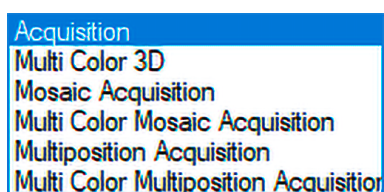


Figure 6.3: Measurement Mode list

All measurement modes are predefined by Miltenyi Biotec. If the settings need to be changed, contact Miltenyi Biotec Technical Support. See **Technical support on page 67**.

Deactivating the **Autosave** function for the **xyz-Table Z** is recommended.

### 6.3.1 Acquiring an image



1 To start the live preview, click the **Start Live Preview** button in the **Measurement Wizard**.

2 Move in the z-direction to the start position.

3 Click **Current Position > Set as Zero** in the **Settings 2** pane.



4 Click the **Set current position as start** button under **Scan Range > Start**.

5 Move in the z-direction to the end position.



6 Click the **Set current position as end** button under **Scan Range > End**.

7 Define the step size, corresponding to the number of images, by entering it under **Scan Range > Stepsize**.



8 To stop the live preview, click the **Stop Live Preview** button. The laser is switched off.

9 Check **Settings** in the **Autosave** area of the **Experiment Manager**.




10 To start the acquisition, click the **Start Measurement** button.



### 6.3.2 Acquiring a multicolor z-stack


1 In the **Measurement Wizard**, select **Measurement Mode > Multi Color 3D**.

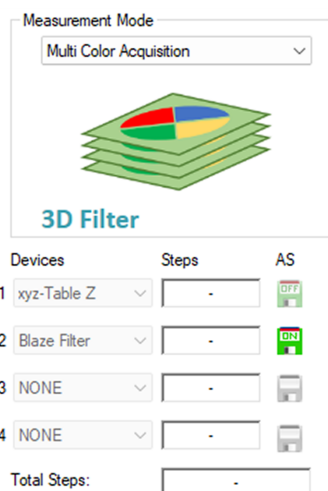

2 To activate the **Autosave** function for the devices, click the **Autosave** button at the end of the row under **AS**. The button turns green.

Optional: If you use **Fast Tiling Scan**, check if the **Autosave** function for the **xyz-Table Z** is deactivated. The button must appear dimmed.

3 Select the correct excitation and emission wavelengths for your sample under **Channel Setup** in the **Settings 1** pane. The selected checkboxes are used.

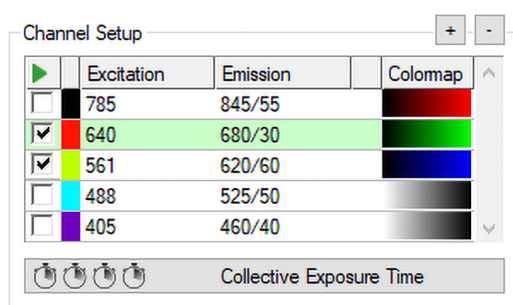

4 Move the slider under **Laser Transmission Control** while the filter is highlighted to set the laser power for each individual filter.

**NOTICE** Higher laser power bleaches the sample faster.


5 Select the shortest emission wavelength, and create a sharp picture.



6 Correct the different focus lengths for different emission wavelengths using the **Chromatic Correction**, starting with the next longer wavelength.



7 If all emission wavelengths have been corrected, click **Save Settings** on the toolbar.

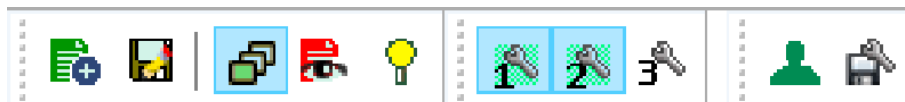




- 8 To start the acquisition, proceed as described under **Acquiring an image on page 35**.

### 6.3.3 Acquiring a 3D mosaic image

- 1 In the **Measurement Wizard**, select **Measurement Mode > 3D Mosaic**.




- 2 In the **xyz-Table Visual XY** view, click the **Mosaic** tab.

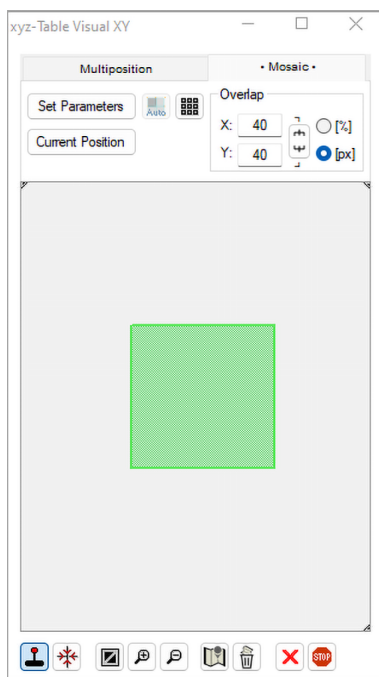

- 3 Delete all positions of this **Measurement Mode**.
- 4 To set an **Overlap**, enter a value and select percentage (%) or pixel (px).

The overlap must be at least 40 px.

- To draw a ROI you would like to cover with a mosaic, click-drag with the right mouse button. From the pop-up menu, select **Define mosaic here**.

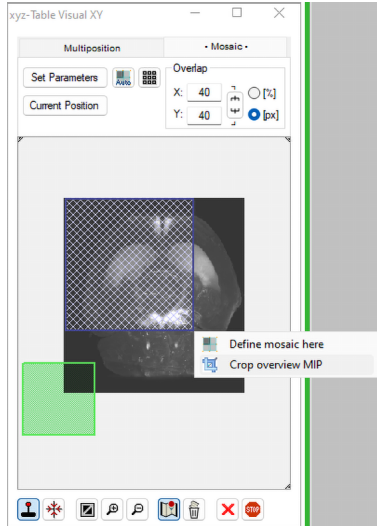


- To activate the **Autosave** function for the devices, click the **Autosave** button at the end of the row under **AS**. The button turns green.  
Optional: If you use **Fast Tiling Scan**, check if the **Autosave** function for the **xyz-Table Z** is deactivated. The button must appear dimmed.
- To start the acquisition, proceed as described under **Acquiring an image on page 35**.

## 6


### 6.3.4 Acquiring a multicolor mosaic image

- In the **Measurement Wizard**, select **Measurement Mode > Multi Color Mosaic Acquisition**.

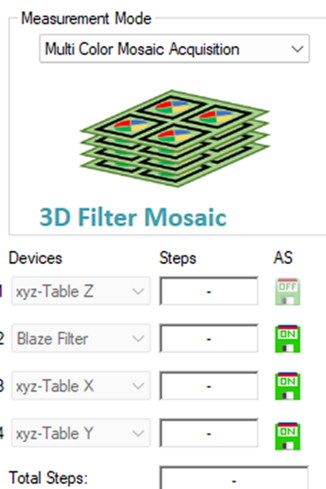

- To activate the **Autosave** function for the devices, click the **Autosave** button at the end of the row under **AS**. The button turns green.  
Optional: If you use **Fast Tiling Scan**, check if the **Autosave** function for the **xyz-Table Z** is deactivated. The button must appear dimmed.

- 3 Select the correct excitation and emission wavelengths for your sample under **Channel Setup** in the **Settings 1** pane. The selected checkboxes are used.

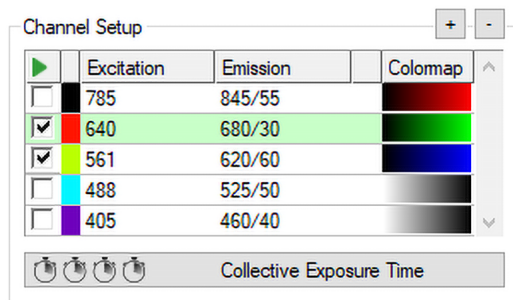

- 4 Move the slider under **Laser Transmission Control** while the filter is highlighted to set the laser power for each individual filter.

**NOTICE** Higher laser power bleaches the sample faster.




- 5 Correct the different focus lengths for different emission wavelengths using the **Chromatic Correction**, starting with the next longer wavelength.



- 6 If all emission wavelengths have been corrected, click **Save Settings** on the toolbar.




- 7 In the **xyz-Table Visual XY** view, click the **Mosaic** tab.

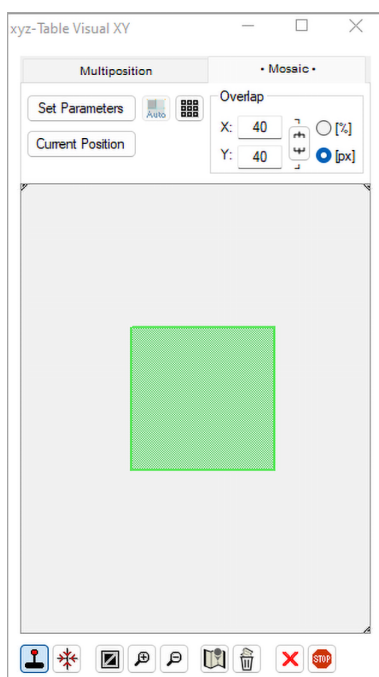

- 8 To set an **Overlap**, enter a value and select percentage (%) or pixel (px).

The overlap must be at least 40 px.

- 9 To draw a ROI you would like to cover with a mosaic, click-drag with the right mouse button. From the pop-up menu, select **Define mosaic here**.

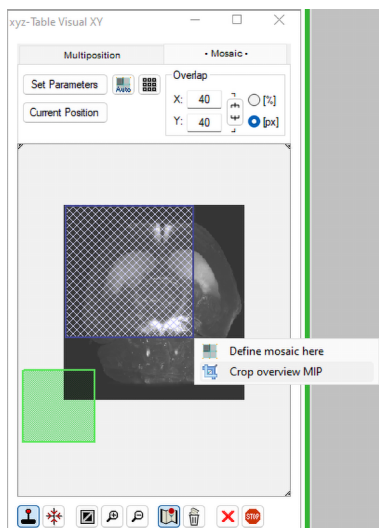


- 10 To start the acquisition, proceed as described under **Acquiring an image on page 35**.

### 6.3.5 Acquiring a 3D multiposition dataset


- 1 In the **Measurement Wizard**, select **Measurement Mode > Multiposition Acquisition**.



- 2 In the **xyz-Table Visual XY** view, click the **Multiposition** tab.

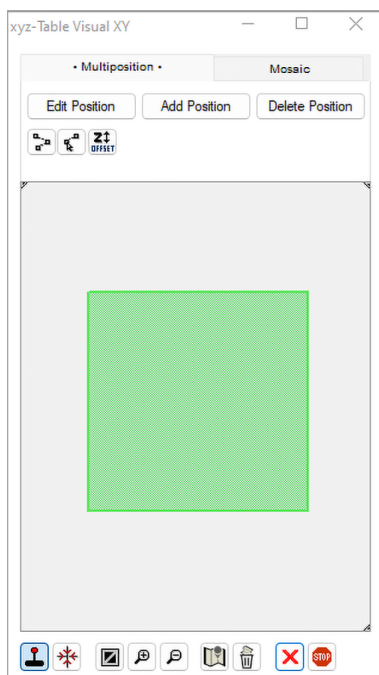

- 3 Delete all positions of this measurement mode.
- 4 Move the measurement positions to your region of interest.
- 5 Click **Add Position** for each position you want to add.

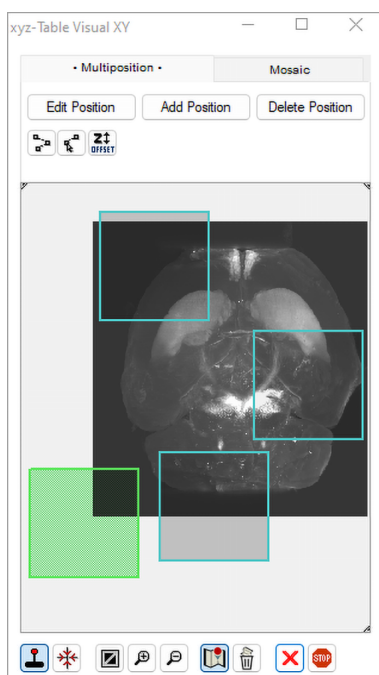

- 6 To activate the **Autosave** function for the devices, click the **Autosave** button at the end of the row under **AS**. The button turns green.  
Optional: If you use **Fast Tiling Scan**, check if the **Autosave** function for the **xyz-Table Z** is deactivated. The button must appear dimmed.
- 7 To start the acquisition, proceed as described under **Acquiring an image on page 35**.

## 6.4 Experiment Manager

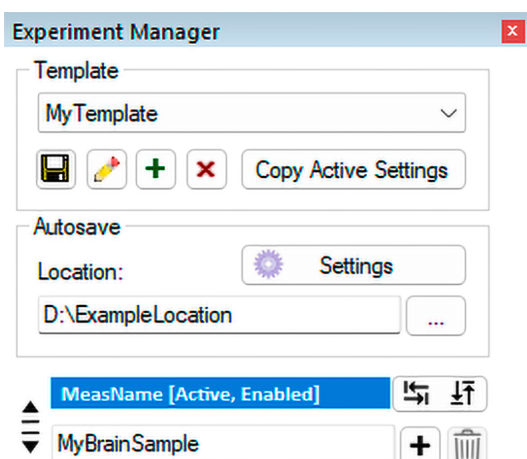


Figure 6.4: Experiment Manager

With the **Experiment Manager**, you can save, load, or configure predefined templates and set up experiments. Default templates are provided by Miltenyi Biotec. These are specialized templates for the use in combination with the multisample holders. If a template is changed and not saved, ImSpector asks if the settings shall be updated for the current template before closing. In batch measurements, which consist of many single measurements, different status are possible. **Active** means that this measurement has focus. **Enabled** means that this measurement will be performed if a batch measurement is started. The status are indicated by a color code. See **Table 6.2**.

You can copy settings from the active measurement by clicking **Copy Active Settings**. The following window is shown:

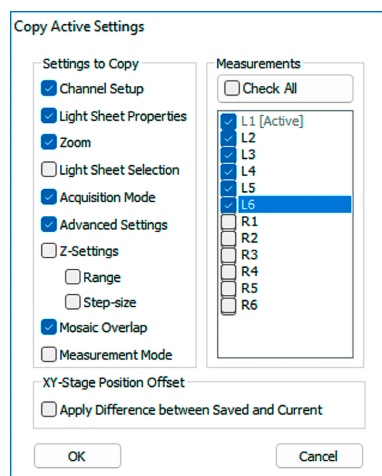


Figure 6.5: Copy Active Settings window

You can select the settings you would like to copy to other measurements on the left and select the measurements which those settings should be applied to on the right. To select all available measurements, use the **Check All** checkbox.

The **Settings** button in the **Autosave** group box provides more advanced options to modify your data output format.

| Color      | Status                             |
|------------|------------------------------------|
| Dark blue  | <b>Active, Enabled</b>             |
| Light blue | <b>Inactive, Enabled</b>           |
| Yellow     | <b>Active, Disabled</b>            |
| Grey       | <b>Inactive, Disabled</b>          |
| Red        | <b>Active, Enabled</b> current xyz |

**Table 6.2:** Color code of measurement status

## 6.5 Setting up a batch measurement

With the batch measurement feature, you can perform several measurements one after the other without any further intervention. These measurements can be saved as a template.



- 1 In the **Experiment Manager**, click the **Create New Measurement** button.
- 2 Name your measurement by right-clicking the colored top bar. An autosave prefix has also to be set, for example, MyBrainSample.
- 3 Define the measurement settings for your measurement. See **Setting up an experiment on page 33**.
- 4 Click the **Create New Measurement** button repeatedly to add further measurements, if required. By dragging the indicator on the left-hand side, the order is changed.
- 5 Click the **Use Current Position** button under **Batch Measurement Position**, or define a mosaic.
- 6 Click the **Go to Position** button to move the stage to the measurement position.
- 7 Optional: By right-clicking, the selected measurement can be made inactive.
- 8 To start the acquisition, click the **Start Batch Measurement** button.



For safety reasons, it is not possible to change the objective lenses during a batch measurement.





# 7

## Data visualization and analysis

The acquired data is shown in the stack window of the ImSpector software. Once the data is saved, you can open it via the **Image Series Viewer**. The data can be broken down into its constituent parts, examined, and analyzed.

### 7.1 Stack window

The stack window in the ImSpector software shows the acquired data:

- The file name is shown in the header.
- The toolbar located on the right of the image contains buttons for frequently used data handling functions. See **Table 7.1**.
- Some options can be accessed by right-clicking the image. See **Table 7.2**.

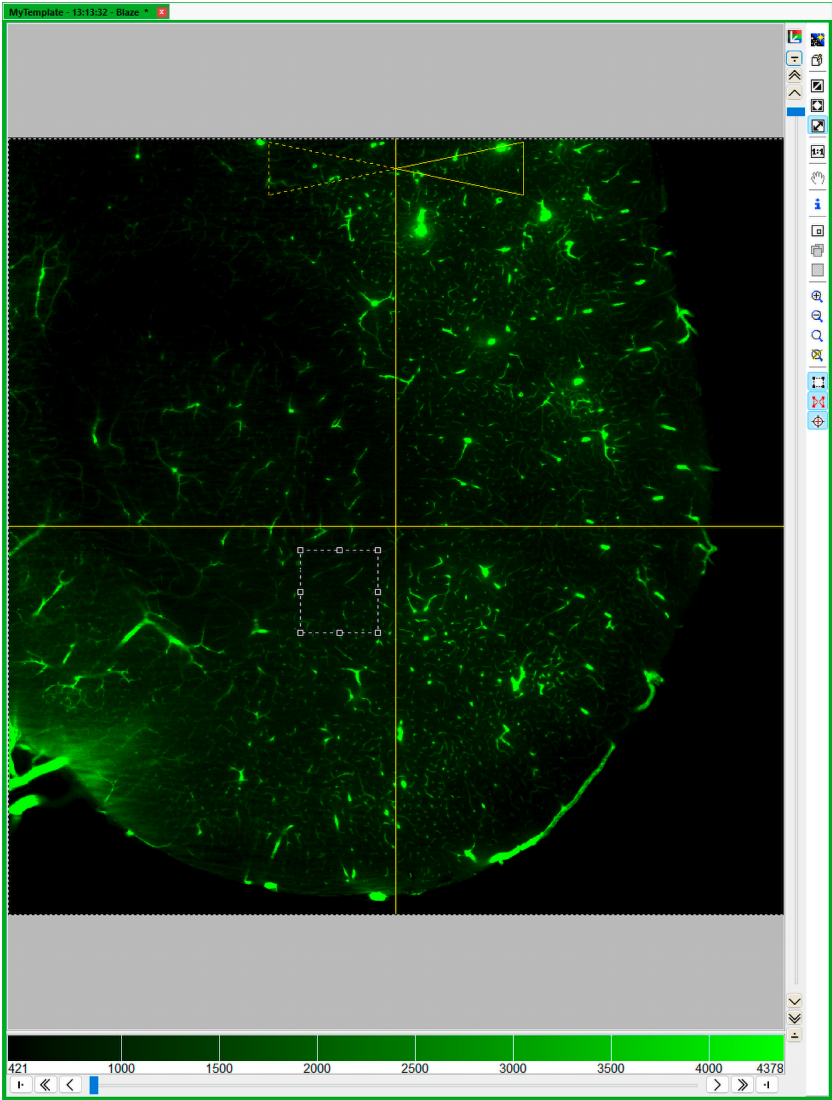



Figure 7.1: Stack window

The row below the image shows the colormap, which is used to display the different pixel values in different colors. You can get additional information for each pixel by moving the cursor over it.

|                                                                                     |                                    |                                                                                                                                                      |
|-------------------------------------------------------------------------------------|------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------|
|  | <b>Display Settings</b>            | Opens the <b>Display Settings</b> dialog below the stack window                                                                                      |
|  | <b>Auto Contrast</b>               | Adjusts the contrast of the colormap automatically to the minimum and maximum pixel value within the whole image or a shown region of interest (ROI) |
|  | <b>Change Size and Dim</b>         | Opens the <b>Stack properties</b> window to enable a change in size and dimension of the data stack                                                  |
|  | <b>Fit Frame to Image</b>          | Fits the frame to the image. The aspect ratio is kept.                                                                                               |
|  | <b>No Aspect Ratio</b>             | Fits the image to the frame. The aspect ratio is NOT kept.                                                                                           |
|  | <b>Fit Image to Frame</b>          | Maximizes the image in a frame of constant size                                                                                                      |
|  | <b>Show Image in Original Size</b> | Shows the image in its original size                                                                                                                 |

|                                                                                     |                                     |                                                                                                                                                              |
|-------------------------------------------------------------------------------------|-------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------|
|    | <b>Move the Image</b>               | Grabs the image and moves it inside the frame; can only be used if the image is bigger than the frame                                                        |
|    | <b>Show Measurement Information</b> | Opens an additional window with metadata of the image                                                                                                        |
|    | <b>Set as New ROI</b>               | After drawing a ROI, this can be set as new camera or scan ROI. If the ROI is not shown, the button changes to <b>Reset ROI to full FOV</b> (field of view). |
|    | <b>Show Merged Channel</b>          | Shows an additional merged channel                                                                                                                           |
|    | <b>Show Merged Channel Only</b>     | Shows only the merged image                                                                                                                                  |
|    | <b>Zoom In</b>                      | Makes the image appear nearer                                                                                                                                |
|    | <b>Zoom Out</b>                     | Makes the image appear farther away                                                                                                                          |
|    | <b>Zoom to Selection</b>            | Zooms into the ROI                                                                                                                                           |
|    | <b>Reset Zoom</b>                   | Resets the zoom settings                                                                                                                                     |
|    | <b>ROI for Active Profile</b>       | Shows or hides the ROI                                                                                                                                       |
|  | <b>Show Crosshair</b>               | Shows or hides the crosshair                                                                                                                                 |

**Table 7.1:** Buttons of the toolbar in the stack window

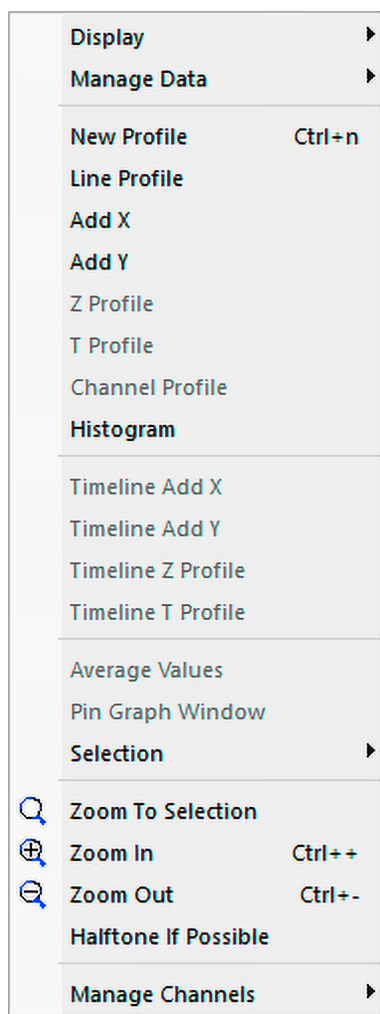

**Figure 7.2:** Right-click image options

| Name                 | Function                                                                                                                        |
|----------------------|---------------------------------------------------------------------------------------------------------------------------------|
| Display              | Changes display options                                                                                                         |
| Manage Data          | Copies, pastes, or edits data stacks                                                                                            |
| New Profile          | Adds a new profile                                                                                                              |
| Line Profile         | Shows an intensity profile along a drawn line                                                                                   |
| Add X                | Sums up all line functions in an interval x if a ROI was set                                                                    |
| Add Y                | Sums up all line functions in an interval y if a ROI was set                                                                    |
| Z Profile            | Shows an intensity profile around the z-axis (third data axis) if a ROI was set                                                 |
| T Profile            | Shows an intensity profile around the t-axis (fourth data axis) if a ROI was set                                                |
| Channel Profile      | Selects the channel profile                                                                                                     |
| Histogram            | Shows a histogram if a ROI was set; the x-axis represents the pixel value and the y-axis how often the pixel values are present |
| Average Values       | Shows the average values using <b>Add X</b> or <b>Add Y</b>                                                                     |
| Group Plot Window    | Pins a graph window to the stack window                                                                                         |
| Selection            | Allows to change various ROI values                                                                                             |
| Halftone If Possible | Interpolates between pixels                                                                                                     |
| Manage Channels      | Allows to adjust various settings of the different color channels                                                               |

**Table 7.2:** Right-click image options

## 7.2 Image Series Viewer

The Image Series Viewer allows you to quickly review data that was saved directly to the hard drive via the **Autosave** function. The Image Series Viewer consists of three separate windows:

- **Image Series Browser**
- **Image Series Axes**
- **Image Series Manager**

### 7.2.1 Image Series Browser

The **Image Series Browser** provides an overview of all recent measurements of the current user.

|                                                                                     |                                            |
|-------------------------------------------------------------------------------------|--------------------------------------------|
|  | Measurement was completed                  |
|  | Measurement was canceled                   |
|  | Measurement data was successfully imported |

**Table 7.3:** Status of measurement

## 7.2.2 Image Series Manager


With the **Image Series Manager**, you can import and export data, run a processing workflow, or reuse measurement settings.

## 7.2.3 Image Series Axes

In the **Image Series Axes** window, you can view a slice of a multidimensional data stack. In addition, you can define operations and destinations for individual axes. The selected operations will be applied when you import or export the data. Some operations, for example **Maximum intensity projection**, can be restricted to a range on the axis. The range is specified by additional indicators on the slider.

## 7.3 Stack Properties

In the **Stack Properties** dialog, you can see information about the stack itself.

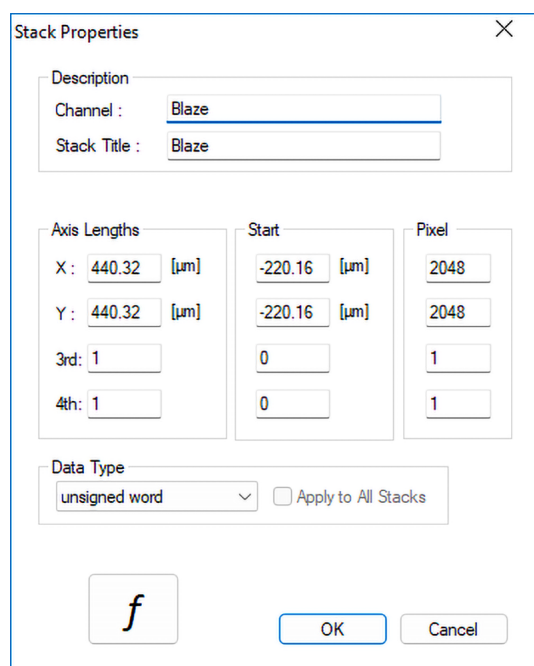

**Figure 7.3:** Stack Properties dialog


- The **Channel** and the **Stack Title** are shown in the **Description** box.
- The dimensions of the stack in  $\mu\text{m}$  are shown in the **Axis Lengths** box.
- The starting values in  $\mu\text{m}$  are shown in the **Start** box.
- The number of pixels and slices in all dimensions are shown in the **Pixel** box.
- The data type can be selected in the **Data Type** menu.
- The **f** button can be used to open the last used analysis tool again.

## 7.4 Changing the view

The ImSpector software can handle up to four-dimensional stacks. You can adjust the way acquired data is presented via the **Display Settings** dialog:



- 1 Click the **Display Settings** button to open the **Display Settings** dialog below the stack window.

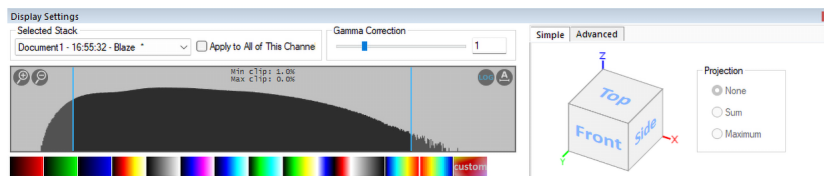

The **Display Settings** dialog consists of two parts:

- Colormap dialog
- Axis handling dialog with **Simple** and **Advanced** settings

### 7.4.1 Colormap dialog

The ImSpector software uses colormaps to distinguish between different pixel values. The current range is shown below the image in the stack window. Click a channel or a stack to change the corresponding settings. Alternatively, you can select the channel or the stack directly under **Selected Stack**. If you select the **Apply to All of This Channel** checkbox, you change the colormap for all stacks instead of only the selected one.


- Click the **Adjust Colormap Range** button to automatically adjust the colormap range.
- Move the blue bars with the mouse to adjust the colormap range manually.
- Click the **Zoom In** button to zoom into the selected histogram.
- Click the **Zoom Reset** button to reset the zoom.
- Press the **LOG** button to switch between linear and natural logarithmic scaling.
- Move the slider under **Gamma Correction** to adjust the gamma correction.

|                                                                                     |                              |                                                             |
|-------------------------------------------------------------------------------------|------------------------------|-------------------------------------------------------------|
|  | <b>Zoom In</b>               | Zooms into the selected histogram range                     |
|  | <b>Zoom Reset</b>            | Resets the histogram zoom                                   |
|  | <b>LOG</b>                   | Switches the scaling between linear and natural logarithmic |
|  | <b>Adjust Colormap Range</b> | Adjusts the colormap range                                  |

**Table 7.4:** Buttons of the **Display Settings** dialog

### 7.4.2 Simple view


The **Simple** tab shows a cube with the sides **Top**, **Front**, and **Side**. Additionally, the axes x, y, and z are shown.

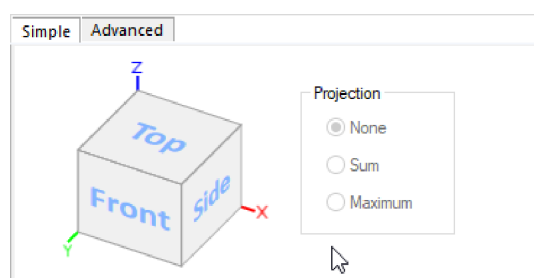

**Figure 7.4:** Simple view

You can click one of the cube sides to change the view.


To view the data from the top:

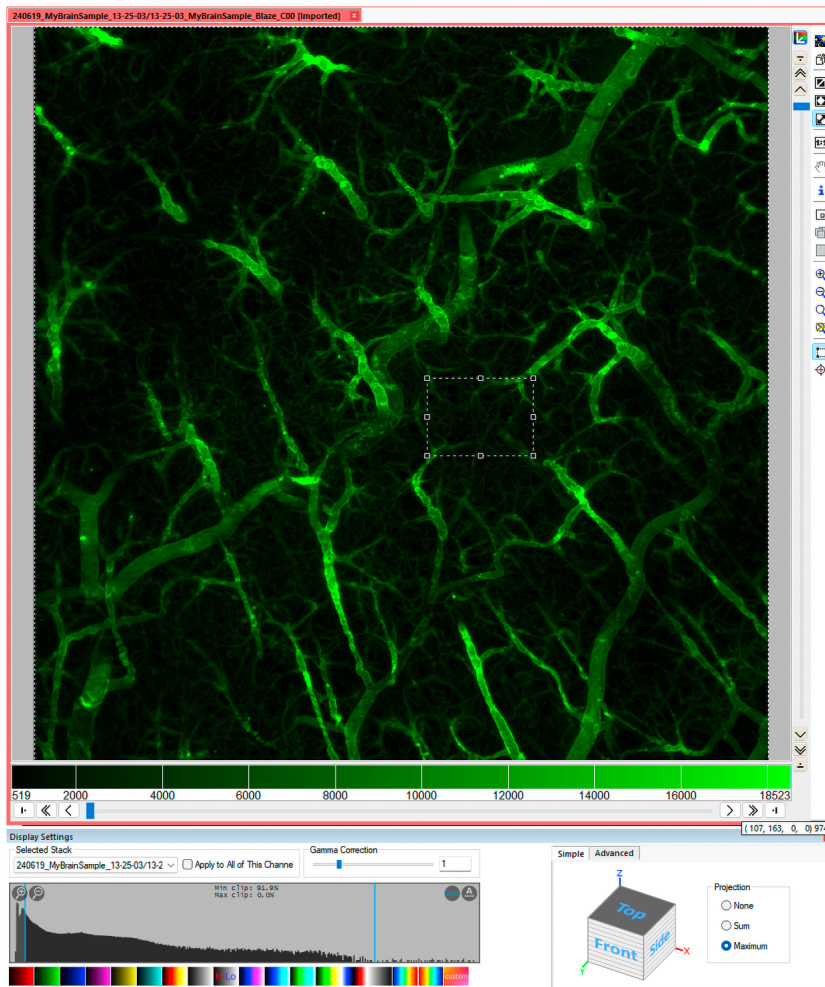

- 1 Click **Top**. Use the scroll bar on the right in the stack window to scroll along the z-axis from top to bottom through the data.




To view the data from the front:

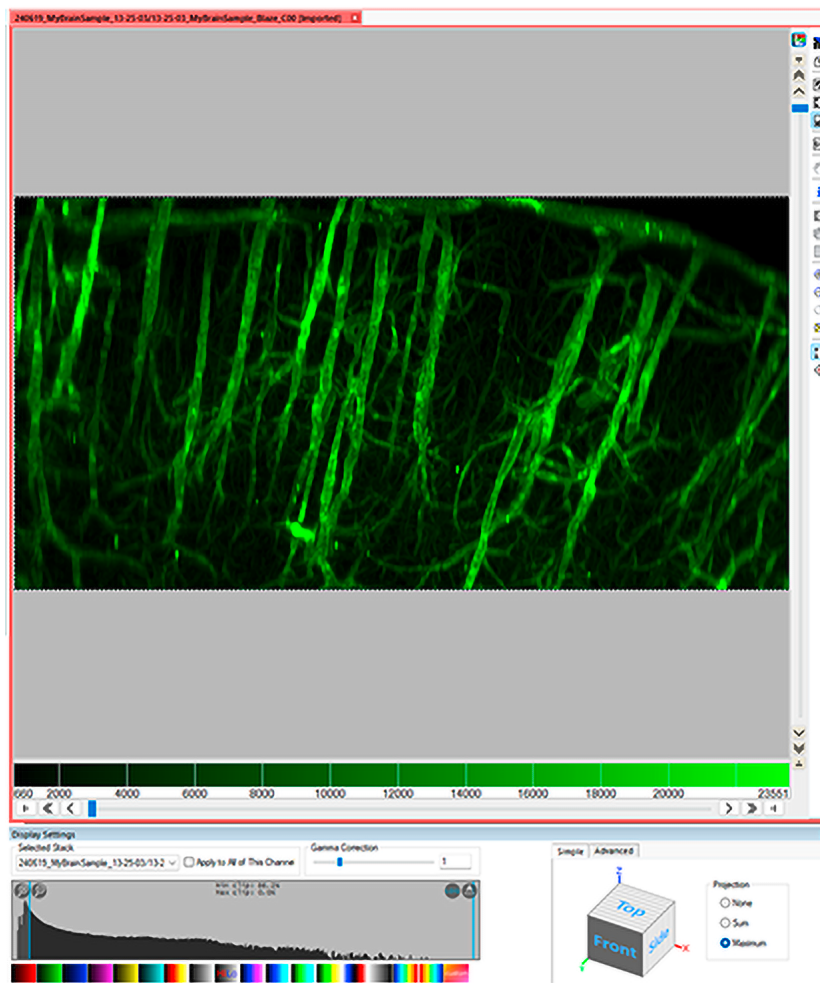

- 1 Click **Front**. Use the scroll bar on the right in the stack window to scroll along the y-axis from front to back through the data.

To view the data from the side:

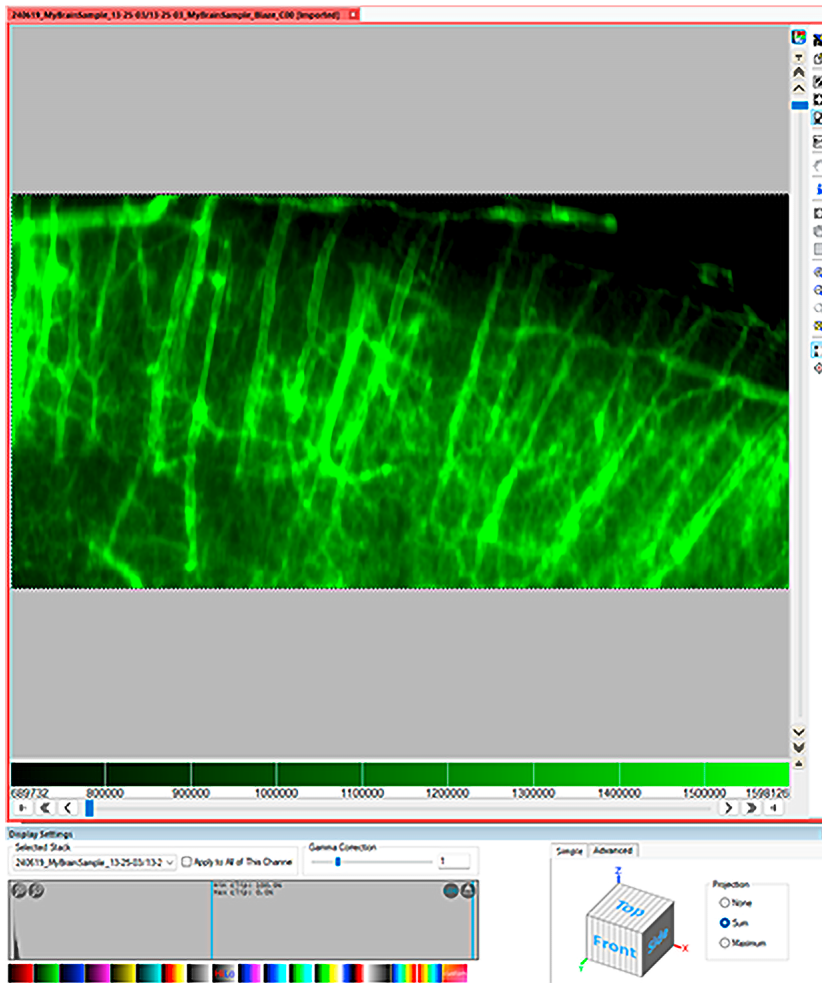

- 1 Click **Side**. Use the scroll bar on the right in the document window to scroll along the x-axis through the data.

A projection can be calculated for each view. By default, **Projection > None** is set, so not any projection is calculated:

- Click **Projection > Sum** to calculate a projection that sums up all pixels along an axis.
- or
- Click **Projection > Maximum** to calculate a projection that uses the maximum value for a pixel along the axis.

### 7.4.3 Advanced view

The **Advanced** tab shows a summary of all used projections. You can change the projections by clicking them in the list.

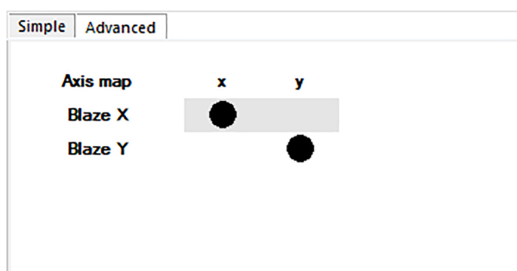


Figure 7.5: Advanced view

## 7.5 Merging and splitting

In the default view, the ImSpector software shows data stacks in **Gallery Mode**. In this mode, data stacks are shown side by side.

To create an overlay of multiple data stacks:

- 1 Right-click the image and select **Manage Channels > Merge Stacks**.

To split a data stack along an axis into a single layer stack:

- 1 Right-click the image and select **Manage Channels > Split Projected Axis**.

## 7.6 Defining a region of interest (ROI)

The ImSpector software allows you to select parts of the data in the stack window. This ROI can be used for various tasks such as copying data, creating profiles, or adjusting colormaps. The pixel positions and the physical size of the ROI are shown.

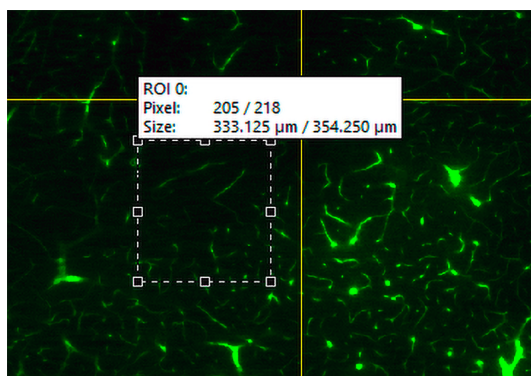



Figure 7.6: Region of interest

To define a ROI:



- 1 Ensure that the **ROI for Active Profile** button is activated.

- 2 Right-click the handles on the ROI to open a menu that allows you to customize the shape between rectangle, circle, freeform, and polygon.

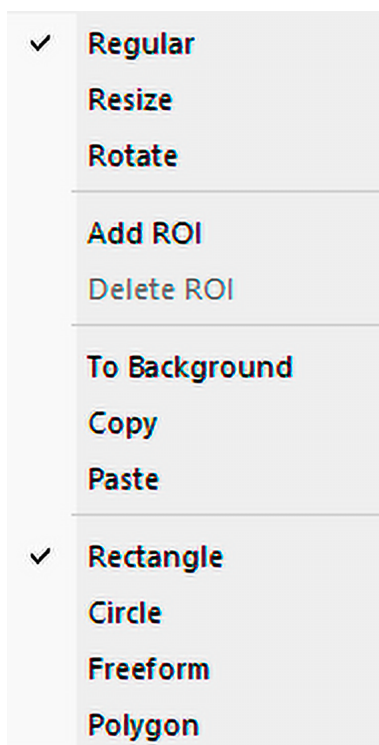

- 3 Click the **Set as new ROI** button to apply the settings.

## 7.7 Showing measurement information

The ImSpector software can show various metadata of your measurement.

The following metadata can be processed:

- Instrument information (laser or camera settings, objectives used, scanner, ...)
- Measurement information (creation date, instrument mode, ...)
- Stack information (pixel number, physical size, ...)
- User-specific data (name, institution, ...)
- Comments of the measurement

To show measurement information:



- 1 Click the **Show Measurement Information** button on the right in the stack window.

To edit comments:

- 1 Click **Edit comment**. A free text field opens.
- 2 Enter any text.

# 8

## Data management

This chapter describes how to open, save, and export data.

### 8.1 Types of windows

Depending on the type of data, the ImSpector software uses different types of windows:

- Stack windows are used for measurements, viewing and analyzing data. A stack window can have several document windows in a tab structure as follows: a **green** tab with access to the instrument or a **gray** tab without. Also see **Stack window on page 45**.
- Profile windows are used to show graphs from evaluation functions.
- Stack windows opened by the **Image Series Viewer** are similar to generic stack windows. Yet all data axis handling actions can be performed in the **Display Settings** dialog. This window has a **red** tab.

### 8.2 Managing files

The stack window is always located on the right-hand side of the ImSpector software.

If you have more than one document open, you can either have them in separate windows or in a tab structure. To switch between these views, activate or deactivate **Window > Show Document Windows in Tabs** in the menu bar.

If you right-click the tab of the stack window, you can select from the following options:

- Open a new stack window or a new graph window
- Close a document or a window
- Save a document
- Export a document

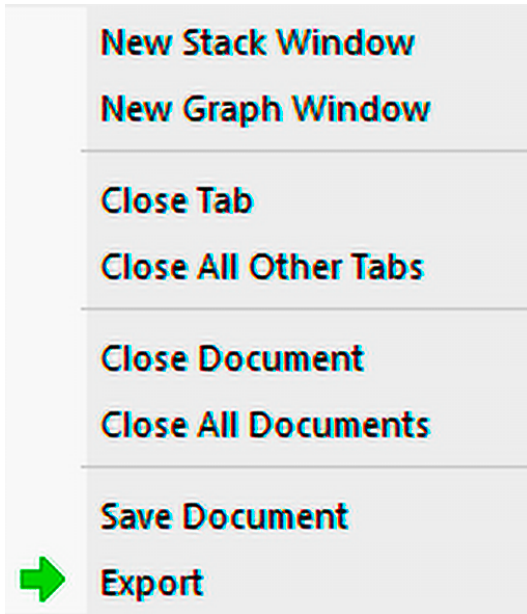


Figure 8.1: Stack window options

If you right-click a document that has no hardware access (gray tab), the option **Activate Hardware** is also available.

The **Experiment Manager** is located on the far right of the ImSpector software. If not, deactivate hiding it under **Window > Hide Experiment Manager**. For the function of the individual buttons, see **Table 8.1**.

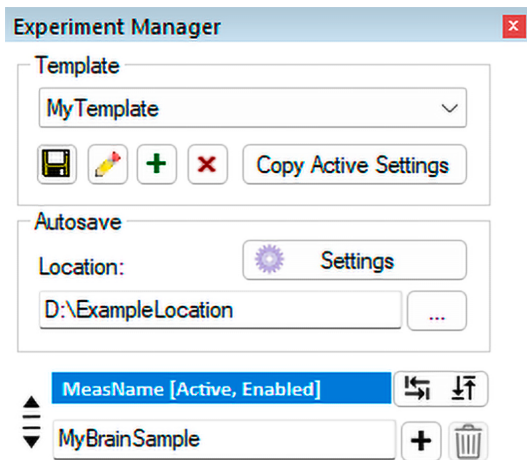



Figure 8.2: Experiment Manager

|                                                                                   |                                                                |                                                                                         |
|-----------------------------------------------------------------------------------|----------------------------------------------------------------|-----------------------------------------------------------------------------------------|
|  | <b>Save current settings</b>                                   | Saves an open template                                                                  |
|  | <b>Rename selected custom template</b>                         | Renames the selected custom template. <b>Note:</b> Default templates cannot be renamed. |
|  | <b>Save current settings as new template</b>                   | Creates a new template                                                                  |
|  | <b>Delete selected custom template</b>                         | Deletes the selected template                                                           |
|  | <b>Create new measurement</b>                                  | Adds a further measurement to the experiment                                            |
|  | <b>Remove measurement</b>                                      | Removes the selected measurement from the experiment                                    |
|  | <b>Activate and move stage to defined measurement position</b> | Activates the selected measurement                                                      |

**Table 8.1:** Buttons of the **Experiment Manager**

### 8.2.1 New documents


New documents can be created via **File > New** in the menu bar. See **File menu on page 14**. The new documents can be created with the current settings and the default settings.

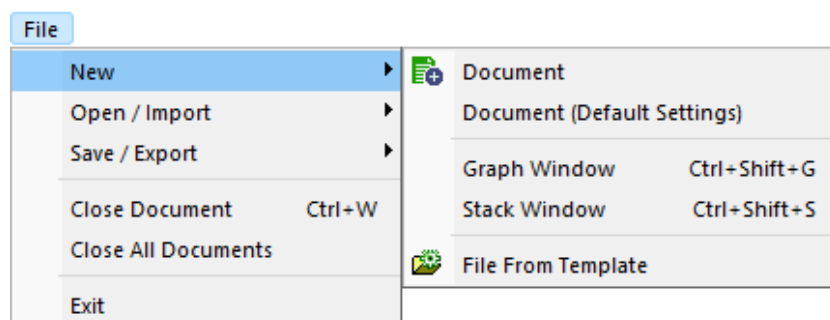

**Figure 8.3:** File menu

### 8.2.2 New graph or stack windows

New graph windows or new stack windows can be created via **File > New** in the menu bar. See **File menu on page 14**. The new windows show TIF files.

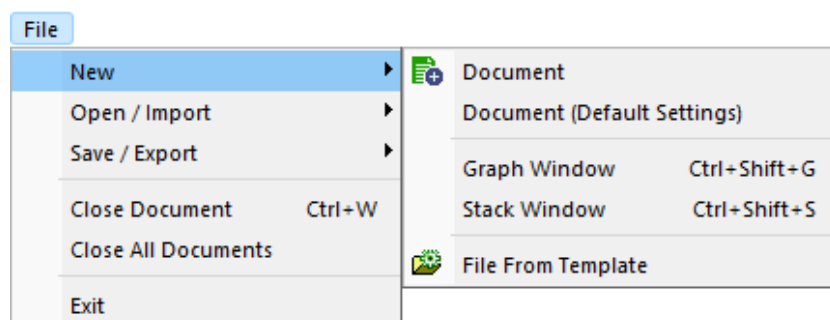

**Figure 8.4:** File menu

### 8.2.3 Open documents

Existing documents can be opened or imported via **File > Open / Import** in the menu bar. See **File menu on page 14**. Files with the extension .msr can be opened. Data saved via the **Autosave** function can be imported with the **Image Series Viewer**.

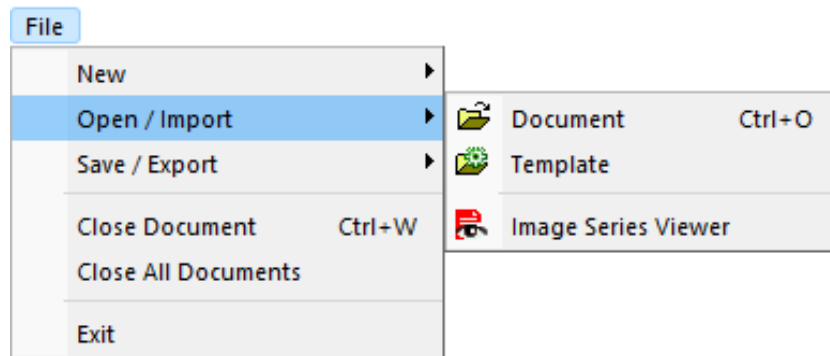


Figure 8.5: File menu

### 8.2.4 Save documents

Documents can be saved or exported via **File > Save / Export** in the menu bar. See **File menu on page 14**.

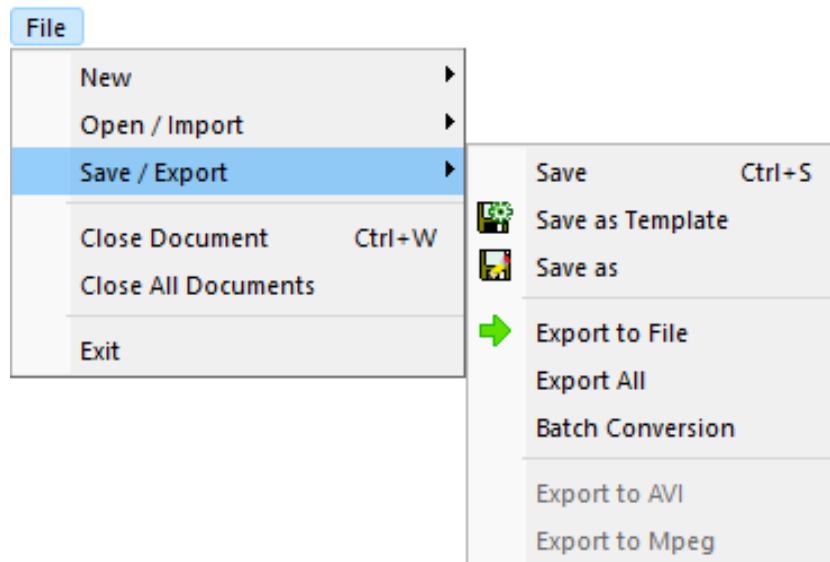


Figure 8.6: File menu

The ImSpector software saves data as SPV files and measurements as MSR files.

MSR files can be exported to different file formats, for example:

- TIFF, OME-TIFF
- RGB TIFF
- Binary Double Files (image data without any other file information)
- ASCII Data Files
- Postscript Files

To export an open file as a TIFF file:

- 1 Select the document to export.
- 2 Select **File > Save / Export > Export to File** or **File > Save / Export > Export All**.



- 3** Select the color depth **8 Bit** or **16 Bit** under **Data Type**.
- 4** Optional: Select the **Scale data** checkbox to avoid data truncation errors.
- 5** Optional: Clear the **Use .ome.tiff Extension** checkbox. If OME-TIFF is not used, the third and/or fourth dimension data axes can be split into single TIFF files.
- 6** Select either **Multi Channel** or **Single Channel** to save the TIFF files for all channels or for the selected channel only.



# 9

## Settings

### 9.1 User settings

To add, remove, or modify a user, you can access the user settings by clicking the corresponding icon on the toolbar.

### 9.2 Hardware settings

The access to the hardware settings under **Settings > Hardware** is restricted.

Changing hardware settings must only be performed by Miltenyi Biotec authorized service providers.



# 10

## Important keyboard shortcuts

The ImSpector software has keyboard shortcuts to help the experienced user to use the software as efficiently as possible. The most important shortcuts are summarized in the following table.

| Function                                                                                      | Keyboard shortcut                        |
|-----------------------------------------------------------------------------------------------|------------------------------------------|
| <b>Main Application</b>                                                                       |                                          |
| Shows or hides <b>Live Dialogs</b>                                                            | Tab                                      |
| Toggles Overview mode                                                                         | Shift + Tab                              |
| Shows or hides the <b>Image Series Viewer</b>                                                 | Ctrl + I                                 |
| Shows or hides <b>Display Settings</b>                                                        | Alt + C                                  |
| Creates a new stack window                                                                    | Ctrl + Shift + S                         |
| Creates a new graph window                                                                    | Ctrl + Shift + G                         |
| Opens the software documentation                                                              | F1                                       |
| Full View mode                                                                                | F5                                       |
| <b>Image Data</b>                                                                             |                                          |
| Adjusts the contrast of the colormap automatically to the minimum pixel value in stack window | F9                                       |
| Adjusts the contrast of the colormap automatically to the maximum pixel value in stack window | F10                                      |
| Applies colormap to all tiles                                                                 | Ctrl + Q                                 |
| Toggles fixed colormap scaling                                                                | Ctrl + K                                 |
| Toggles overview MIP                                                                          | F4                                       |
| Zooms in                                                                                      | Ctrl + plus sign key (+) or mouse wheel  |
| Zooms out                                                                                     | Ctrl + minus sign key (-) or mouse wheel |
| Copies colormap settings                                                                      | Ctrl + Alt + C                           |
| Pastes colormap settings                                                                      | Ctrl + Alt + V                           |
| Toggles channel overlay                                                                       | Ctrl + G                                 |
| Shows data stacks in one to nine gallery columns                                              | Ctrl + 1–9                               |
| Moves through ROI history                                                                     | Ctrl + ↑/↓                               |
| Shows data stack properties                                                                   | Ctrl + T                                 |

| Function                                             | Keyboard shortcut            |
|------------------------------------------------------|------------------------------|
| Pans the image                                       | Spacebar + left mouse button |
| Adjusts the maximum level in graph window (Fit best) | F9                           |
| Adjusts the minimum level in graph window (Grow)     | F10                          |

**Table 10.1:** Important keyboard shortcuts

# 11

## Technical support

For technical support, contact your local Miltenyi Biotec representative or Miltenyi Biotec Technical Support at Miltenyi Biotec headquarters:

Miltenyi Biotec B.V. & Co. KG  
Friedrich-Ebert-Straße 68  
51429 Bergisch Gladbach  
Germany  
Phone: +49 2204 8306 3803  
technicalsupport@miltenyi.com

Visit **[www.miltenyibiotec.com](http://www.miltenyibiotec.com)** for local Miltenyi Biotec Technical Support contact information.





# 12

## Legal notes

### 12.1 Limited warranty

Except as stated in a specific warranty statement which may accompany this product or as otherwise agreed in writing by an authorized representative of Miltenyi Biotec, Miltenyi Biotec's warranty to you, the original purchaser and end user ("you" or "your"), with respect to the product accompanied by this limited warranty shall be subject to the following provisions and the general terms and conditions of sale of the company within the Miltenyi Biotec group which supplied the product in effect at the date of purchase. Those terms and conditions of sale may vary by country and region. Nothing in this document should be construed as constituting an additional warranty.

Miltenyi Biotec warrants that this product will operate or perform substantially in conformance with Miltenyi Biotec's published specifications and be free from material defects in material and workmanship, when subjected to normal, proper, and intended usage by properly trained personnel, for the period of time set forth in the product documentation or package inserts accompanying the product (the "Warranty Period").

Miltenyi Biotec agrees, during the Warranty Period, to repair or replace, at Miltenyi Biotec's option, the defective product so as to cause the same to operate in substantial conformance with said published specifications; provided that you shall (a) promptly notify Miltenyi Biotec in writing upon the discovery of any nonconformity or defect, which notice shall include the product model and serial number (if applicable) and details of the warranty claim; and (b) return the nonconforming or defective product to Miltenyi Biotec, freight prepaid, only after receipt of a Return Material Authorization ("RMA") from Miltenyi Biotec, which may include biohazard decontamination procedures and other product-specific handling instructions, if applicable.

Miltenyi Biotec shall have no obligation to make repairs, replacements, or corrections to the product or any component thereof required, in whole or in part, as the result of (i) normal wear and tear, (ii) improper handling, installation, operation, storage, service, maintenance, or repair, (iii) failure to follow the instructions, cautions, warnings, and notes set forth in the product documentation provided with the product or provided by Miltenyi Biotec from time to time, (iv) abnormal use, misuse, neglect, abuse, mishandling, misapplication, modification, or alteration of the product, (v) use of the product in a manner for which it was not designed, (vi) causes external to the product such as, but not limited to, power failure or electrical power surges, (vii) use of the product in combination with equipment, accessories, consumables, or software not supplied or approved by Miltenyi Biotec, or (viii) accident, disaster, or acts of God. ANY INSTALLATION, MAINTENANCE, REPAIR, SERVICE, OR ALTERATION TO OR OF, OR OTHER TAMPERING WITH, THE PRODUCT PERFORMED BY ANY PERSON OR ENTITY OTHER THAN MILTENYI BIOTEC AUTHORIZED PERSONNEL WITHOUT MILTENYI BIOTEC'S PRIOR WRITTEN APPROVAL, OR ANY USE OF REPLACEMENT PARTS NOT SUPPLIED BY MILTENYI BIOTEC, SHALL IMMEDIATELY VOID AND CANCEL ALL WARRANTIES WITH RESPECT TO THE AFFECTED PRODUCT.

Miltenyi Biotec's warranty does not cover products sold AS IS or WITH ALL FAULTS, or which had its serial number defaced, altered, or removed, or any consumables or parts identified as being supplied by a third party.

Miltenyi Biotec must be informed promptly if a claim is made under this limited warranty. If a material or manufacturing defect occurs within the Warranty Period, Miltenyi Biotec will take the appropriate steps, at Miltenyi Biotec's option, to make repairs, replacements, or corrections to the product or any component thereof as may be required to restore the full usability of your product. If Miltenyi Biotec determines that a product for which you have requested warranty services is not covered by the warranty hereunder, you shall pay or reimburse Miltenyi Biotec for all costs of investigating and responding to such request at Miltenyi Biotec's then prevailing time and materials rates. If Miltenyi Biotec provides repair services or replacement parts that are not covered by this warranty, you shall pay Miltenyi Biotec therefor at Miltenyi Biotec's then prevailing time and materials rates.

THE OBLIGATIONS CREATED BY THIS WARRANTY STATEMENT TO REPAIR OR REPLACE A DEFECTIVE PRODUCT ARE EXCLUSIVE AND SHALL BE THE SOLE REMEDY OF BUYER IN THE EVENT OF A DEFECTIVE PRODUCT. EXCEPT AS EXPRESSLY PROVIDED IN THIS WARRANTY STATEMENT, MILTENYI BIOTEC HEREBY DISCLAIMS ALL OTHER WARRANTIES, WHETHER EXPRESS OR IMPLIED, ORAL OR WRITTEN, WITH RESPECT TO THE PRODUCT, INCLUDING WITHOUT LIMITATION ALL IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE. MILTENYI BIOTEC DOES NOT WARRANT THAT THE PRODUCT IS ERROR-FREE OR WILL ACCOMPLISH ANY PARTICULAR RESULT.

UNDER NO CIRCUMSTANCES WILL MILTENYI BIOTEC BE LIABLE FOR ANY LOSS OF USE, INTERRUPTION OF BUSINESS, OR ANY INDIRECT, SPECIAL, INCIDENTAL, PUNITIVE, OR CONSEQUENTIAL DAMAGES OF ANY KIND (INCLUDING LOST PROFITS) REGARDLESS OF THE FORM OF ACTION WHETHER IN CONTRACT, TORT (INCLUDING NEGLIGENCE), STRICT PRODUCT LIABILITY, OR OTHERWISE, EVEN IF MILTENYI BIOTEC HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. IN NO EVENT SHALL THE TOTAL LIABILITY OF MILTENYI BIOTEC HEREUNDER EXCEED THE GREATER OF 100.00 EUROS OR THE AMOUNT YOU ACTUALLY PAID FOR THE PRODUCT GIVING RISE TO SUCH LIABILITY, REGARDLESS OF THE CAUSE OF ACTION, IN CONTRACT, TORT, STRICT LIABILITY, OR OTHERWISE.

NOT ALL JURISDICTIONS ALLOW SUCH LIMITATIONS OF DAMAGES SO THE FOREGOING LIMITATIONS MAY NOT APPLY TO YOU. This warranty statement gives you specific legal rights and you may have other rights, which may vary from jurisdiction to jurisdiction.

## 12.2 Trademarks

Blaze, Choros, MACS, and the Miltenyi Biotec logo are registered trademarks or trademarks of Miltenyi Biotec B.V. & Co. KG and/or its affiliates in various countries worldwide.







## Miltenyi Biotec

### Germany/Austria

Miltenyi Biotec B.V. & Co. KG  
Friedrich-Ebert-Straße 68  
51429 Bergisch Gladbach  
Germany  
Phone +49 2204 8306-0  
Fax +49 2204 85197  
macsde@miltenyi.com

### USA/Canada

Miltenyi Biotec  
North America, Inc.  
1201 Clopper Road  
Gaithersburg,  
MD 20878, USA  
Phone 800 FOR MACS  
Phone +1 866 811 4466  
Fax +1 877 591 1060  
macsau@miltenyi.com

### Australia

Miltenyi Biotec  
Australia Pty. Ltd.  
Unit 11, 2 Eden Park Drive  
Macquarie Park NSW 2113  
Australia  
Phone +61 2 8877 7400  
Fax +61 2 9889 5044  
macsau@miltenyi.com

### Benelux

Miltenyi Biotec B.V.  
Dellaertweg 9C  
2316 WZ Leiden  
The Netherlands  
macsnl@miltenyi.com

### Customer service for:

#### The Netherlands

Phone 0800 4020120  
Fax 0800 4020100

#### Belgium

Phone 0800 94016  
Fax 0800 99626

#### Luxembourg

Phone 800 24971  
Fax 800 24984

### China

Miltenyi Biotec Technology &  
Trading (Shanghai) Co., Ltd.  
Room A401, 4/F  
No. 1077, Zhangheng Road  
Pudong New Area  
201203 Shanghai, P.R. China  
Phone +86 21 6084 0210  
Fax +86 21 6235 0953  
macscn@miltenyi.com.cn

### France

Miltenyi Biotec SAS  
10 rue Mercœur  
75011 Paris, France  
Phone +33 1 56 98 16 16  
macsfr@miltenyi.com

### Hong Kong

Miltenyi Biotec Hong Kong Ltd.  
Unit 301, Lakeside 1  
No. 8 Science Park West Avenue  
Hong Kong Science Park  
Pak Shek Kok, New Territories  
Hong Kong  
Phone +852 3751 6698  
Fax +852 3619 5772  
macshk@miltenyi.com.hk

### India

Miltenyi Biotec India Pvt. Ltd.  
Vatika Business Centre,  
Floor No. 6  
Divya Sree Omega  
Kondapur, Serilingampally  
K.V., Rangareddy  
Telangana 500084, India  
Phone +91 040 45175910  
macsin@miltenyi.com

### Italy

Miltenyi Biotec S.r.l.  
Via Paolo Nanni Costa, 30  
40133 Bologna, Italy  
Phone +39 051 6 460 411  
Fax +39 051 6 460 499  
macsit@miltenyi.com

### Japan

Miltenyi Biotec K.K.  
NEX-Eitai Building 5F  
16-10 Fuyuki, Koto-ku,  
Tokyo 135-0041, Japan  
Phone +81 3 5646 8910  
Fax +81 3 5646 8911  
macsjp@miltenyi.com

### Nordics and Baltics

Miltenyi Biotec Norden AB  
Medicon Village  
Scheeleorget 1  
223 81 Lund, Sweden  
macsse@miltenyi.com

### Customer service for:

#### Sweden

Phone 0200 111 800  
Fax +46 280 72 99

#### Denmark

Phone 80 20 30 10  
Fax +46 46 280 72 99

### Norway, Finland, Iceland, and Baltic countries

Phone +46 46 280 72 80  
Fax +46 46 280 72 99

### Singapore

Miltenyi Biotec Asia Pacific Pte Ltd.  
438B Alexandra Road, Block B  
Alexandra Technopark  
#06-01  
Singapore 119968  
Phone +65 6238 8183  
Fax +65 6238 0302  
macssg@miltenyi.com

### South Korea

Miltenyi Biotec Korea Co., Ltd.  
Donggeuk 7F,  
562 Nonhyeon-ro  
Gangnam-gu  
Seoul 06136, South Korea  
Phone +82 2 555 1988  
Fax +82 2 555 8890  
macskr@miltenyi.com

### Spain

Miltenyi Biotec S.L.  
C/Virgilio 2, Edificio II, Planta -1  
28223 Pozuelo de Alarcón  
Madrid, Spain  
Phone +34 91 512 12 90  
Fax +34 91 512 12 91  
macses@miltenyi.com

### Switzerland

Miltenyi Biotec Swiss AG  
Soodstrasse 52  
8134 Adliswil  
Switzerland  
Phone +41 32 623 08 47  
Fax +49 2204 85197  
macsch@miltenyi.com

### United Kingdom

Miltenyi Biotec Ltd.  
Almac House, Church Lane  
Bisley, Surrey GU24 9DR, UK  
Phone +44 1483 799 800  
Fax +44 1483 799 811  
macsuk@miltenyi.com

[www.miltenyibiotec.com](http://www.miltenyibiotec.com)

Miltenyi Biotec provides products and services worldwide. Visit [www.miltenyibiotec.com/local](http://www.miltenyibiotec.com/local) to find your nearest Miltenyi Biotec contact.

Unless otherwise specifically indicated, Miltenyi Biotec products and services are for research use only and not for therapeutic or diagnostic use.  
Copyright © 2025 Miltenyi Biotec and/or its affiliates. All rights reserved.