

Components	2 mL CD90.1 MicroBeads, mouse and rat: MicroBeads conjugated to monoclonal anti-mouse/rat CD90.1 antibodies (isotype: mouse IgG2a).
Capacity	For 2×10^9 total cells.
Product format	All components are supplied in buffer containing stabilizer and 0.05% sodium azide.
Storage	Store protected from light at 2–8 °C. Do not freeze. The expiration date is indicated on the vial labels.

Safety information

For research use only. Not intended for any animal or human therapeutic or diagnostic use.

Before use, please consult the Safety Data Sheet for information regarding hazards and safe handling practices.

Cell separation protocols






Fully automated cell labeling and separation using the autoMACS® Pro Separator

Alternatively:

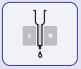

Manual magnetic labeling

↓  Subsequent manual separation	or	↓  Subsequent semi-automated cell separation using the MultiMACS Cell24 Separator Plus	or	↓  Subsequent automated cell separation using the autoMACS® Pro Separator
---	----	--	----	---

General notes

▲ For tips concerning sample preparation (e.g. with the gentleMACS™ Dissociator), magnetic labeling and separation, visit www.miltenyibiotec.com/faq and www.miltenyibiotec.com/protocols.

▲ For product-specific background information and applications of this product, refer to the respective product page at www.miltenyibiotec.com/130-121-273.

Reagent and instrument requirements

- Buffer: Prepare a solution containing phosphate-buffered saline (PBS), pH 7.2, 0.5% bovine serum albumin (BSA), and 2 mM EDTA by diluting MACS® BSA Stock Solution (# 130-091-376) 1:20 with autoMACS® Rinsing Solution (# 130-091-222). Degas buffer before use, as air bubbles could block the column.

- (Optional) Pre-Separation Filters (30 µm) (# 130-041-407) to remove cell clumps.
- Choose the appropriate MACS Separator and MACS Columns:

Column	Max. number of labeled cells	Max. number of total cells	Separator
Positive selection			
MS	10^7	2×10^8	MiniMACS, OctoMACS
LS	10^8	2×10^9	MidiMACS, QuadroMACS
Multi-24 Column Block (per column)	10^8	10^9	MultiMACS Cell24 Separator Plus
Depletion			
LD	10^8	5×10^8	MidiMACS, QuadroMACS, VarioMACS, SuperMACS II
Positive selection or depletion			
autoMACS	2×10^8	4×10^9	autoMACS Pro

▲ **Note:** For depletion the use of LD Columns is strongly recommended. For details refer to the LD Column data sheet.



Fully automated cell labeling and separation using the autoMACS® Pro Separator

▲ Refer to the user manual for instructions on how to use the autoMACS® Pro Separator.

▲ All buffer temperatures should be ≥ 10 °C.

▲ Place tubes in the following Chill Rack positions:

position A = sample, **position B** = negative fraction, **position C** = positive fraction.

1. For appropriate resuspension volumes and cell concentrations, please visit www.automacspro.com/autolabeling.
2. Switch on the instrument for automatic initialization.
3. Go to the **Reagent** menu and select **Read Reagent**. Scan the 2D barcode of each reagent vial with the barcode scanner on the autoMACS Pro Separator. Place the reagent into the appropriate position on the reagent rack.
4. Place sample and collection tubes into the Chill Rack.
5. Go to the **Separation** menu and select the reagent name for each sample from the **Labeling** submenu (the correct labeling, separation (**Possel**), and wash protocols will be selected automatically for **positive selection** by default).
 - ▲ **Note:** If depletion of target cells is desired, please select program **Depletes** from the menu.
 - ▲ **Note:** Choose program **Possels** if higher recovery is required at costs of purity or contact technical support macstec@miltenyibiotec.de.
6. Enter sample volume into the **Volume** submenu. Press **Enter**.
7. Select **Run**.
8. Collect enriched CD90.1⁺ T cell fraction at position C = positive fraction.

9. (Optional) Collect negative cell fraction at position B containing unlabeled cells that are depleted from CD90.1⁺ T cells.



Manual magnetic labeling

- ▲ Work fast, keep cells cold, and use pre-cooled solutions (2–8 °C).
- ▲ Volumes for magnetic labeling given below are for up to 10⁷ total cells. When working with fewer cells, use the same volumes as indicated. When working with higher cell numbers, scale up all reagent volumes and total volumes accordingly.
- ▲ For optimal performance it is important to obtain a single-cell suspension before magnetic labeling.

1. Prepare cells and determine cell number.
2. Resuspend cell pellet in 90 µL of buffer per 10⁷ total cells.
3. Add 10 µL of CD90.1 MicroBeads per 10⁷ total cells.
4. Mix well and incubate for 10 minutes in the refrigerator (2–8 °C).
5. Proceed to subsequent magnetic cell separation.
 - ▲ **Note:** A minimum of 500 µL is required for magnetic separation. If necessary, add buffer to the cell suspension.



Subsequent semi-automated cell separation using the MultiMACS™ Cell24 Separator Plus

The MultiMACS™ Cell24 Separator Plus can be used with the Multi-24 Column Block.

6. Prepare and prime the instrument.
7. Follow instructions given on the Touch Screen Display and in the respective user manual.
8. Choose program “**Possel2**”. Collect enriched cells according to respective user manual.



Subsequent manual cell separation

- ▲ Always wait until the column reservoir is empty before proceeding to the next step.

6. Place column in the magnetic field of a suitable MACS Separator. For details refer to the respective MACS Column data sheet.
7. Prepare column by rinsing with the appropriate amount of buffer:

MS: 500 µL	LS: 3 mL
------------	----------
8. Apply cell suspension onto the column. Collect flow-through containing unlabeled cells, representing the CD90.1[−] cell fraction.
9. Wash column with the appropriate amount of buffer. Collect unlabeled cells that pass through and combine with the effluent from step 8.

MS: 2×500 µL	LS: 1×3 mL
--------------	------------
10. Remove column from the separator and place it on a suitable collection tube. Pipette the appropriate amount of buffer onto the column. Immediately flush out the magnetically labeled CD90.1⁺ T cells by firmly pushing the plunger into the column.

MS: 1 mL	LS: 5 mL
----------	----------

11. (Optional) To increase the purity of CD90.1⁺ T cells, the eluted fraction can be enriched over a second MS or LS Column. Repeat the magnetic separation procedure as described in steps 6 to 10 by using a new column.



Subsequent automated cell separation using the autoMACS® Pro Separator

- ▲ Refer to the user manual for instructions on how to use the autoMACS® Pro Separator.

- ▲ All buffer temperatures should be ≥10 °C.
- ▲ Place tubes in the following Chill Rack positions:
position A = sample, **position B** = negative fraction,
position C = positive fraction.

6. Prepare and prime the instrument.
7. Follow the instructions that are given in the user manual.
8. For a standard separation choose one of the following programs:

Positive selection: Possel

Collect positive fraction from row C of the tube rack.


- ▲ **Note:** Choose program **Possels** if higher recovery is required at costs of purity or contact technical support macstec@miltenyibiotec.de.

Depletion: Depletes

Collect negative fraction in row B of the tube rack.

Example of a separation using CD90.1 MicroBeads

A single-cell suspension from mouse spleen was prepared using the gentleMACS™ Dissociator. CD90.1⁺ T cells were isolated from this single-cell suspension using the CD90.1 MicroBeads and the autoMACS® Pro Separator. Cells were fluorescently stained with REAfinity™ Antibodies CD45-VioGreen™, CD45R-APC-Vio® 770 and CD90.1-VioBlue® and analyzed by flow cytometry using the MACSQuant® Analyzer. Viable leukocytes were gated for analysis based on scatter signals, 7-AAD Staining Solution fluorescence, and CD45 expression.

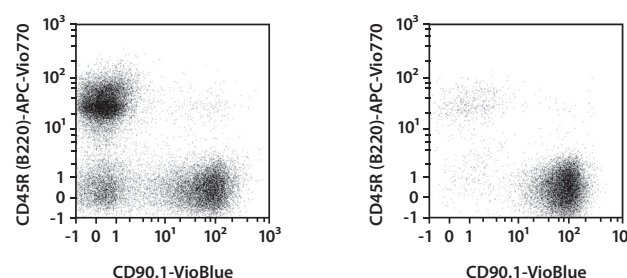

For more information or assistance refer to our technical support.

Check out Miltenyi Biotec's flow cytometry solutions at www.miltenyibiotec.com/MACSQuant and explore the extensive antibody portfolio at www.miltenyibiotec.com/antibodies.

Refer to www.miltenyibiotec.com for all data sheets and protocols. Miltenyi Biotec provides technical support worldwide. Visit www.miltenyibiotec.com/local to find your nearest Miltenyi Biotec contact.

Legal notices

Limited product warranty

Miltenyi Biotec B.V. & Co. KG and/or its affiliate(s) warrant this product to be free from material defects in workmanship and materials and to conform substantially with Miltenyi Biotec's published specifications for the product at the time of order, under normal use and conditions in accordance with its applicable documentation, for a period beginning on the date of delivery of the product by Miltenyi Biotec or its authorized distributor and ending on the expiration date of the product's applicable shelf life stated on the product label, packaging or documentation (as applicable) or, in the absence thereof, ONE (1) YEAR from date of delivery ("Product Warranty"). Miltenyi Biotec's Product Warranty is provided subject to the warranty terms as set forth in Miltenyi Biotec's General Terms and Conditions for the Sale of Products and Services available on Miltenyi Biotec's website at www.miltenyibiotec.com, as in effect at the time of order ("Product Warranty"). Additional terms may apply. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS. THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING IF A PRODUCT IS SUITABLE FOR CUSTOMER'S PARTICULAR PURPOSE AND APPLICATION METHODS.

Technical information

The technical information, data, protocols, and other statements provided by Miltenyi Biotec in this document are based on information, tests, or experience which Miltenyi Biotec believes to be reliable, but the accuracy or completeness of such information is not guaranteed. Such technical information and data are intended for persons with knowledge and technical skills sufficient to assess and apply their own informed judgment to the information. Miltenyi Biotec shall not be liable for any technical or editorial errors or omissions contained herein.

All information and specifications are subject to change without prior notice. Please contact Miltenyi Biotec Technical Support or visit www.miltenyibiotec.com for the most up-to-date information on Miltenyi Biotec products.

Licenses

This product and/or its use may be covered by one or more pending or issued patents and/or may have certain limitations. Certain uses may be excluded by separate terms and conditions. Please contact your local Miltenyi Biotec representative or visit Miltenyi Biotec's website at www.miltenyibiotec.com for more information.

The purchase of this product conveys to the customer the non-transferable right to use the purchased amount of the product in research conducted by the customer (whether the customer is an academic or for-profit entity). This product may not be further sold. Additional terms and conditions (including the terms of a Limited Use Label License) may apply.

CUSTOMER'S USE OF THIS PRODUCT MAY REQUIRE ADDITIONAL LICENSES DEPENDING ON THE SPECIFIC APPLICATION. THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING FOR ITSELF WHETHER IT HAS ALL APPROPRIATE LICENSES IN PLACE. Miltenyi Biotec provides no warranty that customer's use of this product does not and will not infringe intellectual property rights owned by a third party. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

Trademarks

autoMACS, gentleMACS, MACS, MACSQuant, MidiMACS, the Miltenyi Biotec logo, MiniMACS, MultiMACS, OctoMACS, QuadroMACS, REAfinity, SuperMACS, VarioMACS, Vio, VioBlue, and VioGreen are registered trademarks or trademarks of Miltenyi Biotec and/or its affiliates in various countries worldwide.

Copyright © 2020 Miltenyi Biotec and/or its affiliates. All rights reserved.