



Miltenyi Biotec

# Pre-Separation Filters

20  $\mu\text{m}$   
30  $\mu\text{m}$   
70  $\mu\text{m}$

130-101-812  
130-041-407  
130-095-823

## Contents

### 1. Description

- 1.1 Important safety information
- 1.2 Background information
- 1.3 Applications
- 1.4 Technical specifications


## 1. Description

This product is for research use only.

For use by professional laboratory personnel only.

Components	<b>50 Pre-Separation Filters (20 <math>\mu\text{m}</math>)</b> Green colored, sterile packed or
	<b>50 Pre-Separation Filters (30 <math>\mu\text{m}</math>)</b> Yellow colored, sterile packed or
	<b>50 Pre-Separation Filters (70 <math>\mu\text{m}</math>)</b> Purple colored, sterile packed

**Storage** Store Pre-Separation Filters dry. The expiration date is indicated on the box label.

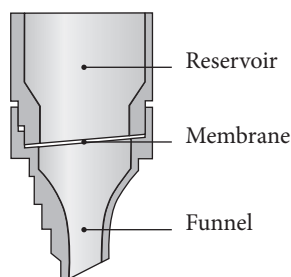

### 1.1 Important safety information

#### Warning

All biological material must be considered potentially infectious. Regulations for the treatment and disposal of infectious materials must be observed.

### 1.2 Background information

Pre-Separation Filters are designed for the easy removal of cell aggregates or large particles from single-cell suspensions to ensure effective magnetic or fluorescent cell labeling, and thus optimal flow within cell separation columns and flow cytometer. They fit onto MACS® MS, LS, LD, CS, Large Cell, and Whole Blood Columns, as well as onto tubes with a capacity of 5 mL and 13–15 mL.

One set contains 50 sterile and ready-to-use filters.

### 1.3 Applications

- Removal of large particles from blood cell samples and from single-cell suspensions prior to fluorescent or magnetic cell labeling, cell separation, flow cytometric analysis, or cell sorting.

### 1.4 Technical specifications

- For single use only.
- Mesh size: 20  $\mu\text{m}$ , 30  $\mu\text{m}$ , or 70  $\mu\text{m}$ , respectively.
- Recommended sample volume: 500–3500  $\mu\text{L}$

▲ The membrane of the Pre-Separation Filter is slightly hydrophobic. Moisturize Pre-Separation Filter by pipetting 500  $\mu\text{L}$  of buffer into the reservoir and discard effluent. Use same buffer as used to resuspend cells.

▲ Avoid drops trapped in the funnel of the filter. Tap the rim of the filter to elute the suspension completely.

Refer to [www.miltenyibiotec.com](http://www.miltenyibiotec.com) for all data sheets and protocols. Miltenyi Biotec provides technical support worldwide. Visit [www.miltenyibiotec.com](http://www.miltenyibiotec.com) for local Miltenyi Biotec Technical Support contact information.

## Legal notices

### Limited product warranty

Miltenyi Biotec B.V. & Co. KG and/or its affiliate(s) warrant this product to be free from material defects in workmanship and materials and to conform substantially with Miltenyi Biotec's published specifications for the product at the time of order, under normal use and conditions in accordance with its applicable documentation, for a period beginning on the date of delivery of the product by Miltenyi Biotec or its authorized distributor and ending on the expiration date of the product's applicable shelf life stated on the product label, packaging or documentation (as applicable) or, in the absence thereof, ONE (1) YEAR from date of delivery ("Product Warranty"). Miltenyi Biotec's Product Warranty is provided subject to the warranty terms as set forth in Miltenyi Biotec's General Terms and Conditions for the Sale of Products and Services available on Miltenyi Biotec's website at [www.miltenyibiotec.com](http://www.miltenyibiotec.com), as in effect at the time of order ("Product Warranty"). Additional terms may apply. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS. THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING IF A PRODUCT IS SUITABLE FOR CUSTOMER'S PARTICULAR PURPOSE AND APPLICATION METHODS.

### Technical information

The technical information, data, protocols, and other statements provided by Miltenyi Biotec in this document are based on information, tests, or experience which Miltenyi Biotec believes to be reliable, but the accuracy or completeness of such information is not guaranteed. Such technical information and data are intended for persons with knowledge and technical skills sufficient to assess and apply their own informed judgment to the information. Miltenyi Biotec shall not be liable for any technical or editorial errors or omissions contained herein.

All information and specifications are subject to change without prior notice. Please contact Miltenyi Biotec Technical Support or visit [www.miltenyibiotec.com](http://www.miltenyibiotec.com) for the most up-to-date information on Miltenyi Biotec products.

### Licenses

This product and/or its use may be covered by one or more pending or issued patents and/or may have certain limitations. Certain uses may be excluded by separate terms and conditions. Please contact your local Miltenyi Biotec representative or visit Miltenyi Biotec's website at [www.miltenyibiotec.com](http://www.miltenyibiotec.com) for more information.

The purchase of this product conveys to the customer the non-transferable right to use the purchased amount of the product in research conducted by the customer (whether the customer is an academic or for-profit entity). This product may not be further sold. Additional terms and conditions (including the terms of a Limited Use Label License) may apply.

CUSTOMER'S USE OF THIS PRODUCT MAY REQUIRE ADDITIONAL LICENSES DEPENDING ON THE SPECIFIC APPLICATION. THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING FOR ITSELF WHETHER IT HAS ALL APPROPRIATE LICENSES IN PLACE. Miltenyi Biotec provides no warranty that customer's use of this product does not and will not infringe intellectual property rights owned by a third party. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

### Trademarks

MACS and the Miltenyi Biotec logo are registered trademarks or trademarks of Miltenyi Biotec and/or its affiliates in various countries worldwide.

Copyright © 2022 Miltenyi Biotec and/or its affiliates. All rights reserved.