

Contents

1. Description
 - 1.1 Principle of the MACS® Separation
 - 1.2 Background information
 - 1.3 Applications
 - 1.4 Reagent and instrument requirements
2. Protocol
 - 2.1 Sample preparation
 - 2.2 Magnetic labeling
 - 2.3 Magnetic separation
3. Example of a separation using the Anti-Ly-6G MicroBeads UltraPure

Warnings

Reagents contain sodium azide. Under acidic conditions sodium azide yields hydrazoic acid, which is extremely toxic. Azide compounds should be diluted with running water before discarding. These precautions are recommended to avoid deposits in plumbing where explosive conditions may develop.

1. Description

This product is for research use only.

Components	2 mL Anti-Ly-6G MicroBeads UltraPure, mouse: MicroBeads conjugated to monoclonal anti-mouse Ly-6G antibodies (isotype: human IgG1).
Capacity	For 2×10^9 total cells, up to 200 separations.
Product format	Anti-Ly-6G MicroBeads UltraPure are supplied in buffer containing stabilizer and 0.05% sodium azide.
Storage	Store protected from light at 2–8 °C. Do not freeze. The expiration date is indicated on the vial label.

1.1 Principle of the MACS® Separation

First, the Ly-6G⁺ cells are magnetically labeled with Anti-Ly-6G MicroBeads UltraPure. Then, the cell suspension is loaded onto a MACS® Column, which is placed in the magnetic field of a MACS Separator. The magnetically labeled Ly-6G⁺ cells are retained within the column. The unlabeled cells run through; this cell fraction is thus depleted of Ly-6G⁺ cells. After removing the column from the magnetic field, the magnetically retained Ly-6G⁺ cells can be eluted as the positively selected cell fraction. To increase the purity of Ly-6G⁺ cells from lung and blood, the positively selected cell fraction containing the Ly-6G⁺ cells must be separated over a second column.

1.2 Background information

The Anti-Ly-6G MicroBeads Ultrapure were developed for direct, positive selection or depletion of mouse neutrophils from lymphoid tissues. Ly-6G is highly expressed on neutrophils, at lower level on a subset of eosinophils, and transiently during developmental stages on monocytes. This kit contains Anti-Ly-6G clone REA526, which displays negligible binding to Fc receptors.

1.3 Applications

- Positive selection or depletion of cells expressing human Ly-6G antigen from mouse bone marrow, lung, and blood.

1.4 Reagent and instrument requirements

- Buffer: Prepare a solution containing phosphate-buffered saline (PBS), pH 7.2, 0.5% bovine serum albumin (BSA), and 2 mM EDTA by diluting MACS® BSA Stock Solution (# 130-091-376) 1:20 with autoMACS® Rinsing Solution (# 130-091-222). Keep buffer cold (2–8 °C). Degas buffer before use, as air bubbles could block the column.

▲ **Note:** EDTA can be replaced by other supplements such as anticoagulant citrate dextrose formula-A (ACD-A) or citrate phosphate dextrose (CPD). BSA can be replaced by other proteins such as mouse serum albumin, mouse serum, or fetal bovine serum (FBS). Buffers or media containing Ca²⁺ or Mg²⁺ are not recommended for use.

- MACS Columns and MACS Separators: Ly-6G⁺ cells can be enriched by using MS or LS Columns or depleted with the use of LD Columns. Positive selection can also be performed by using the autoMACS Pro Separator or the MultiMACS™ Cell24 Separator Plus.

Column	Max. number of labeled cells	Max. number of total cells	Separator
Positive selection			
MS	10 ⁷	2 × 10 ⁸	MiniMACS, OctoMACS, VarioMACS, SuperMACS II
LS	10 ⁸	2 × 10 ⁹	MidiMACS, QuadroMACS, VarioMACS, SuperMACS II,
	10 ⁸	10 ⁹	MultiMACS Cell24 Separator Plus
autoMACS	2 × 10 ⁸	4 × 10 ⁹	autoMACS Pro
Multi-24 Column Block (per column)	10 ⁸	10 ⁹	MultiMACS Cell24 Separator Plus
Depletion			
LD	10 ⁸	5 × 10 ⁸	MidiMACS, QuadroMACS, VarioMACS, SuperMACS II, MultiMACS Cell24 Separator Plus

▲ **Note:** Column adapters are required to insert certain columns into the VarioMACS™ or SuperMACS™ II Separators. For details refer to the respective MACS Separator data sheet.

▲ **Note:** If separating with LS Columns and the MultiMACS Cell24 Separator Plus use the Single-Column Adapter. Refer to the user manual for details.

- (Optional) Fluorochrome-conjugated antibodies for flow cytometric analysis, e.g., Anti-Ly-6G-VioBlue® or Anti-Ly-6G-APC. For more information about antibodies refer to www.miltenyibiotec.com/antibodies.
- (Optional) Propidium Iodide Solution (# 130-093-233) or 7-AAD Staining Solution (# 130-111-568) for flow cytometric exclusion of dead cells.
- (Optional) Dead Cell Removal Kit (# 130-090-101) for the depletion of dead cells.
- (Optional) Pre-Separation Filters (30 µm) (# 130-041-407) to remove cell clumps.

2. Protocol

2.1 Sample preparation

▲ For highest recovery and purity of Ly-6G⁺ neutrophils from mouse lung, it is recommended to perform enzymatic disaggregation using the Lung Dissociation Kit, mouse (# 130-095-927).

▲ When working with solid tissues, prepare a single-cell suspension using the gentleMACS™ Dissociator and an appropriate tissue dissociation kit.

▲ If working with bone marrow, blood or lung single-cell suspension please refer to Magnetic labeling 2.2. When working with blood, erythrocyte lysis must be performed prior to magnetic labeling by using the Red Blood Cell Lysis Solution (10× (# 130-094-183)).

▲ Dead cells may bind non-specifically to MACS MicroBeads. To remove dead cells, we recommend using density gradient centrifugation or the Dead Cell Removal Kit (# 130-090-101).



2.2 Magnetic labeling

▲ Work fast, keep cells cold, and use pre-cooled solutions. This will prevent capping of antibodies on the cell surface and non-specific cell labeling.

▲ Volumes for magnetic labeling given below are for up to 10⁷ total cells. When working with fewer than 10⁷ cells, use the same volumes as indicated. When working with higher cell numbers, scale up all reagent volumes and total volumes accordingly (e.g. for 2×10⁷ total cells, use twice the volume of all indicated reagent volumes and total volumes).

▲ For optimal performance it is important to obtain a single-cell suspension before magnetic labeling. Pass cells through 30 µm nylon mesh (Pre-Separation Filters (30 µm), # 130-041-407) to remove cell clumps which may clog the column. Moisten filter with buffer before use.

▲ The recommended incubation temperature is 2–8 °C. Higher temperatures and/or longer incubation times may lead to non-specific cell labeling. Working on ice may require increased incubation times.

1. Determine cell number.
2. Centrifuge cell suspension at 300×g for 10 minutes. Aspirate supernatant completely.
3. Resuspend cell pellet in 90 µL of buffer per 10⁷ total cells.
4. Add 10 µL of Anti-Ly-6G MicroBeads UltraPure per 10⁷ total cells.
5. Mix well and incubate for 10 minutes in the refrigerator (2–8 °C).
6. Wash cells by adding 1–2 mL of buffer per 10⁷ cells and centrifuge at 300×g for 10 minutes. Aspirate supernatant completely.
7. Resuspend up to 10⁸ cells in 500 µL of buffer.
▲ **Note:** For higher cell numbers, scale up buffer volume accordingly.
8. Proceed to magnetic separation (2.3).



2.3 Magnetic separation

▲ Choose an appropriate MACS Column and MACS Separator according to the number of total cells and the number of Ly-6G⁺ cells. For details refer to the table in section 1.4.

▲ Always wait until the column reservoir is empty before proceeding to the next step.

Magnetic separation with MS or LS Columns

1. Place column in the magnetic field of a suitable MACS Separator. For details refer to the respective MACS Column data sheet.
2. Prepare column by rinsing with the appropriate amount of buffer:

MS: 500 µL LS: 3 mL

3. Apply cell suspension onto the column. Collect flow-through containing unlabeled cells.
4. Wash column with the appropriate amount of buffer. Collect unlabeled cells that pass through and combine with the flow-through from step 3.

MS: 3×500 µL LS: 3×3 mL

▲ **Note:** Perform washing steps by adding buffer aliquots as soon as the column reservoir is empty.

5. Remove column from the separator and place it on a suitable collection tube.
6. Pipette the appropriate amount of buffer onto the column. Immediately flush out the magnetically labeled cells by firmly pushing the plunger into the column.

MS: 1 mL LS: 5 mL

7. (Optional) To increase the purity of Ly-6G⁺ cells isolated from lung, the eluted fraction must be enriched over a second MS or LS Column. Repeat the magnetic separation procedure as described in steps 1 to 6 by using a new column.

Magnetic separation with the MultiMACS™ Cell24 Separator

Refer to the the MultiMACS™ Cell Separator user manual for instructions on how to use the MultiMACS Cell24 Separator.

Magnetic separation with the autoMACS® Pro Separator

▲ Refer to the respective user manual for instructions on how to use the autoMACS® Pro Separator.

▲ Buffers used for operating the autoMACS Pro Separator should have a temperature of $\geq 10^\circ\text{C}$.

1. Prepare and prime the instrument.
2. Apply tube containing the sample and provide tubes for collecting the labeled and unlabeled cell fractions. Place sample tube in row A of the tube rack and the fraction collection tubes in rows B and C.
3. For a standard separation choose the following program:

Positive selection from blood or bone marrow: Possel

Collect positive fraction in row C of the tube rack.


Positive selection from lung: Posseld2

Collect positive fraction in row C of the tube rack.


3. Example of a separation using Anti-Ly-6G MicroBeads UltraPure

Ly-6G⁺ neutrophils were isolated from mouse blood (A), lung (B), or bone marrow (C) using Anti-Ly-6G MicroBeads UltraPure, one or two MS Columns (for bone marrow and blood or lung, respectively), and a MiniMACS™ Separator. Cells were fluorescently stained with Anti-Ly-6G-VioBlue® and Anti-Ly-6C-APC and analyzed by flow cytometry using the MACSQuant® X. Cell debris and dead cells were excluded from the analysis based on scatter signals and propidium iodide fluorescence.


A) Blood

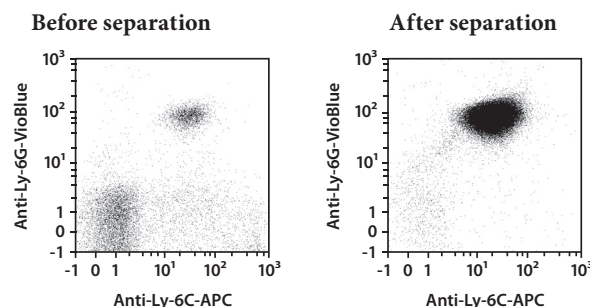

B) Bone marrow

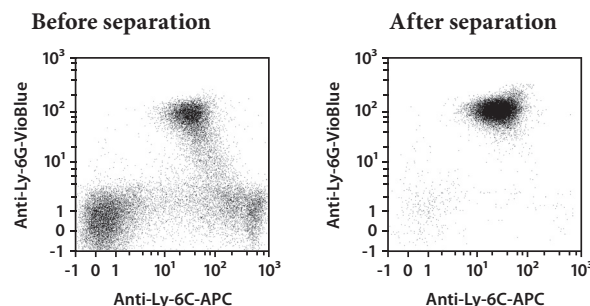

C) Lung

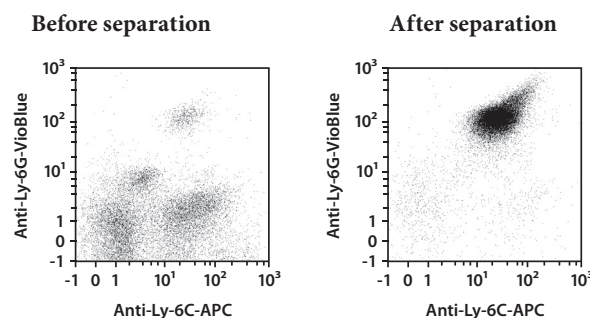

Refer to www.miltenyibiotec.com for all data sheets and protocols. Miltenyi Biotec provides technical support worldwide. Visit www.miltenyibiotec.com/local to find your nearest Miltenyi Biotec contact.

Legal notices

Limited product warranty

Miltenyi Biotec B.V. & Co. KG and/or its affiliate(s) warrant this product to be free from material defects in workmanship and materials and to conform substantially with Miltenyi Biotec's published specifications for the product at the time of order, under normal use and conditions in accordance with its applicable documentation, for a period beginning on the date of delivery of the product by Miltenyi Biotec or its authorized distributor and ending on the expiration date of the product's applicable shelf life stated on the product label, packaging or documentation (as applicable) or, in the absence thereof, ONE (1) YEAR from date of delivery ("Product Warranty"). Miltenyi Biotec's Product Warranty is provided subject to the warranty terms as set forth in Miltenyi Biotec's General Terms and Conditions for the Sale of Products and Services available on Miltenyi Biotec's website at www.miltenyibiotec.com, as in effect at the time of order ("Product Warranty"). Additional terms may apply. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING IF A PRODUCT IS SUITABLE FOR CUSTOMER'S PARTICULAR PURPOSE AND APPLICATION METHODS.

Technical information

The technical information, data, protocols, and other statements provided by Miltenyi Biotec in this document are based on information, tests, or experience which Miltenyi Biotec believes to be reliable, but the accuracy or completeness of such information is not guaranteed. Such technical information and data are intended for persons with knowledge and technical skills sufficient to assess and apply their own informed judgment to the information. Miltenyi Biotec shall not be liable for any technical or editorial errors or omissions contained herein.

All information and specifications are subject to change without prior notice. Please contact Miltenyi Biotec Technical Support or visit www.miltenyibiotec.com for the most up-to-date information on Miltenyi Biotec products.

Licenses

This product and/or its use may be covered by one or more pending or issued patents and/or may have certain limitations. Certain uses may be excluded by separate terms and conditions. Please contact your local Miltenyi Biotec representative or visit Miltenyi Biotec's website at www.miltenyibiotec.com for more information.

The purchase of this product conveys to the customer the non-transferable right to use the purchased amount of the product in research conducted by the customer (whether the customer is an academic or for-profit entity). This product may not be further sold. Additional terms and conditions (including the terms of a Limited Use Label License) may apply.

CUSTOMER'S USE OF THIS PRODUCT MAY REQUIRE ADDITIONAL LICENSES DEPENDING ON THE SPECIFIC APPLICATION. THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING FOR ITSELF WHETHER IT HAS ALL APPROPRIATE LICENSES IN PLACE. Miltenyi Biotec provides no warranty that customer's use of this product does not and will not infringe intellectual property rights owned by a third party. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

Trademarks

autoMACS, gentleMACS, MACS, MACSQuant, MidiMACS, the Miltenyi Biotec logo, MiniMACS, MultiMACS, OctoMACS, QuadroMACS, SuperMACS, VarioMACS, and VioBlue are registered trademarks or trademarks of Miltenyi Biotec and/or its affiliates in various countries worldwide.

Copyright © 2020 Miltenyi Biotec and/or its affiliates. All rights reserved.