

2 mL Streptavidin MicroBeads	Order no. 130-048-101
1 mL Streptavidin MicroBeads	Order no. 130-048-102

## Index

1. Description
  - 1.1 Principle of MACS® Separation
  - 1.2 Background and product applications
  - 1.3 Reagent and instrument requirements
2. Protocol
  - 2.1 Sample preparation
  - 2.2 Magnetic labeling
  - 2.3 Magnetic separation
3. Example of a separation using Streptavidin MicroBeads
4. References

## 1. Description

<b>Components</b>	2 mL Streptavidin MicroBeads (# 130-048-101) or 1 mL Streptavidin MicroBeads (# 130-048-102): MicroBeads conjugated to streptavidin.
<b>Size</b>	2 mL for 2×10 <sup>9</sup> total cells, up to 200 separations or 1 mL for 10 <sup>9</sup> total cells, up to 100 separations.
<b>Product format</b>	Streptavidin MicroBeads are supplied in a solution containing stabilizer and 0.05% sodium azide.
<b>Storage</b>	Store protected from light at 2–8 °C. Do not freeze. The expiration date is indicated on the vial label.

### 1.1 Principle of MACS® Separation

First, the cells are stained with a biotinylated primary antibody or ligand. Subsequently, the cells are magnetically labeled with Streptavidin MicroBeads. Then the cell suspension is loaded onto a MACS® Column which is placed in the magnetic field of a MACS Separator. The magnetically labeled cells are retained in the column while the unlabeled cells run through. After removing the column from the magnetic field, the magnetically retained cells can be eluted as the positively selected cell fraction.

### 1.2 Background and product applications

Streptavidin MicroBeads are developed for the separation of cells according to surface markers labeled with biotinylated antibodies or ligands.

Streptavidin MicroBeads have been utilized for the isolation of human cells such as isotype specific<sup>1</sup> or allergen specific<sup>2</sup> B cells, epithelial cells<sup>3</sup>, follicular dendritic cells<sup>4,5,6</sup>, or stromal cells<sup>7</sup>. They have also been used for the isolation of mouse cells such as MHC class I<sup>+</sup> cells<sup>8</sup>, isotype specific B cells<sup>9</sup>, liver sinusoidal endothelial cells<sup>10</sup>, or for the sequential isolation of mouse dendritic cells,

macrophages, and medullary and cortical epithelial cells<sup>11</sup>. Further, Streptavidin MicroBeads have been used for the isolation of rat T cells<sup>12</sup> and rat retinal ganglion cells<sup>13</sup>, rhesus monkey T cells<sup>14</sup>, sheep B cell subpopulations<sup>15</sup>, fish (carp) B cells<sup>16</sup> or plant (potato) protoplasts<sup>17</sup>.

### Example applications

- Positive selection or depletion of cells labeled with biotinylated antibodies.
- Positive selection or depletion of cells labeled with biotinylated ligands.

### 1.3 Reagent and instrument requirements

- **Labeling buffer:** Prepare a solution containing phosphate buffered saline (PBS) pH 7.2, supplemented with 2 mM EDTA, optionally also supplemented with 0.5% **biotin-free** bovine serum albumin (BSA). Keep buffer cold (4–8 °C). Degas buffer before use, as air bubbles could block the column.

▲ **Note:** For magnetic labeling, use **biotin-free** labeling buffer only, for example, buffer without BSA or buffer supplemented with biotin-free BSA. If the BSA preparation is contaminated with biotin, free biotin competes with biotinylated antibody for the binding of the Streptavidin MicroBeads. The use of cell culture medium is not recommended as it may contain free biotin. Alternatively, use Anti-Biotin MicroBeads (# 130-090-485) which do not bind to free biotin.

- **Separation buffer:** Prepare a solution containing PBS pH 7.2, 0.5% BSA, and 2 mM EDTA by diluting MACS BSA Stock Solution (# 130-091-376) 1:20 with autoMACS™ Rinsing Solution (# 130-091-222). Keep buffer cold (4–8 °C). Degas buffer before use, as air bubbles could block the column.

▲ **Note:** EDTA can be replaced by other supplements such as anticoagulant citrate dextrose formula-A (ACD-A) or citrate phosphate dextrose (CPD). BSA can be replaced by other proteins such as human serum albumin, human serum, or fetal calf serum. Buffers or media containing Ca<sup>2+</sup> or Mg<sup>2+</sup> are not recommended for use.

- MACS® Columns and MACS Separators: Cells labeled with Streptavidin MicroBeads can be enriched by using MS, LS or XS Columns (positive selection). Streptavidin MicroBeads can be used for depletion of cells on LD, CS or D Columns. Cells which strongly express the biotin labeled antigen can also be depleted using MS, LS or XS Columns. Positive selection or depletion can also be performed by using the autoMACS™ Separator.

Column	Max. number of labeled cells	Max. number of total cells	Separator
<b>Positive selection</b>			
MS	10 <sup>7</sup>	2×10 <sup>8</sup>	MiniMACS, OctoMACS, VarioMACS, SuperMACS
LS	10 <sup>8</sup>	2×10 <sup>9</sup>	MidiMACS, QuadroMACS, VarioMACS, SuperMACS
XS	10 <sup>9</sup>	2×10 <sup>10</sup>	SuperMACS
<b>Depletion</b>			
LD	10 <sup>8</sup>	5×10 <sup>8</sup>	MidiMACS, QuadroMACS, VarioMACS, SuperMACS
CS	2×10 <sup>8</sup>		VarioMACS, SuperMACS
D	10 <sup>9</sup>		SuperMACS
<b>Positive selection or depletion</b>			
autoMACS	2×10 <sup>8</sup>	4×10 <sup>9</sup>	autoMACS

▲ **Note:** Column adapters are required to insert certain columns into VarioMACS™ or SuperMACS™ Separators. For details see the respective MACS Separator data sheet.

- Biotinylated primary antibody or ligand.
- (Optional) Anti-Biotin-FITC (# 130-090-857), Anti-Biotin-PE (# 130-090-756), Anti-Biotin-APC (# 130-090-856) or fluorochrome conjugated streptavidin for flow cytometric analysis.
- (Optional) Propidium iodide (PI) or 7-AAD for flow cytometric exclusion of dead cells.
- (Optional) Pre-Separation Filters (# 130-041-407) to remove cell clumps.
- (Optional) Dead Cell Removal Kit (# 130-090-101) for the labeling and discrimination of dead cells by flow cytometry.

## 2. Protocol


### 2.1 Sample preparation

When working with anticoagulated peripheral blood or buffy coat, peripheral blood mononuclear cells (PBMCs) should be isolated by density gradient centrifugation, e.g., using Ficoll-Paque™. For details see section General Protocols in the user manuals or visit [www.miltenyibiotec.com/protocols](http://www.miltenyibiotec.com/protocols).

▲ **Note:** To remove platelets after density gradient separation, resuspend cell pellet in buffer and centrifuge at 200×g for 10–15 minutes at 20 °C. Carefully aspirate supernatant. Repeat washing step.

When working with tissues, prepare a single-cell suspension by a standard preparation method. For details see section General Protocols in the user manuals or visit [www.miltenyibiotec.com/protocols](http://www.miltenyibiotec.com/protocols).

▲ **Note:** Dead cells may bind non-specifically to MACS MicroBeads. To remove dead cells, we recommend using density gradient centrifugation or the Dead Cell Removal Kit (# 130-090-101).

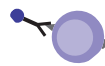

### 2.2 Magnetic labeling

▲ Work fast, keep cells cold, and use pre-cooled solutions. This will prevent capping of antibodies on the cell surface and non-specific cell labeling.

▲ Volumes for magnetic labeling given below are for up to 10<sup>7</sup> total cells. When working with fewer than 10<sup>7</sup> cells, use the same volumes as indicated. When working with higher cell numbers, scale up all reagent volumes and total volumes accordingly (e.g. for 2×10<sup>7</sup> total cells, use twice the volume of all indicated reagent volumes and total volumes).

▲ For optimal performance it is important to obtain a single-cell suspension before magnetic separation. Pass cells through 30 μm nylon mesh (Pre-Separation Filters, # 130-041-407) to remove cell clumps which may clog the column.

▲ The optimal relative centrifugal force (RCF) and centrifugation time may be different depending on the cell sample.

1. Determine cell number.
2. Centrifuge cell suspension at 300×g for 10 minutes. Aspirate supernatant completely.
3. Resuspend cell pellet and label cells with biotinylated antibody at time and titer recommended by the manufacturer. Typically, labeling for 10 minutes is sufficient.
  - ▲ **Note:** The biotinylated antibody should be used at its optimal titer, i.e., with optimal labeling intensity and no background labeling.
4. Wash cells to remove unbound primary antibody by adding 1–2 mL of **labeling buffer** per 10<sup>7</sup> cells and centrifuge at 300×g for 10 minutes. Aspirate supernatant completely.
5. (Optional) Repeat washing step.
6. Resuspend cell pellet in 90 μL of **labeling buffer** per 10<sup>7</sup> total cells.
7. Add 10 μL of Streptavidin MicroBeads per 10<sup>7</sup> total cells.
8. Mix well and refrigerate for 15 minutes (4–8 °C).
  - ▲ **Note:** Working on ice requires increased incubation time. Higher temperatures and/or longer incubation times lead to non-specific cell labeling.
9. (Optional) Add 10 μL of Anti-Biotin-fluorochrome per 10<sup>7</sup> cells and refrigerate for 5 minutes in the dark (4–8 °C).
10. Wash cells by adding 1–2 mL of buffer per 10<sup>7</sup> cells and centrifuge at 300×g for 10 minutes. Aspirate supernatant completely.
11. Resuspend up to 10<sup>8</sup> cells in 500 μL of **separation buffer**.
  - ▲ **Note:** For higher cell numbers, scale up buffer volume accordingly.
  - ▲ **Note:** For depletion with LD columns, resuspend cell pellet in 500 μL of buffer for up to 1.25×10<sup>8</sup> cells.
12. Proceed to magnetic separation (2.3).



## 2.3 Magnetic separation

▲ Choose an appropriate MACS® Column and MACS Separator according to the number of total cells and the number of magnetically labeled cells. For details see table in section 1.3.

### Magnetic separation with MS or LS Columns

1. Place column in the magnetic field of a suitable MACS Separator. For details see respective MACS Column data sheet.
2. Prepare column by rinsing with appropriate amount of **separation buffer**:

MS: 500 µL      LS: 3 mL

3. Apply cell suspension onto the column.
4. Collect unlabeled cells that pass through and wash column with appropriate amount of **separation buffer**. Perform washing steps by adding **separation buffer** three times. Only add new buffer when the column reservoir is empty.

MS: 3×500 µL      LS: 3×3 mL

Collect total effluent. This is the unlabeled cell fraction.

5. Remove column from the separator and place it on a suitable collection tube.
6. Pipette an appropriate amount of **separation buffer** onto the column. Immediately flush out the magnetically labeled cells by firmly pushing the plunger into the column.

MS: 1 mL      LS: 5 mL

▲ **Note:** To increase the purity of the magnetically labeled fraction pass the cells over a new, freshly prepared column.

### Magnetic separation with XS Columns

For instructions on the column assembly and the separation refer to the XS Column data sheet.

### Depletion with LD Columns

1. Place LD Column in the magnetic field of a suitable MACS Separator. For details see LD Column data sheet.
2. Prepare column by rinsing with 2 mL of **separation buffer**.
3. Apply cell suspension onto the column.
4. Collect unlabeled cells that pass through and wash with 2×1 mL of **separation buffer**. Collect total effluent. This is the unlabeled cell fraction.

### Depletion with CS Columns

1. Assemble CS Column and place it in the magnetic field of a suitable MACS Separator. For details see CS Column data sheet.
2. Prepare column by filling and rinsing with 60 mL of **separation buffer**. Attach a 22G flow resistor to the 3-way stopcock of the assembled column. For details see CS Column data sheet.
3. Apply cell suspension onto the column.
4. Collect unlabeled cells that pass through and wash with 30 mL of **separation buffer** from top. Collect total effluent. This is the unlabeled cell fraction.

### Depletion with D Columns

For instructions on column assembly and separation refer to the D Column data sheet.

## Magnetic separation with autoMACS™ Separator

▲ Refer to the autoMACS™ user manual for instructions on how to use the autoMACS Separator.


1. Prepare and prime autoMACS Separator.
2. Place tube containing magnetically labeled cells in the autoMACS Separator. Choose a separation program according to the recommendations in the autoMACS user manual.

▲ **Note:** Program choice depends on the isolation strategy, the strength of magnetic labeling, and the frequency of magnetically labeled cells. For details see autoMACS user manual, section autoMACS Cell Separation Programs.


## 3. Example of a separation using Streptavidin MicroBeads

Separation of human peripheral blood mononuclear cells (PBMCs) using biotinylated anti-human CD4 antibody, Streptavidin MicroBeads, an MS Column, and a MiniMACS™ Separator. Cells are fluorescently stained with Anti-Biotin-PE (# 130-090-756). Cell debris and dead cells were excluded from the analysis based on scatter signals and PI fluorescence.


Before separation

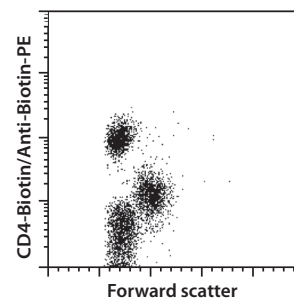

Negative fraction

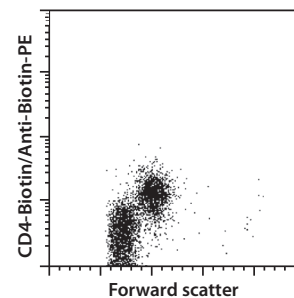

Positive fraction

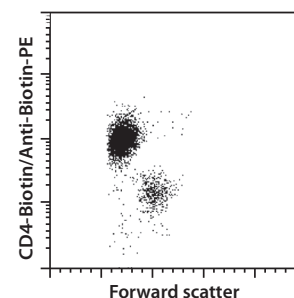

## 4. References

1. Callard, R. E. *et al.* (1995) CD40 Cross-Linking Inhibits Specific Antibody Production by Human B Cells. *Inter. Immunol.* 7: 1809–1815. [258]
2. Zürcher, A.W. *et al.* (1995) IgE-Producing Hybridomas established after B-cell culture in the CD40 system. *Immunol. Letters* 46: 49–57. [226]
3. Clarke, C. *et al.* (1994) An Immunomagnetic Separation Method Using Superparamagnetic (MACS) Beads for Large-Scale Purification of Human Mammary Luminal and Myoepithelial Cells. *Epith. Cell Biol.* 3: 38–46. [26]
4. Schmitz, J. *et al.* (1993) Optimizing Follicular Dendritic Cell Isolation by Discontinuous Gradient Centrifugation and Use of the Magnetic Cell Sorter (MACS). *J. Immunol. Methods* 159: 189–196. [106]
5. Schmitz, J. *et al.* (1994) Follicular Dendritic Cells Retain HIV-1 Particles on Their Plasma Membrane, but are not Productively Infected in Asymptomatic Patients with Follicular Hyperplasia. *J. Immunol.* 153: 1352–1359. [110]
6. Heath, S. *et al.* (1995) Follicular Dendritic Cells and Human Immunodeficiency Virus Infectivity. *Nature* 377: 740–744. [177]
7. Gronthos, S. and Simmons, P. J. (1995) The Growth Factor Requirements of STRO-1-Positive Human Bone Marrow Stromal Precursors Under Serum-Deprived Conditions *In Vitro*. *Blood* 15: 929–940. [149]
8. Klarmann, B. *et al.* (1995) Mutant Frequency at the H-2K Class I and HPRT Genes in T Lymphocytes from the X-Ray-Exposed Mouse. *Int. J. Radiat. Biol.* 67: 421–430. [140]
9. Serwe, M. and Sablitzky, F. (1993) V(D)J Recombination in B Cells is Impaired but not Blocked by Targeted Deletion of the Immunoglobulin Heavy Chain Intron Enhancer. *EMBO J.* 12: 2321–2327. [112]
10. Do H. *et al.* (1999) Expression of Factor VIII by Murine Liver Sinusoidal Endothelial Cells. *J. Biol. Chem.* 274: 19587–19592. [679]
11. Volkman, A. *et al.* (1997) Antigen-Presenting Cells in the Thymus That Can Negatively Select MHC Class II-Restricted T Cells Recognizing a Circulating Self Antigen. *J. Immunol.* 158: 693–706. [317]
12. Hsu, C. H. *et al.* (1996) Immunophylaxis of Allergen-induced Immunoglobulin E Synthesis and Airway Hyperresponsiveness in vivo by Genetic Immunization. *Nature Med.* 2: 540–544. [220]
13. Shoge, K. *et al.* (1999) Rat retinal ganglion cells culture enriched with the magnetic cell sorter. *Neuroscience Letters* 259: 111–114. [843]
14. Wallace, M. *et al.* (1994) Antiviral Activity of Primate  $\gamma\delta$  T Lymphocytes Isolated by Magnetic Cell Sorting. *J. Med. Primatol.* 23: 131–135. [276]
15. Gupta, V. K. *et al.* (1998) Two B cell subpopulations have distinct recirculation characteristics. *Eur. J. Immunol.* 28: 1597–1603. [476]
16. Koumans-van Diepen, J. C. E. *et al.* (1994) Immunocytochemical Analysis of Mitogen Responses of Carp (*Cyprinus Carpio* L.) Peripheral Blood Leucocytes. *Vet. Immunol. Immunopathol.* 42: 209–219. [79]
17. Dörr, I. *et al.* (1994) Selecting Somatic Hybrid Plants Using Magnetic Protoplast Sorting. *BIO/TECHNOLOGY* 12: 511–515. [32]

Refer to [www.miltenyibiotec.com](http://www.miltenyibiotec.com) for all data sheets and protocols. Miltenyi Biotec provides technical support worldwide. Visit [www.miltenyibiotec.com/local](http://www.miltenyibiotec.com/local) to find your nearest Miltenyi Biotec contact.

## Warnings

Reagents contain sodium azide. Under acidic conditions sodium azide yields hydrazoic acid, which is extremely toxic. Azide compounds should be diluted with running water before discarding. These precautions are recommended to avoid deposits in plumbing where explosive conditions may develop.

## Legal notices

### Limited product warranty

Miltenyi Biotec B.V. & Co. KG and/or its affiliate(s) warrant this product to be free from material defects in workmanship and materials and to conform substantially with Miltenyi Biotec's published specifications for the product at the time of order, under normal use and conditions in accordance with its applicable documentation, for a period beginning on the date of delivery of the product by Miltenyi Biotec or its authorized distributor and ending on the expiration date of the product's applicable shelf life stated on the product label, packaging or documentation (as applicable) or, in the absence thereof, ONE (1) YEAR from date of delivery ("Product Warranty"). Miltenyi Biotec's Product Warranty is provided subject to the warranty terms as set forth in Miltenyi Biotec's General Terms and Conditions for the Sale of Products and Services available on Miltenyi Biotec's website at [www.miltenyibiotec.com](http://www.miltenyibiotec.com), as in effect at the time of order ("Product Warranty"). Additional terms may apply. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS. THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING IF A PRODUCT IS SUITABLE FOR CUSTOMER'S PARTICULAR PURPOSE AND APPLICATION METHODS.

### Technical information

The technical information, data, protocols, and other statements provided by Miltenyi Biotec in this document are based on information, tests, or experience which Miltenyi Biotec believes to be reliable, but the accuracy or completeness of such information is not guaranteed. Such technical information and data are intended for persons with knowledge and technical skills sufficient to assess and apply their own informed judgment to the information. Miltenyi Biotec shall not be liable for any technical or editorial errors or omissions contained herein.

All information and specifications are subject to change without prior notice. Please contact Miltenyi Biotec Technical Support or visit [www.miltenyibiotec.com](http://www.miltenyibiotec.com) for the most up-to-date information on Miltenyi Biotec products.

### Licenses

This product and/or its use may be covered by one or more pending or issued patents and/or may have certain limitations. Certain uses may be excluded by separate terms and conditions. Please contact your local Miltenyi Biotec representative or visit Miltenyi Biotec's website at [www.miltenyibiotec.com](http://www.miltenyibiotec.com) for more information.

The purchase of this product conveys to the customer the non-transferable right to use the purchased amount of the product in research conducted by the customer (whether the customer is an academic or for-profit entity). This product may not be further sold. Additional terms and conditions (including the terms of a Limited Use Label License) may apply.

CUSTOMER'S USE OF THIS PRODUCT MAY REQUIRE ADDITIONAL LICENSES DEPENDING ON THE SPECIFIC APPLICATION. THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING FOR ITSELF WHETHER IT HAS ALL APPROPRIATE LICENSES IN PLACE. Miltenyi Biotec provides no warranty that customer's use of this product does not and will not infringe intellectual property rights owned by a third party. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

### Trademarks

autoMACS, MACS, MidiMACS, the Miltenyi Biotec logo, MiniMACS, OctoMACS, QuadroMACS, SuperMACS, and VarioMACS are registered trademarks or trademarks of Miltenyi Biotec and/or its affiliates in various countries worldwide. All other trademarks mentioned in this publication are the property of their respective owners and are used for identification purposes only.

Ficoll-Paque is a trademark of GE Healthcare companies.

Copyright © 2020 Miltenyi Biotec and/or its affiliates. All rights reserved.