



Miltenyi Biotec



CliniMACS® Formulation Unit

Automated sampling and final formulation

The CliniMACS Formulation Unit expands the processing capabilities of the CliniMACS Prodigy® Platform to include automated sampling of cells and automated formulation to deliver fresh and frozen final cell products. As an accessory of the CliniMACS Prodigy it further reduces hands-on time from sample to final cell product.

- Closed, end-to-end automated process from cell sample to cell product
- Expands automation to include fresh or frozen final formulation and fill in desired cell concentrations
- Customized solutions available to meet your needs

► miltenyibiotec.com/prodigy



The game changer in cell product formulation and final fill

Full end-to-end automation with the CliniMACS Formulation Unit avoids risky manual cell formulation steps. Uninterrupted cell processing and final formulation guarantee safety for your cells and short cuts quality control procedures.

Our Customized Application Service and process development teams are available for you to enable the translation of individual cell product manufacturing and formulation procedures into automated GMP-compliant processes. Contact us and let us bring your final cell product formulation to the next level.

CliniMACS Formulation Set

The CliniMACS Formulation Set is a single-use, biocompatible tubing module that must be connected to the CliniMACS Prodigy Tubing Set. It allows for automated sampling and contains eight EVA product bags, ready for both fresh or frozen product strategies.



Figure 1: The installation of the CliniMACS Formulation Set is guided by the user-friendly CliniMACS Prodigy touchscreen.

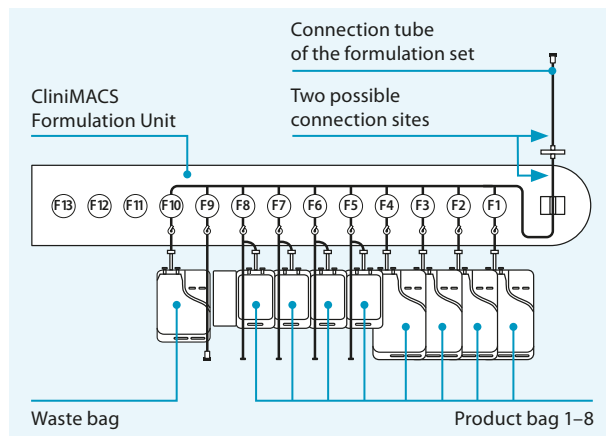


Figure 2: Schematic drawing of the CliniMACS Formulation Set, installed on the CliniMACS Formulation Unit.

CliniMACS Formulation Solution and Cryo Supplement

The CliniMACS Formulation Solution has been developed for *ex vivo* cell-based drug product formulation. As it already contains human serum albumin (HSA), it serves as a ready-to-use solution and can be directly welded to a CliniMACS Prodigy Tubing Set.

CliniMACS Cryo Supplement, 3x can be sterile-connected to the CliniMACS Formulation Set and is the perfect choice for the cryopreservation of your cells. In combination with the CliniMACS Formulation Solution it serves as a freezing solution with 7.5% DMSO final concentration.



Figure 3: CliniMACS Formulation Solution and Cryo Supplement come in one liter bags and can be sterile-connected to your process setup.



Miltenyi Biotec

Miltenyi Biotec B.V. & Co. KG | Friedrich-Ebert-Straße 68 | 51429 Bergisch Gladbach | Germany | Phone +49 2204 8306-0 | Fax +49 2204 85197
macsde@miltenyi.com | www.miltenyibiotec.com

Miltenyi Biotec provides products and services worldwide. Visit www.miltenyibiotec.com/local to find your nearest Miltenyi Biotec contact.

Unless otherwise specifically indicated, Miltenyi Biotec products and services are for research use only and not for therapeutic or diagnostic use. CliniMACS, CliniMACS Prodigy and the Miltenyi Biotec logo are registered trademarks or trademarks of Miltenyi Biotec and/or its affiliates in various countries worldwide. Copyright © 2021 Miltenyi Biotec and/or its affiliates. All rights reserved.



Ask your local Miltenyi Biotec representative for more information or visit

► miltenyibiotec.com/prodigy