Miltenyi Biotec

CD19 MicroBeads
human

Order no. 130-050-301

Contents

1. Description
   - 1.1 Principle of the MACS® Separation
   - 1.2 Background information
   - 1.3 Applications
   - 1.4 Reagent and instrument requirements
2. Protocol
   - 2.1 Sample preparation
3. Example of a separation using the CD19 MicroBeads

1. Description

This product is for research use only.

Components	2 mL CD19 MicroBeads, human : MicroBeads conjugated to monoclonal anti-human CD19 antibodies (isotype: mouse IgG1).
Capacity	For 10^9 total cells
Product format	CD19 MicroBeads are supplied in buffer containing stabilizer and 0.05% sodium azide.
Storage	Store protected from light at +2 to +8 °C. Do not freeze. The expiration date is indicated on the vial label.

1.1 Principle of the MACS Separation

First, the $CD19^+$ cells are magnetically labeled with CD19MicroBeads, human. Then, the cell suspension is loaded onto a MACS Column, which is placed in the magnetic field of a MACS Separator. The magnetically labeled $CD19^+$ cells are retained within the column. The unlabeled cells run through; this cell fraction is thus depleted of $CD19^+$ cells. After removing the column from the magnetic field, the magnetically retained $CD19^+$ cells can be eluted as the positively selected cell fraction. To increase the purity, the positively selected cell fraction containing the $CD19^+$ cells must be separated over a second column.

1.2 Background information

CD19 MicroBeads have been developed for the separation of human B cells based on the expression of the CD19 antigen. CD19 is a 95–120 kDa glycosylated transmembrane protein that is critically involved in signal transduction processes that regulate development, activation, and differentiation of B lymphocytes. CD19 is expressed on lineage B cells from the early lineage-committed pro-B cell stage to the B cell blast stage including most malignant B cells. Expression is down-regulated during the differentiation into plasma cells. Furthermore, CD19 is expressed on follicular dendritic cells.

1.3 Applications

- Positive selection or depletion of cells expressing human CD19 antigen.
- Isolation or depletion of B cells from peripheral blood mononuclear cells (PBMCs), bodily fluids (e.g. bronchial lavage), or single-cell suspensions from tissue (e.g. lymphoid and tumor tissue).

1.4 Reagent and instrument requirements

- Buffer: Prepare a solution containing phosphate-buffered saline (PBS), pH 7.2, 0.5% bovine serum albumin (BSA), and 2 mM EDTA by diluting MACS BSA Stock Solution (# 130-091-376) 1:20 with autoMACS® Rinsing Solution (# 130-091-222). Keep buffer cold (+2 to +8 °C). Degas buffer before use, as air bubbles could block the column.

  ▲ **Note:** EDTA can be replaced by other supplements such as anticoagulant citrate dextrose formula-A (ACD-A) or citrate phosphate dextrose (CPD). BSA can be replaced by other proteins such as human serum albumin, human serum, or fetal bovine serum (FBS). Buffers or media containing Ca^{2+} or Mg^{2+} are not recommended for use.

- Table 1 provides an overview of MACS Columns and MACS Separators with links to the corresponding complete protocols. Use the link matching the MACS Column and MACS Separator used.

Table 1: MACS Columns and MACS Separators overview.

MACS Column (link to protocol)	Max. number of labeled cells	Max. number of total cells	MACS Separator
Columns for manual separation			
MS Column	10^7	2×10^8	MiniMACS, OctoMACS
LS Column	10^8	2×10^9	MidiMACS, QuadroMACS
XLS Column	5×10^8	10^{10}	QuadroMACS XL
LD Column	10^8	5×10^8	MidiMACS, QuadroMACS
Columns for instruments			
LS Column	10^8	10^9	MultiMACS Cell24 Separator Plus
autoMACS Column	2×10^8	4×10^9	autoMACS NEO
autoMACS Column	2×10^8	4×10^9	autoMACS Pro Separation program: Positive selection: Possel Depletion: Depletes

Miltenyi Biotec B.V. & Co. KG
Friedrich-Ebert-Straße 68, 51429 Bergisch Gladbach, Germany
Phone +49 2204 8306-0, Fax +49 2204 85197
macsde@miltenyi.com
www.miltenyibiotec.com

140-000-068.09

MACS Column (link to protocol)	Max. number of labeled cells	Max. number of total cells	MACS Separator
Multi-24 Column Block	10^8 per column 2.4×10^9 per block	10^9 per column 2.4×10^{10} per block	MultiMACS Cell24 Separator Plus

▲ **Note:** If separating with LS or LD Columns and the MultiMACS Cell24 Separator Plus, use the Single-Column Adapter. Refer to the user manual for details.

▲ **Note:** Cells which strongly express the CD19 antigen can also be depleted using MS, LS, or XLS Columns.

▲ **Note:** For autolabeling compatibilty for the autoMACS NEO Seperator, please see here.

- (Optional) Fluorochrome-conjugated antibodies for flow cytometric analysis, e.g., CD19 Antibody, anti-human. For more information about antibodies refer to www.miltenyibiotec.com/antibodies.
- (Optional) Propidium Iodide Solution (# 130-093-233) or 7-AAD Staining Solution (# 130-111-568) for flow cytometric exclusion of dead cells.
- (Optional) Dead Cell Removal Kit (# 130-090-101) for the depletion of dead cells.
- (Optional) Pre-Separation Filters (30 µm) (# 130-041-407) to remove cell clumps.

2. Protocol

2.1 Sample preparation

When working with anticoagulated peripheral blood or buffy coat, peripheral blood mononuclear cells (PBMCs) should be isolated by using a MACS PBMC Isolation Kit or by density gradient centrifugation, for example, using Ficoll-Paque™.


▲ **Note:** To remove platelets after density gradient separation, resuspend cell pellet in buffer and centrifuge at 200 × g for 10–15 minutes at +20 °C. Carefully aspirate supernatant. Repeat washing step.

For details refer to the data sheet or the protocols section at www.miltenyibiotec.com/protocol.

▲ Dead cells may bind non-specifically to MACS MicroBeads. To remove dead cells, it is recommended to use density gradient centrifugation or the Dead Cell Removal Kit (# 130-090-101).

3. Example of a separation using CD19 MicroBeads, human

$CD19^+$ cells were isolated from human PBMCs using CD19 MicroBeads, human, an MS Column, and a MiniMACS™ Separator. Cells were fluorescently stained with CD19-PE and analyzed by flow cytometry. Cell debris and dead cells were excluded from the analysis based on scatter signals and propidium iodide fluorescence.

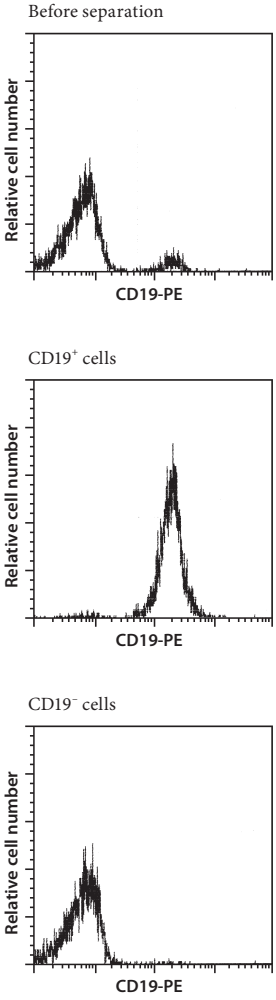

Refer to **www.miltenyibiotec.com** for all data sheets and protocols. Miltenyi Biotec provides technical support worldwide. Visit www.miltenyibiotec.com for local Miltenyi Biotec Technical Support contact information.

Legal notices

Limited product warranty

Miltenyi Biotec B.V. & Co. KG and/or its affiliate(s) warrant this product to be free from material defects in workmanship and materials and to conform substantially with Miltenyi Biotec's published specifications for the product at the time of order, under normal use and conditions in accordance with its applicable documentation, for a period beginning on the date of delivery of the product by Miltenyi Biotec or its authorized distributor and ending on the expiration date of the product's applicable shelf life stated on the product label, packaging or documentation (as applicable) or, in the absence thereof, ONE (1) YEAR from date of delivery ("Product Warranty"). Miltenyi Biotec's Product Warranty is provided subject to the warranty terms as set forth in Miltenyi Biotec's General Terms and Conditions for the Sale of Products and Services available on Miltenyi Biotec's website at www.miltenyibiotec.com, as in effect at the time of order ("Product Warranty"). Additional terms may apply. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING IF A PRODUCT IS SUITABLE FOR CUSTOMER'S PARTICULAR PURPOSE AND APPLICATION METHODS.

Technical information

The technical information, data, protocols, and other statements provided by Miltenyi Biotec in this document are based on information, tests, or experience which Miltenyi Biotec believes to be reliable, but the accuracy or completeness of such information is not guaranteed. Such technical information and data are intended for persons with knowledge and technical skills sufficient to assess and apply their own informed judgment to the information. Miltenyi Biotec shall not be liable for any technical or editorial errors or omissions contained herein.

All information and specifications are subject to change without prior notice. Please contact Miltenyi Biotec Technical Support or visit www.miltenyibiotec.com for the most up-to-date information on Miltenyi Biotec products.

Licenses

This product and/or its use may be covered by one or more pending or issued patents and/or may have certain limitations. Certain uses may be excluded by separate terms and conditions. Please contact your local Miltenyi Biotec representative or visit Miltenyi Biotec's website at www.miltenyibiotec.com for more information.

The purchase of this product conveys to the customer the non-transferable right to use the purchased amount of the product in research conducted by the customer (whether the customer is an academic or for-profit entity). This product may not be further sold. Additional terms and conditions (including the terms of a Limited Use Label License) may apply.

CUSTOMER'S USE OF THIS PRODUCT MAY REQUIRE ADDITIONAL LICENSES DEPENDING ON THE SPECIFIC APPLICATION. THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING FOR ITSELF WHETHER IT HAS ALL APPROPRIATE LICENSES IN PLACE. Miltenyi Biotec provides no warranty that customer's use of this product does not and will not infringe intellectual property rights owned by a third party. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

Trademarks

autoMACS, MACS, MidiMACS, the Miltenyi Biotec logo, MidiMACS, MultiMACS, OctoMACS, and QuadroMACS are registered trademarks or trademarks of Miltenyi Biotec B.V. & Co. KG and/or its affiliates in various countries worldwide.

Ficoll-Paque is a trademark of GE Healthcare companies

Copyright © 2026 Miltenyi Biotec and/or its affiliates. All rights reserved.