

Contents

1. Description
 - 1.1 Principle of the REAl ease MACS[®] Separation
 - 1.2 Background information
 - 1.3 Applications
 - 1.4 Reagent and instrument requirements
2. Protocol
 - 2.1 Protocol overview
 - 2.2 Sample preparation
 - 2.3 Magnetic labeling
 - 2.4 Magnetic separation and removal of magnetic labeling
 - 2.4.1 Highest recovery
 - 2.4.2 Highest purity
 - 2.5 Magnetic labeling and separation using autoMACS[®] Separators
 - 2.6 (Optional) Removal of the REAl ease Complex and second magnetic labeling with REAl ease MicroBeads
 - 2.6.1 Removal of the REAl ease Complex
 - 2.6.2 Second magnetic labeling with REAl ease MicroBeads
 - 2.7 (Optional) Second magnetic labeling with MACS MicroBeads
3. Example of a separation using the REAl ease CD90.1 MicroBead Kit

1. Description

This product is for research use only.

Components	2×1 mL REAl ease CD90.1-Biotin, mouse 2×5 mL REAl ease Anti-Biotin MicroBeads (CD90.1, mouse) 2×4 mL REAl ease Bead Release Reagent (50×) 4 mL REAl ease Release Reagent 2×4 mL REAl ease Stop Reagent
Capacity	For 2×10 ⁹ total cells, up to 100 separations.
Product format	REAl ease Stop Reagent is supplied in buffer containing 0.05% sodium azide. All other reagents are supplied in buffer containing stabilizer and 0.05% sodium azide.
Storage	Store protected from light at 2–8 °C. Do not freeze. The expiration date is indicated on the vial label.

1.1 Principle of the REAl ease MACS Separation

The REAl ease Technology relies on recombinantly engineered antibody fragments to label specific cell surface markers. The fragments are developed to have low affinity for epitopes. However, when fragments are multimerized as a REAl ease Biotin Complex (i.e., REAl ease CD90.1-Biotin, mouse) they bind to epitopes with high avidity. REAl ease Technology can control the multimer/monomer state of antibody fragments. With this technology a controlled label release is possible where monomerized antibody fragments dissociate from the cell surface, enabling users to obtain bead- and label-free cells.

First, the target cells are labeled with REAl ease CD90.1-Biotin (REAl ease Biotin Complex). Subsequently, REAl ease Anti-Biotin MicroBeads (CD90.1, mouse) bind to the REAl ease Biotin Complex. Then, the cell suspension is loaded onto a MACS Column, which is placed in the magnetic field of a MACS Separator. The magnetically labeled cells are retained within the column. The unlabeled non-target cells flow through; this cell fraction is thus depleted of CD90.1⁺ cells. Target cells are eluted from the column using the REAl ease Bead Release Reagent, which simultaneously removes the MicroBeads from the cells. For highest recovery, the column is removed from the magnetic field and target cells are eluted. For highest purity, target cells are eluted within the magnetic field of the column. Finally, during the subsequent incubation with the REAl ease Release Reagent, the REAl ease Biotin Complex monomerizes and dissociates from the cell surface leaving the cells free of all labels.

1.2 Background information

The REAl ease CD90.1 MicroBead Kit, mouse has been developed for positive selection of mouse CD90.1⁺ cells from single-cell suspensions of lymphoid and non-lymphoid tissues. The mouse monoclonal antibody reacts with mouse CD90.1 (Thy1.1), a GPI-anchored conserved membrane glycoprotein. In the mouse strains AKR/J, PL, and FVB/N, CD90.1 is a pan T cell marker and can be found on thymocytes, hematopoietic stem cells in the bone marrow, intraepithelial cells (dendritic epidermal cells) in skin, and on neurons, such as retinal ganglion cells. The antibody does not cross-react with mouse CD90.2 (Thy1.2), but cross-reacts with rat CD90 antigen (expressed on thymocytes, hematopoietic stem cells, neurons, and other cell types).

1.3 Applications

- Positive selection of mouse CD90.1⁺ cells from single-cell suspensions, e.g., from spleen or lung preparations.
- Isolation of cells for the adoptive transfer of CD90.1⁺ T cells into CD90.2⁺ inbred strains or immunodeficient mice.
- Isolation of CD90.1⁺ label-free cells.

1.4 Reagent and instrument requirements

- Separation buffer: Prepare a solution containing phosphate-buffered saline (PBS), pH 7.2, 0.5% bovine serum albumin (BSA), and 2 mM EDTA by diluting MACS BSA Stock Solution (#130-091-376) 1:20 with autoMACS Rinsing Solution (#130-091-222). Use buffer at room temperature (19–25 °C). Store buffer cold (2–8 °C). Degas buffer before use, as air bubbles could block the column.

▲ **Note:** BSA can be replaced by other proteins such as mouse serum albumin, mouse serum, or fetal bovine serum (FBS). Buffers or media containing Ca^{2+} or Mg^{2+} are not recommended for use.

- REAlase Bead Release buffer: Prepare a 1:50 dilution of REAlase Bead Release Reagent (50×), e.g., for 1 mL add 20 µL of REAlase Bead Release Reagent to 980 µL of separation buffer.

▲ **Note:** Use freshly prepared buffer the same day. Store at room temperature.

▲ **Note:** For highest recovery protocol, prepare 1 mL per MS Column and 5 mL per LS Column. For highest purity protocol, prepare 3 mL per MS Column and 14 mL per LS Column.

- MACS Columns and MACS Separators: CD90.1⁺ cells can be enriched by using MS or LS Columns. Positive selection can also be performed by using autoMACS Columns on the autoMACS NEO Separator.

Column	Max. number of labeled cells	Max. number of total cells	Separator
Positive selection			
MS	10 ⁷	2×10 ⁸	MiniMACS, OctoMACS, SuperMACS II
LS	10 ⁸	2×10 ⁹	MidiMACS, QuadroMACS, SuperMACS II
autoMACS	2×10 ⁸	4×10 ⁹	autoMACS NEO Separator


▲ **Note:** Column adapters are required to insert certain columns into SuperMACS™ II Separators. For details refer to the respective MACS Separator data sheet.

- (Optional) Fluorochrome-conjugated antibodies for flow cytometric analysis, e.g., CD45 Antibody, anti-mouse, VioGreen™ and CD90.1 Antibody, anti-mouse, FITC. For more information about antibodies refer to www.miltenyibiotec.com/antibodies.
- (Optional) Propidium Iodide Solution (#130-093-233), 7-AAD Staining Solution (#130-111-568), or DAPI Staining Solution (#130-111-570) for flow cytometric exclusion of dead cells.
- (Optional) Dead Cell Removal Kit (#130-090-101) for the depletion of dead cells.
- (Optional) Pre-Separation Filters (30 µm) (#130-041-407) to remove cell clumps.


2. Protocol

2.1 Protocol overview

A) Highest recovery

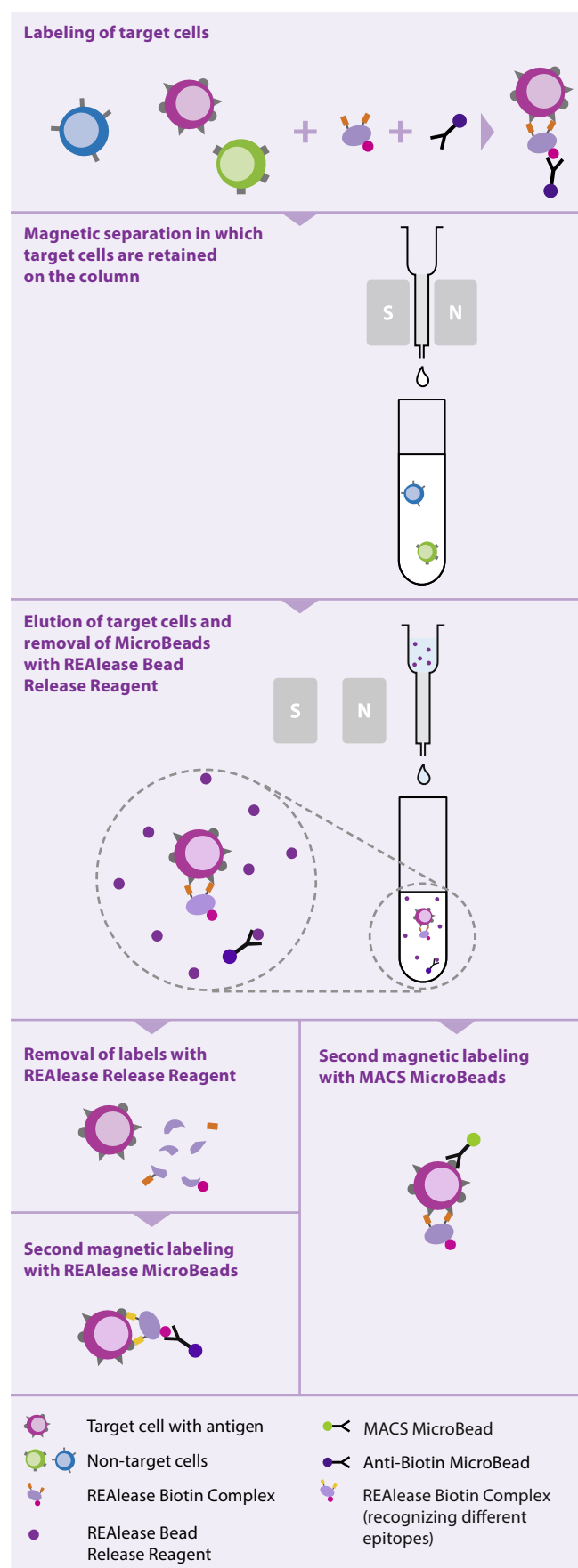

B) Highest purity

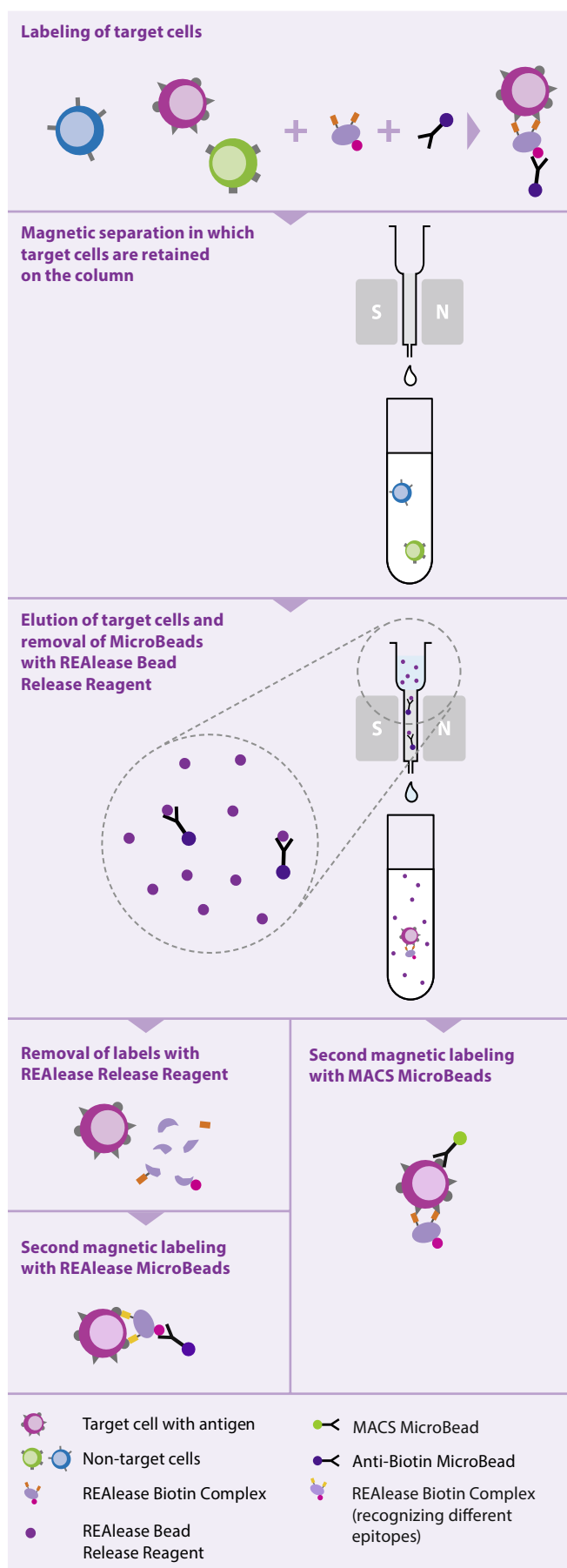

2.2 Sample preparation

When working with tissues, prepare a single-cell suspension using the gentleMACS™ Dissociator.

For details refer to www.miltenyibiotec.com/gentlemacs.



2.3 Magnetic labeling

▲ Dead cells may bind non-specifically to MACS MicroBeads. To remove dead cells, it is recommended to use density gradient centrifugation or the Dead Cell Removal Kit (# 130-090-101).

▲ The recommended incubation temperature is at room temperature (19–25 °C).

▲ Volumes for magnetic labeling given below are for up to 2×10^7 total cells. When working with fewer than 2×10^7 cells, use the same volumes as indicated. When working with higher cell numbers, scale up all reagent volumes and total volumes accordingly (e.g. for 4×10^7 total cells, use twice the volume of all indicated reagent volumes and total volumes).

▲ For optimal performance it is important to obtain a single-cell suspension before magnetic labeling. Pass cells through 30 μ m nylon mesh (Pre-Separation Filters (30 μ m), # 130-041-407) to remove cell clumps which may clog the column. Moisten filter with buffer before use.

1. Determine cell number.
2. Centrifuge cell suspension at $300 \times g$ for 10 minutes. Aspirate supernatant completely.
3. Resuspend cell pellet in 80 μ L of separation buffer per 2×10^7 total cells.
4. Add 20 μ L of REAlease CD90.1-Biotin per 2×10^7 total cells.
5. Mix well and incubate for 5 minutes.
6. Add 100 μ L of REAlease Anti-Biotin MicroBeads (CD90.1, mouse) per 2×10^7 total cells.
7. Mix well and incubate for 5 minutes.
8. (Optional) Add staining antibodies, e.g., CD45 Antibody, anti-mouse, VioGreen and CD90.1 Antibody, anti-mouse, FITC, and incubate according to manufacturer's recommendation.
▲ **Note:** These staining antibodies cannot be removed from the cells.
9. Dilute up to 5×10^7 cells in a total volume of 500 μ L with separation buffer.
▲ **Note:** For volumes larger than 500 μ L a dilution is not needed.
10. Proceed to magnetic separation (2.4).



2.4 Magnetic separation and removal of magnetic labeling

▲ Choose an appropriate MACS Column and MACS Separator according to the number of total cells and the number of CD90.1⁺ cells. For details refer to the table in section 1.4.

▲ Always wait until the column reservoir is empty before proceeding to the next step.

▲ The recommended incubation temperature is at room temperature (19–25 °C).

Magnetic separation with MS or LS Columns

- Place column in the magnetic field of a suitable MACS Separator. For details refer to the respective MACS Column data sheet.
- Prepare column by rinsing with the appropriate amount of separation buffer:

MS: 500 µL LS: 3 mL

- Apply cell suspension onto the column. Collect flow-through containing unlabeled cells.
- Wash column with the appropriate amount of separation buffer. Collect unlabeled cells that pass through and combine with the flow-through from step 3.

MS: 3×500 µL LS: 3×3 mL

▲ **Note:** Perform washing steps by adding buffer aliquots as soon as the column reservoir is empty.

- Proceed either to 2.4.1 if highest recovery is needed or 2.4.2 if highest purity is needed.

2.4.1 Highest recovery

- Remove column from the separator and place it on a suitable collection tube.
- (Optional) If removal of MicroBeads is not required, pipette appropriate amount of separation buffer. Immediately flush out target cells by firmly pushing the plunger into the column. Eluted cells are ready for downstream applications, e.g. flow cytometry analysis.

MS: 1 mL LS: 5 mL

- For removal of MicroBeads proceed with step 9.
- Pipette the appropriate amount of REAlease Bead Release buffer (prepared by REAlease Bead Release Reagent (50×), refer to chapter 1.4) onto the column. Immediately flush out target cells by firmly pushing the plunger into the column.

MS: 1 mL LS: 5 mL

- Mix well and incubate for 10 minutes.
- Cells are now free from MicroBeads and ready for analysis and downstream applications.
- (Optional) Proceed either to
 - 2.6 Removal of REAlease Complex and second magnetic labeling with REAlease MicroBeads or proceed to
 - 2.7 Second magnetic labeling with MACS MicroBeads.

2.4.2 Highest purity

- Place a new tube under the column to collect target cells
- Pipette the appropriate amount of REAlease Bead Release buffer (prepared by REAlease Bead Release Reagent (50×), refer to chapter 1.4) onto the column. Collect flow-through containing target cells.

MS: 3×1 mL LS: 2×4 mL + 1×6 mL

- Mix well and incubate for 10 minutes.
- Cells are now free from MicroBeads and ready for analysis and downstream applications.
- (Optional) Remove column from the separator and place it on a suitable collection tube. Apply separation buffer onto the column. Immediately flush out cells sticking on the column by firmly pushing the plunger into the column.

MS: 1 mL LS: 5 mL

- (Optional) Proceed either to
 - 2.6 Removal of REAlease Complex and second magnetic labeling with REAlease MicroBeads or proceed to
 - 2.7 Second magnetic labeling with MACS MicroBeads.

2.5 Magnetic labeling and separation using autoMACS Separators

▲ Refer to the user manual and the short instructions for instructions on how to use the autoMACS Separators.

▲ Buffers used for operating the autoMACS Separators should have a temperature of ≥10 °C.

▲ Place tubes in the following Chill Rack positions:

position A = sample, position B = unlabeled (negative) fraction,

position C = labeled (positive) fraction.

2.5.1 Magnetic labeling and separation using the autoMACS NEO Separator

▲ The autoMACS NEO Separator enables stage loading to extend column capacity for selected reagents, minimizing the need to divide larger samples.

▲ For more information on selecting alternative separation programs, stage loading-compatible reagents, autolabeling-compatible reagents, and the minimal and maximal volumes for each reagent and Chill Rack, refer to www.miltenyibiotec.com/automacs-neo-sample-processing.

2.6 (Optional) Removal of the REAlease Complex and second magnetic labeling with REAlease MicroBeads

▲ The recommended incubation temperature is at room temperature (19–25 °C).

▲ For second magnetic labeling with MACS Anti-Biotin MicroBeads proceed through all steps of chapter 2.6 to remove the REAlease Biotin Complex.

2.6.1 Removal of the REAlease Complex

1. Centrifuge cell suspension at 300×g for 10 minutes. Aspirate supernatant completely.
2. Resuspend cell pellet in appropriate amount of separation buffer:
MS: 1 mL LS: 5 mL
3. Add an appropriate amount of REAlease Release Reagent:
MS: 20 µL LS: 100 µL
4. Mix well and incubate for 5 minutes.
5. Cells are now free from REAlease Complex and MicroBeads and are ready for analysis or downstream applications.
6. (Optional) For second magnetic labeling with REAlease MicroBeads continue with 2.6.2.

2.6.2 Second magnetic labeling with REAlease MicroBeads

1. Centrifuge cell suspension at 300×g for 10 minutes. Aspirate supernatant completely.
2. Resuspend cell pellet in 40 µL of REAlease Stop Reagent per 10⁷ total cells.
3. Mix well.
4. For a second magnetic labeling follow the labeling protocol in the respective REAlease MicroBead Kit data sheet.

▲ **Note:** For best recovery and purity of cells, the amount of MACS MicroBeads for the second positive labeling may need optimization as the starting frequency of target cells may be different from a PBMC sample.

2.7 (Optional) Second magnetic labeling with MACS MicroBeads

▲ For second magnetic labeling with MACS Anti-Biotin MicroBeads proceed through all steps of chapter 2.6 to remove the REAlease Biotin Complex.


1. Centrifuge cell suspension at 300×g for 10 minutes. Aspirate supernatant completely.
2. Add the recommended amount of MACS MicroBeads to label the cells magnetically for the second marker. For details refer to the respective MACS MicroBeads data sheet.

▲ **Note:** For best recovery and purity of cells, the amount of MACS MicroBeads for the second positive labeling may need optimization as the starting frequency of target cells may be different from a PBMC sample.

3. Example of a separation using the REAlease CD90.1 MicroBead Kit

A perfused mouse lung sample was dissociated using the gentleMACS Octo Dissociator with Heaters in combination with the Lung Dissociation Kit, mouse (# 130-095-927). CD90.1⁺ cells were isolated from the single-cell suspension using the REAlease CD90.1 MicroBead Kit, MS Columns, and an OctoMACS™ Separator using the highest recovery protocol. Cells were fluorescently stained with CD45-VioGreen and CD90.1-FITC and analyzed by flow cytometry using the MACSQuant® Analyzer X. Cell debris and dead cells were excluded from the analysis based on scatter signals and DAPI fluorescence.

A) Cell purity



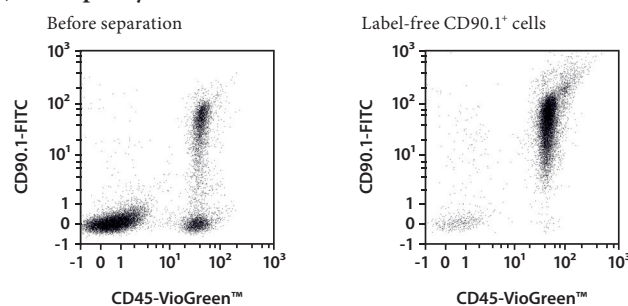


B) Bead-free cells: efficiency of REAlease Anti-Biotin MicroBeads release

Release efficiency was higher than 99% for the REAlease Anti-Biotin MicroBeads (CD90.1). The efficiency was determined by re-applying the isolated cells to a second MACS Column. The ratio between the numbers of cells in the flow-through and the total number of cells applied to the second column allowed us to calculate the efficiency of magnetic labeling removal.

C) Label-free cells: REAlease Biotin Complex release

The efficient removal of all labels was shown by using Anti-Biotin-APC to analyze the cells by flow cytometry for the presence of REAlease Biotin Complex. Directly after isolation, the cells showed staining of biotin ("MicroBead-free CD90.1⁺ cells"), whereas the label-free CD90.1⁺ cells after the REAlease Biotin Complex release were negative for biotin similar to the non-labeled cells before separation.

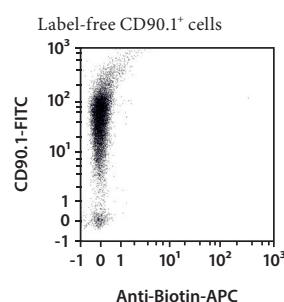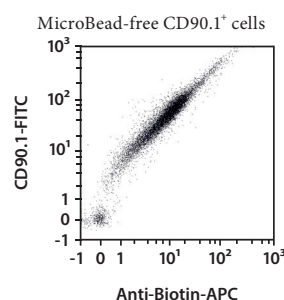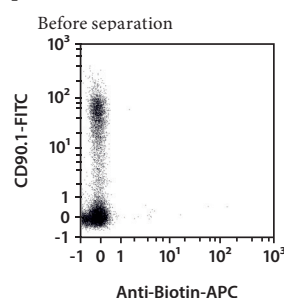

Refer to www.miltenyibiotec.com for all data sheets and protocols. Miltenyi Biotec provides technical support worldwide. Visit www.miltenyibiotec.com for local Miltenyi Biotec Technical Support contact information.

Legal notices

Limited product warranty

Miltenyi Biotec B.V. & Co. KG and/or its affiliate(s) warrant this product to be free from material defects in workmanship and materials and to conform substantially with Miltenyi Biotec's published specifications for the product at the time of order, under normal use and conditions in accordance with its applicable documentation, for a period beginning on the date of delivery of the product by Miltenyi Biotec or its authorized distributor and ending on the expiration date of the product's applicable shelf life stated on the product label, packaging or documentation (as applicable) or, in the absence thereof, ONE (1) YEAR from date of delivery ("Product Warranty"). Miltenyi Biotec's Product Warranty is provided subject to the warranty terms as set forth in Miltenyi Biotec's General Terms and Conditions for the Sale of Products and Services available on Miltenyi Biotec's website at www.miltenyibiotec.com, as in effect at the time of order ("Product Warranty"). Additional terms may apply. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING IF A PRODUCT IS SUITABLE FOR CUSTOMER'S PARTICULAR PURPOSE AND APPLICATION METHODS.

Technical information

The technical information, data, protocols, and other statements provided by Miltenyi Biotec in this document are based on information, tests, or experience which Miltenyi Biotec believes to be reliable, but the accuracy or completeness of such information is not guaranteed. Such technical information and data are intended for persons with knowledge and technical skills sufficient to assess and apply their own informed judgment to the information. Miltenyi Biotec shall not be liable for any technical or editorial errors or omissions contained herein.

All information and specifications are subject to change without prior notice. Please contact Miltenyi Biotec Technical Support or visit www.miltenyibiotec.com for the most up-to-date information on Miltenyi Biotec products.

Licenses

This product and/or its use may be covered by one or more pending or issued patents and/or may have certain limitations. Certain uses may be excluded by separate terms and conditions. Please contact your local Miltenyi Biotec representative or visit Miltenyi Biotec's website at www.miltenyibiotec.com for more information.

The purchase of this product conveys to the customer the non-transferable right to use the purchased amount of the product in research conducted by the customer (whether the customer is an academic or for-profit entity). This product may not be further sold. Additional terms and conditions (including the terms of a Limited Use Label License) may apply.

CUSTOMER'S USE OF THIS PRODUCT MAY REQUIRE ADDITIONAL LICENSES DEPENDING ON THE SPECIFIC APPLICATION. THE CUSTOMER IS SOLELY RESPONSIBLE FOR DETERMINING FOR ITSELF WHETHER IT HAS ALL APPROPRIATE LICENSES IN PLACE. Miltenyi Biotec provides no warranty that customer's use of this product does not and will not infringe intellectual property rights owned by a third party. BY USE OF THIS PRODUCT, THE CUSTOMER AGREES TO BE BOUND BY THESE TERMS.

Trademarks

autoMACS, MACS, MACSQuant, MidiMACS, the Miltenyi Biotec logo, MiniMACS, OctoMACS, QuadroMACS, REAlease, SuperMACS, and Viability are registered trademarks or trademarks of Miltenyi Biotec B.V. & Co. KG and/or its affiliates in various countries worldwide. All other trademarks mentioned in this publication are the property of their respective owners and are used for identification purposes only.

Copyright © 2025 Miltenyi Biotec and/or its affiliates. All rights reserved.